Table of Contents

README	1.1
INSTALL	1.2
Users Guide	1.3
firmware	1.3.1
Decode	1.3.2
GUI	1.3.3
OnlineStatics	1.3.4
MakeEvent	1.3.5
FrontPanel	1.3.6
Logic	1.3.7
MultipleModulesSynchronously	1.3.8
Programmers Guide	1.4
XIA API	1.4.1
PKU code	1.4.2

README

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<http://wuhongyi.cn/PKUXIADAQ/>

- Have any comments or suggestions for this project(function addition or improvement), please send an email to Hongyi Wu(wuhongyi@qq.com). 对本获取程序有任何的意见及建议(功能添加及改进) , 欢迎给吴鸿毅(wuhongyi@qq.com)发邮件。
- We will improve the Chinese/English manuals of the software as soon as possible. Currently, we mainly use the software to explain the use of the software. 我们将会尽快完善软件的中英文使用说明书 , 当前主要以操作演示来讲解软件的使用。

About/关于

This manual applies only to XIA LLC Pixie-16

- This program is developed by the **group of experimental nuclear physics of Peking University**. 本程序由北京大学实验核物理组开发。

- The earliest graphical interface development of this program is based on NSCL DDAS Nscope(<http://docs.nscl.msu.edu/daq/newsite/ddas-1.1/nscope.html>). 最早的图形界面程序是基于 NSCL DDAS Nscope 开发。
- Thanks to Hui Tan's(XIA LLC) support for our development. 特别感谢 谭辉(XIA LLC) 对我们开发的支持。

Technical adviser 技术指导:

- Zhihuan Li 李智焕
- Hui Tan 谭辉(XIA LLC)

Software Developer 软件主要开发者:

- Hongyi Wu 吴鸿毅(wuhongyi@qq.com)
- Jing Li 李晶(lijinger02@126.com)

Principal author of the instruction 说明书主要撰写者:

- Hongyi Wu 吴鸿毅
- Diwen Luo 罗迪雯
- Xiang Wang 王翔

The development of this program is supported by the following 本程序的开发得到以下单位的支持 :

- XIA LLC
 - Institute of Modern Physics, Chinese Academy of Sciences 中国科学院近代物理研究所(IMP)
 - China Institute of Atomic Energy 中国原子能科学研究院(CIAE)
 - The University of Hong Kong 香港大学(HKU)
 - Shandong University, Weihai 山东大学(威海)(SDU)
 - ...
-

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contents / 目录

The following files/folders are included in the user's use package 用户使用程序包中包含以下文件/文件夹:

- Decode(Convert raw binary data to ROOT file / 将原始二进制数据转为 ROOT)
- docs(Web page manual / 使用说明书，网页版)
- firmware(固件)
 - firmware/firmware.md(History of firmware / 历史各版本固件说明)
- GUI(图形软件)
- MakeEvent(事件重构程序，可选)
- NOGUI(Non-graphics software. New version upgrade, temporarily unavailable / 非图形软件。新版本升级中，暂时不可用)
- OnlineStatics(Online monitoring program / 在线监视程序)
- parset(parameter setting file / 参数设置文件)

- PlxSdk.tar.gz(Plx9054 driver)
 - README(markdown version manual)
 - README.md(home page introduction)
 - README.pdf(pdf version manual)
 - software(unofficial standards pixie16 driver API, have been revised by Hongyi Wu)
 - TestTool(Developer testing tool, users don't need / 开发者测试工具，用户不需要！！！)
-

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- The main control interface developed based on the ROOT GUI is highly complex and difficult for users to modify now. Users do not easy to develop their own version based on this program. 当前基于 ROOT GUI 开发的主控制界面复杂度高，用户修改难度大。其它用户不容易基于其发展适合自己的版本。
- We are also developing a web control for online/offline analysis 我们也在开发基于网页控制的获取在线/离线分析程序：
 - ZeroMQ
 - FastCGI
 - JSROOT
 - web
 - ...

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User's Manual Digital Gamma Finder (DGF) PIXIE-16 Version 1.40, October 2009

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Do Not Hot-Swap!

To avoid personal injury, and/or damage to the DGF-Pixie-16, always turn off crate power before removing the DGF-Pixie-16 from the crate!

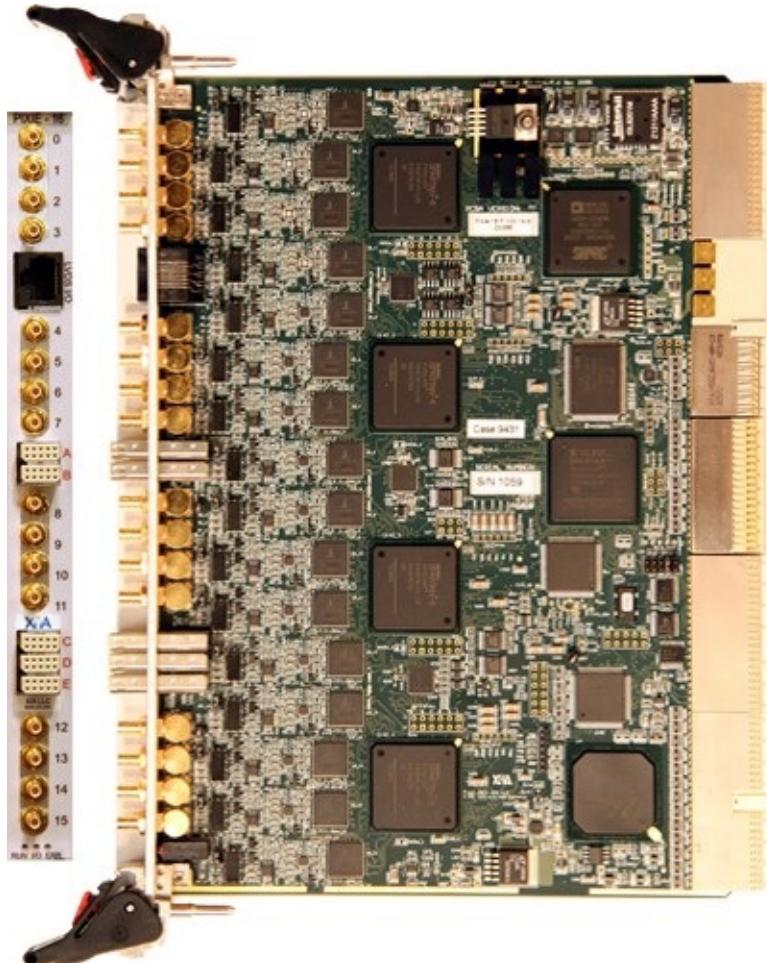
Figure: pixie16

The DGF Pixie-16 is a 16-channel all-digital waveform acquisition and spectrometer card based on the CompactPCI/PXI standard for fast data readout to the host. It combines spectroscopy with waveform digitizing and the option of on-line pulse shape analysis. The Pixie-16 accepts signals from virtually any radiation detector. Incoming signals are digitized by 12/14/16-bit 100/250/500 MSPS ADCs. Waveforms of up to 163.8 μ s in length for each channel can be stored in a FIFO.

The waveforms are available for onboard pulse shape analysis, which can be customized by adding user functions to the core processing software. Waveforms, timestamps, and the results of the pulse shape analysis can be read out by the host system for further off-line processing. Pulse heights are calculated to 16-bit precision and can be binned into spectra with up to 32K channels. The Pixie-16 supports coincidence spectroscopy and can recognize complex hit patterns.

Data readout rates through the CompactPCI/PXI backplane to the host computer can be up to 109 Mbyte/s. The standard PXI backplane, as well as additional custom backplane connections are used to distribute clocks and trigger signals between several Pixie-16 modules for group operation. A complete data acquisition and processing systems can be built by combining Pixie-16 modules with commercially available CompactPCI/PXI processor, controller or I/O modules in the same chassis.

The Pixie-16 is an instrument for waveform acquisition and MCA histogramming for arrays of gamma ray or other radiation detectors such as

- 100 MSPS
 - Segmented HPGe detectors.
 - Scintillator/PMT combinations: NaI, CsI, BGO and many others.
 - Gas detector.
 - Silicon strip detectors.
- 250 MSPS
 - Scintillator
 - LaBr₃
- 500 MSPS
 - Scintillator
 - LaBr₃

The Pixie-16 modules must be operated in a custom 6U CompactPCI/PXI chassis providing high currents at specific voltages not included in the CompactPCI/PXI standard 1 . Currently XIA provides a 14-slot chassis. Put the host computer(or remote PXI controller) in the system slot (slot 1) of your chassis. Put the Pixie-16 modules into any free peripheral slot (slot 2-14) with the chassis still powered down. After modules are installed, power up the chassis (Pixie-16 modules are not hot swappable). If using a remote controller, be sure to boot the host computer after powering up the chassis.

Clock distribution

In a multi-module system there will be one clock master and a number of clock slaves or repeaters. The clock function of a module can be selected by setting shunts on Jumper JP101 at the rear of the board.

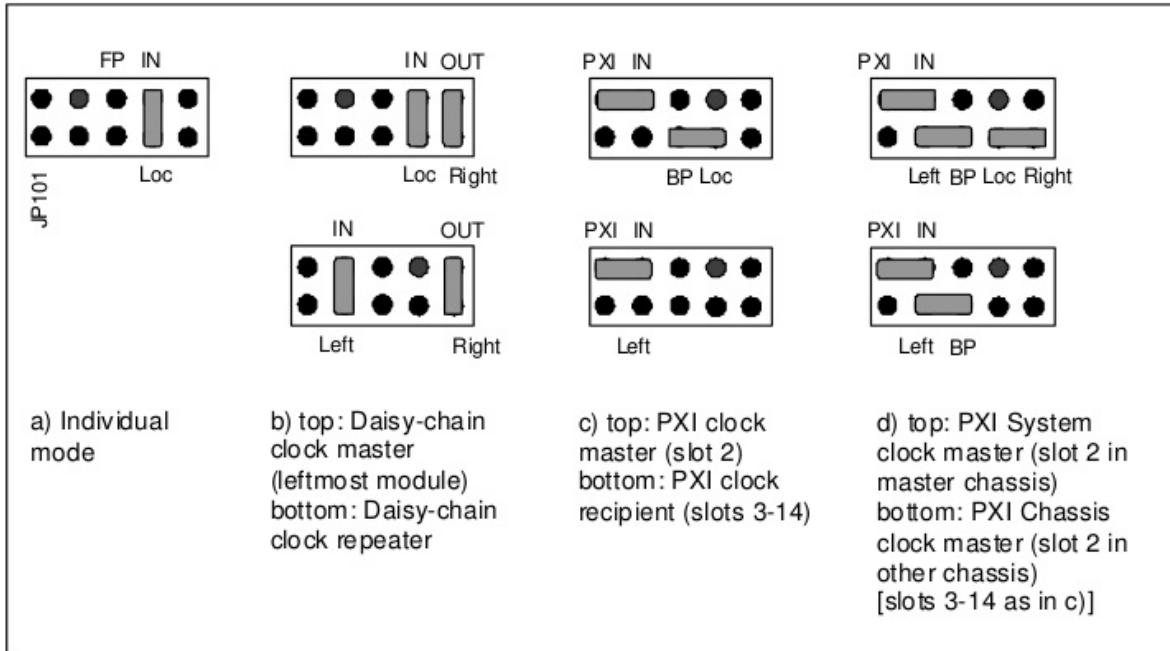
Figure: ClockDistribution

Jumper Settings for different clock distribution modes. An individual module uses its own local clock (a). In a group of modules, there will be one daisy-chained clock master in the leftmost position and several repeaters (b). Alternatively, the PXI clock distribution path can be used, with the module in slot 2 the PXI clock master and the other modules as PXI clock recipients (c). In multi-chassis systems, the module in slot 2 in the clock master chassis should be configured shown in the top picture of (d), the modules in slot 2 in all other chassis as shown in the bottom picture of (d), and modules in any other slot as in the bottom picture of (c). The local clock can be substituted by an LVDS clock input on the front panel by using the “Loc” pin instead of the “FP” pin and setting jumper JP5 to “Clk”.

Individual Clock mode

If only one Pixie-16 module is used in the system, or if clocks between modules do not have to be synchronized, the module should be set into individual clock mode, as shown in Figure a.

Connect pin 7 of JP101 (the clock input) with a shunt to pin 8 (LOC – IN). This will use the local clock crystal as the clock source.

Daisy-chained Clock Mode

The preferred way to distribute clocks between modules is to daisy-chain the clocks from module to module, where each module repeats and amplifies the signal. This requires one master module, located in the leftmost slot of the group of Pixie-16 modules. The master module has the same jumper settings shown in Figure b (top), using its local crystal as the input and sending its output to the right (LOC – IN, OUT – Right). Configure the other modules in the chassis as clock repeaters by setting the jumpers as shown in Figure b (bottom), using the signal from the left neighbor as the input and sending its output to the right (Left – IN, OUT – Right). There must be no gaps between modules.

PXI Clock Mode

A further option for clock distribution is to use the PXI clock distributed on the backplane. This clock is by default generated on the backplane (10MHz), repeated by a fan out buffer and connected to each slot by a dedicated line with minimum skew. Though the 10MHz is too slow to be a useful clock for the Pixie-16, it can be overridden by a signal from a module in slot 2.

The Pixie-16 can be configured to be the PXI clock master in slot 2, by connecting pins 6 and 8 (Loc – BP). All modules, including the clock master, should be set to receive the PXI clock by connecting pin 1 and 3 on JP101 (PXI – IN), see Figure c.

Multi-Chassis Clock Mode

In multi-chassis systems (Figure d), the “left” and “right” connections are taken over by a P16Trigger board on the rear backplane. The module in slot 2 in the clock master chassis should be configured to send its local clock to the P16Trigger board (Loc - Right), receive it back from the P16trigger board and send it to the PXI clock input (Left - BP), and receive the PXI clock as the input for on-board distribution (PXI – IN), see Figure d (top). The modules in slot 2 in all other chassis receive the clock from the P16trigger board and send it to the PXI clock input (Left – BP), and also receive the PXI clock as the input for on-board distribution (PXI – IN), see Figure d (bottom). The modules in any other slot should be set to use the PXI clock as the input (PXI- IN) as in Figure c (bottom).

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Event header as the first 4 words RevD(12-bit,100MHz),RevF(14-bit,100MHz)

Index	Data					
0	[31]	[30:17]	[16:12]	[11:8]	[7:4]	[3:0]
	Finish Code	Event Length	Header Length	CratID	SlotID	Chan#
1	[31:0] EVTTIME_LO[31:0]					
2	[31]	[30:16] CFD Fractional Time[14:0] x 32768			[15:0] EVTTIME_HI[15:0]	
3	[31]	[30:16] Trace Length		[15:0] Event Energy		

(EVTTIME_LO[31:0]+EVTTIME_HI[15:0]x2³²+(CFD_Fractional_Time[14:0]/32768))x10ns

Finish Code: 0-good event,1- pileup event

CFD forced trigger bit: 0- valid ,1-invalid (Threshold was set too high)

Trace Out-of-Range Flag: 0- trace in range, 1- trace is out of range

Event header as the first 4 words RevF(12/14/16-bit,250MHz)

Index	Data					
0	[31]	[30:17]	[16:12]	[11:8]	[7:4]	[3:0]
	Finish Code	Event Length	Header Length	CratID	SlotID	Chan#
1	[31:0] EVTTIME_LO[31:0]					
2	[31]	[30]	[29:16] CFDFractionalTime[13:0]x16384		[15:0] EVTTIME_HI[15:0]	
3	[31]	[30:16] Trace Length		[15:0] Event Energy		

((EVTTIME_LO[31:0]+EVTTIME_HI[15:0]x2³²)x2 - CFD trigger source bit + (CFD_Fractional_Time[13:0]/16384))x4ns

Finish Code: 0-good event,1- pileup event

CFD forced trigger bit: 0- valid ,1-invalid (Threshold was set too high)

CFD trigger source bit:

Trace Out-of-Range Flag: 0- trace in range, 1- trace is out of range

Event header as the first 4 words RevF(12/14-bit,500MHz)

Index	Data					
0	[31]	[30:17]	[16:12]	[11:8]	[7:4]	[3:0]
	Finish Code	Event Length	Header Length	CratID	SlotID	Chan#
1	[31:0] EVTTIME_LO[31:0]					
2	[31:29]	[28:16] CFD Fractional Time[12:0] x 8192		[15:0] EVTTIME_HI[15:0]		
3	[31]	[30:16] Trace Length		[15:0] Event Energy		

(EVTTIME_LO[31:0]+EVTTIME_HI[15:0]x2³²)x10 + ((CFD_Fractional_Time[12:0]/8192)+ CFD trigger source bits[2:0]-1)x2ns

Finish Code: 0-good event,1- pileup event

CFD trigger source bits:

Trace Out-of-Range Flag: 0- trace in range, 1- trace is out of range

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Recording of **raw energy sums and baseline** ENABLED

Index	Data[31:0]
4	EXTTIME_LO
5	EXTTIME_HI[15:0]

Recording of **external timestamp** ENABLED

Index	Data[31:0]
4	QDC sum #0
5	QDC sum #1
6	QDC sum #2
7	QDC sum #3

	QDC sums
8	QDC sum #4
9	QDC sum #5
10	QDC sum #6
11	QDC sum #7

Recording of **QDC sums** ENABLED

Index	Data[31:0]
4	Trailing energy sum
5	Leading energy sum
6	Gap energy sum
7	Baseline value
8	EXTTIME_LO
9	EXTTIME_HI[15:0]

Recording of **raw energy sums and baseline** and **external timestamp** ENABLED

Index	Data[31:0]
4	QDC sum #0
5	QDC sum #1
6	QDC sum #2
7	QDC sum #3
8	QDC sum #4
9	QDC sum #5
10	QDC sum #6
11	QDC sum #7
12	EXTTIME_LO
13	EXTTIME_HI[15:0]

Recording of **QDC sums** and
external timestamp ENABLED

Index	Data [31:0]
4	Trailing energy sum
5	Leading energy sum
6	Gap energy sum
7	Baseline value
8	QDC sum #0
9	QDC sum #1
10	QDC sum #2
11	QDC sum #3
12	QDC sum #4
13	QDC sum #5
14	QDC sum #6
15	QDC sum #7

Recording of raw energy sums and baseline and QDC sums ENABLED

Index	Data [31:0]
4	Trailing energy sum
5	Leading energy sum
6	Gap energy sum
7	Baseline value
8	QDC sum #0
9	QDC sum #1
10	QDC sum #2
11	QDC sum #3
12	QDC sum #4
13	QDC sum #5
14	QDC sum #6
15	QDC sum #7
16	EXTTIME_LO
17	EXTTIME_HI[15:0]

Recording **of raw energy sums and baseline and QDC sums and external timestamp ENABLED**

Recording of trace ENABLED

Index	Data			
n	[31:0]			
	Last word of event header			
n+1	[31:30]	[29:16]	[15:14]	[13:0]
	Not used	ADC Data #1	Not used	ADC Data #0
n+2	[31:30]	[29:16]	[15:14]	[13:0]
	Not used	ADC Data #3	Not used	ADC Data #2
...	...			

If trace recording is enabled, trace data will immediately follow the last word of the event header. Since raw ADC data points are 14-bit number, two 14-bit numbers are packed into one 32-bit word, as shown below. Since the event header could have variable length(4,6,8,10,12,14 ,16 or 18 words) depending on the selection of various output data options, the header length, event length and trace length that are recorded in the first 4 words of the event header should be used to navigate through the output data stream.

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- Good channel. Only channels marked as good will have their events recorded. This setting has no bearing on the channel's capability to issue a trigger. There can be a triggering channel whose data are discarded. Channels not marked as good will be excluded from the automatic offset adjustment.
- Capture trace. When this box is checked, trace will be captured and recorded for each event, along with other list mode information, e.g. timestamp, energy, etc.
- Capture QDC sums. When this box is checked, eight QDC sums will be recorded for each event. QDC sums are consecutive sums of the list mode trace.
- Enable CFD trigger. Check this box to enable this channel's CFD trigger. Otherwise, regular trapezoidal fast trigger will be used.
- Capture raw energy sums and baseline. Check this box to record raw energy sums and baseline value for each event.

The **Base Setup** page controls the analog gain, offset and polarity for each channel. It is useful to click on **Trace & Baseline** in the top control **Monitor** bar to view the signal read from the ADCs while adjusting these parameters. The display shows one or all 16 channels of a module; you can set the sampling interval for each block to capture a longer time frame in **Hist & XDT** page. Click **Draw** to update the graph.

Pulses from the detector should fall in the range from 0 to 12/14/16-bit -1, with the baseline at ~400 to allow for drifts and/or undershoots and no clipping at the upper limit. If there is clipping, adjust the Gain and Offset or click on the Adjust Offsets button to let the software set the DC offsets to proper levels.**TODO** 这段还需要进一步修改

Since the trigger/filter circuits in the FPGA only act on rising pulses, negative pulses are inverted at the input of the FPGA, and the waveforms shown in the ADC trace display include this optional inversion. Thus set the channel's Polarity such that pulses from the detector appear with positive amplitude (rising edge).

In the **Base Setup** tab, you can set the total trace length and the pre-trigger trace delay for the waveforms to be acquired in list mode runs.

The Pixie-16 constantly takes baseline measurements when no pulse is detected and keeps a baseline average to be subtracted from the energy filter output during pulse height reconstruction. Baseline measurements that differ from the average by more than Baseline Cut will be rejected as they are likely contaminated with small pulses below the trigger threshold. You can click on the Acquire Baseline button to view a series of baseline measurements for each channel, and in the single channel view you can build a histogram of baselines to verify that the Baseline Cut does not reject measurements falling into the main (ideally Gaussian) peak in the baseline distribution. Usually, it is sufficient to keep Baseline Cut at its default value.

Note: Since the baseline computation takes into account the exponential decay, no pulses should be noticeable in the baseline display if a) the decay time is set correctly and b) the detector pulses are truly exponential.

Baseline Percent is a parameter used for automatic offset adjustment; by clicking on the Adjust Offsets button, offsets will be set such that the baseline seen in the ADC trace display falls at the Baseline Percent fraction of the full ADC range (e.g. for Baseline Percent = 10% the baseline falls at ADC step 1638 out of 16384 total).

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Assume the digitized waveform stream can be represented by data series $\text{Trace}[i]$, $i=0,1,2,\dots$. First the fast filter response of the digitized waveform is computed as follows:

TODO公式

Where FL is called the fast length and FG is called the fast gap of the digital trapezoidal filter.

- A longer trigger filter rise time averages more samples and thus allows setting lower thresholds without triggering on noise.
- Typically the threshold should be set as low as possible, just above the noise level.
- A longer trigger filter flat top time makes it easier to detect slow rising pulses.

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TODO

The most critical parameter for the energy computation is the signal decay time Tau. It is used to compensate for the falling edge of a previous pulse in the computation of the energy. You can either enter Tau directly for each channel, or enter an approximate value in the right control, select a channel, and click Find it to let the software determine the decay time automatically. Click Accept it to apply the found value to the channel. (If the approximate value is unchanged, the software could not find a better value.) **TODO** 这段还需要进一步修改

- The energy filter flat top time should be about the same as the pulse rise time.
- The energy filter rise time can be varied to balance resolution and throughput. Typically, energy resolution increases with the length of the filter rise time, up to an optimum when longer filters only add more noise into the measurement. The filter dead time is about $TD = 2 \times (T_{rise} + T_{flat})$, and the maximum throughput for Poisson statistics is $1/(TD \cdot e)$. For HPGe detectors, a rise time of 4-6 μ s is usually appropriate.
- Choose the smallest energy filter range that allows setting the optimum energy filter rise time. Larger filter ranges allow

longer filter sums, but increase the granularity of possible values for the energy filter rise time and flat top time and increase the jitter of latching the energy filter output relative to the rising edge of the pulse. This is usually only important for very fast pulses.

Filter range	Filter granularity	Max. SL+SG	Min. SL	Min. SG
1	0.02 us	2.54 us	0.04 us	0.06 us
2	0.04 us	5.08 us	0.08 us	0.12 us
3	0.08 us	10.16 us	0.16 us	0.24 us
4	0.16 us	20.32 us	0.32 us	0.48 us
5	0.32 us	40.64 us	0.64 us	0.96 us
6	0.64 us	81.28 us	1.28 us	1.92 us

Figure: filter range

CFD

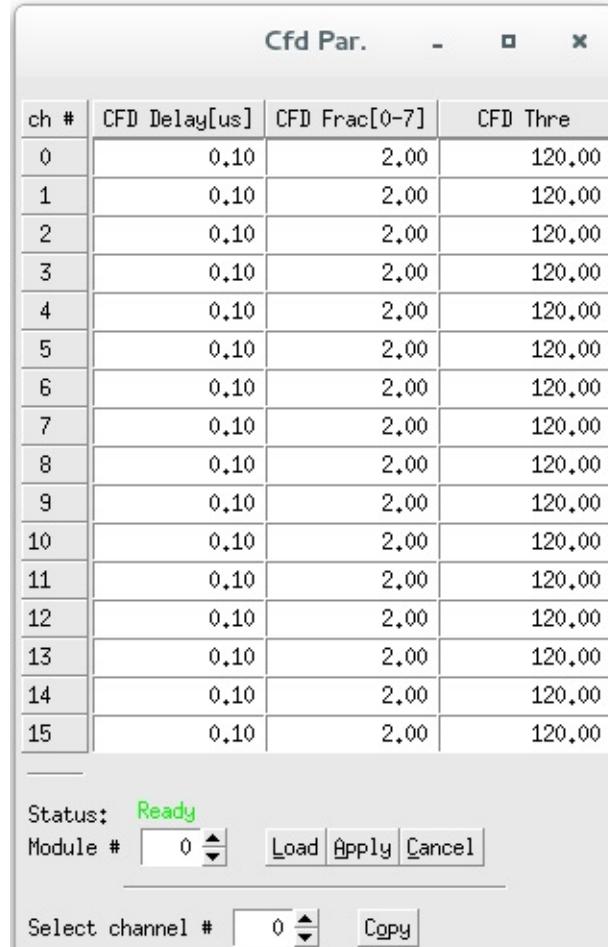
Figure: CFD

.. TODO ..

Then the CFD is computed as follows:

TODO公式

Where D is called the CFD delay length and W is called the CFD scaling factor.

The CFD zero crossing point (ZCP) is then determined when $CFD[i] > 0$ and $CFD[i+1] < 0$. The timestamp is latched at Trace point i, and the fraction time f is given by the ratio of the two CFD response amplitudes right before and after the ZCP.

TODO公式

Where $CFDout1$ is the CFD response amplitude right before the ZCP and $CFDout2$ is the CFD response amplitude (absolute value since $CFDout2$ is negative) right after the ZCP. The Pixie-16 DSP computes the CFD final value as shown below and stored it in the output data stream for online or offline analysis.

Valid CFD Scale values and corresponding CFD scaling factors are

CFD	Scale
0	1
1	0.875
2	0.75
3	0.625
4	0.5
5	0.375
6	0.25
7	0.125

QDC

Qdc Par.								
ch #	QDC len0[us]	QDC len1[us]	QDC len2[us]	QDC len3[us]	QDC len4[us]	QDC len5[us]	QDC len6[us]	QDC len7[us]
0	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
1	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
2	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
3	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
4	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
5	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
6	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
7	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
8	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
9	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
10	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
11	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
12	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
13	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
14	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
15	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13

Status: Ready

Module #

Select channel #

Figure: QDC

.. . TODO. .

Eight QDC sums, each of which can have different length varying from 10 ns to 327.68 μ s, are computed in the signal processing FPGA of a Pixie-16 module for each channel and the sums are written to the output data stream if the user requests so. The recording of QDC sums starts at the waveform point which is Pre-trigger Trace Length earlier than the trigger point, which is either the CFD trigger or local fast trigger depending on whether or not CFD trigger mode is enabled. The eight QDC sums are computed one by one continuously, but they are not overlapping. Both the QDC sum length and the Pre-trigger trace length are set by the user. The recording of QDC sums ends when the eight intervals have passed.

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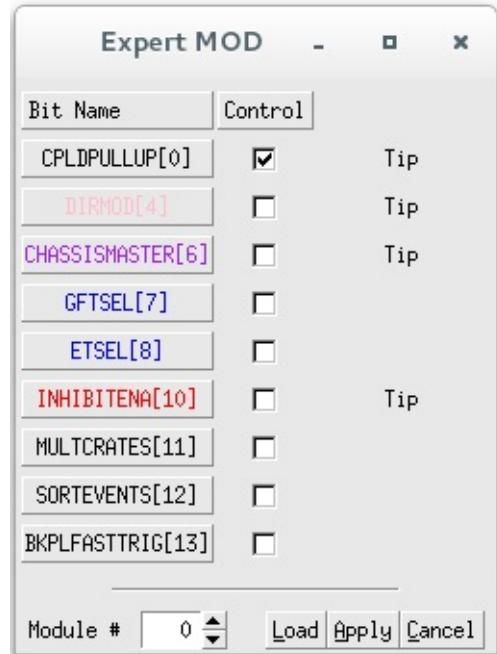
Figure: Module Variables

。。 TODO。 。

Module Control Register B affecting the module as a whole.

- a) Enable pullups for backplane bus lines. This should be enabled for only one module in the crate.
- b) Connect trigger signals to backplane. You can set the module to share triggers over the backplane with other modules, unless you run each module (or channel) independently.
- c) Accept external trigger and run inhibit signals. Enable this option to let this module accept external trigger and run inhibit signals and then put the signals on the backplane so that all modules can see the same signals. This should be enabled for only one module in the crate.
- d) Crate master module (multiple crates only). This option is only used when multiple Pixie-16 crates communicate with each other. By enabling this option, the master module in each crate is responsible for sending synchronization or trigger signals to certain backplane lines. This should be enabled for only one module in the crate.
- e) Enable run inhibit signal input. This option is only applicable to the module which has the “Accept external trigger and run inhibit signals” option enabled. This should be enabled for only one module in the crate.
- f) Multiple crates. This option is only used when multiple Pixie-16 crates communicate with each other.

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。 。 **TODO**。 。

Channel Control Register A affecting each channel individually

- Good channel. Only channels marked as good will have their events recorded. This setting has no bearing on the channel's capability to issue a trigger. There can be a triggering channel whose data are discarded. Channels not marked as good will be excluded from the automatic offset adjustment.

- Histogram energies. When this box is checked, pulse height (energy) computed for each event will be incremented to an energy histogram in the MCA memory.
- Capture trace. When this box is checked, trace will be captured and recorded for each event, along with other list mode information, e.g. timestamp, energy, etc.
- Capture QDC sums. When this box is checked, eight QDC sums will be recorded for each event. QDC sums are consecutive sums of the list mode trace.
- Enable CFD trigger. Check this box to enable this channel's CFD trigger. Otherwise, regular trapezoidal fast trigger will be used.
- Require global external trigger for validation. Check this box to require global external trigger to validate events for this channel.
- Capture raw energy sums and baseline. Check this box to record raw energy sums and baseline value for each event.
- Require channel external trigger for validation. Check this box to require channel external trigger to validate events for this channel.
- Enable pileup rejection. Check this box to enable pileup rejection for this channel. Otherwise, pulses will still be recorded even if they are piled up. **TODO** 这段还需要进一步修改

Logic Set

Logic Trigger

ch #	ExtTrigStr[us]	ExtDeller[us]	ExtTrigoutDel[us]	VetoStr[us]	ChafTrigStr[us]	astTriBalen[us]	Left[0-FFFF]	Itself[0-FFFF]	Right[0-FFFF]	ItselfCoin[0-7]	RightCoin[0-7]	LeftCoin[0-7]	MultiThre[0-31]	Sel Coin/Multi	Sel ChValidTrig
0	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
1	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
2	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
3	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
4	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
5	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
6	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
7	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
8	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
9	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
10	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
11	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
12	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
13	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
14	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0
15	1,50	0,20	0,00	0,30	0,00	0,10	0000	0000	0000	0	0	0	0	0	0

Status: **Ready**

Module # **0**

Select channel # **0**

FastTrigBackplaneLeft[0-FFFF] **0000** FastTrigBackplaneRight[0-FFFF] **0000** Front Panel Monitor Disable Group 000 Ch 00 TestSig 00

Channel Validation Trigger (System FPGA group trigger): (ExtFastTrigGate => 'FT' AND Ext_Fast_Trig_In(module fast trigger)) [Choose 'System FPGA' or 'Front panel channel GATE' in CSRA]

Ch 0-3:	ExtFastTrigCat	GroupTri 0:	GroupTri 0_0	LocalModule Ch	Ch 00	RightModule Ch	Ch 00	LeftModule Ch	Ch 00
Ch 4-7:	ExtFastTrigCat	GroupTri 1:	GroupTri 1_0	LocalModule Ch	Ch 00	RightModule Ch	Ch 00	LeftModule Ch	Ch 00
Ch 8-11:	ExtFastTrigCat	GroupTri 2:	GroupTri 2_0	LocalModule Ch	Ch 00	RightModule Ch	Ch 00	LeftModule Ch	Ch 00
Ch 12-15:	ExtFastTrigCat	GroupTri 3:	GroupTri 3_0	LocalModule Ch	Ch 00	RightModule Ch	Ch 00	LeftModule Ch	Ch 00

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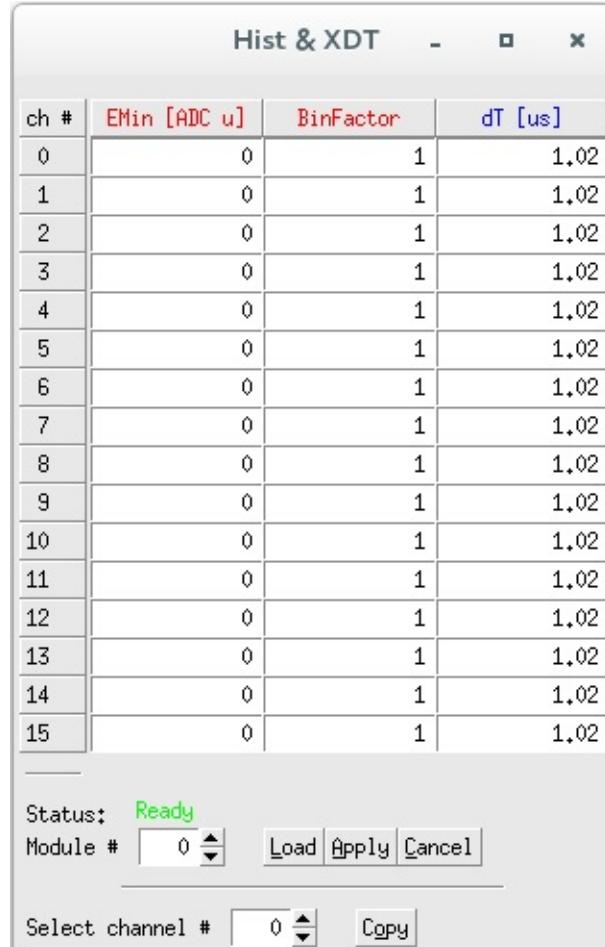
Figure: Hist & XDT

.. TODO ..

The binning factor controls the number of MCA bins in the spectrum. Energies are computed as 16 bit numbers, allowing in principle 64K MCA bins.

However, spectrum memory for each channel is limited to 32K bins, so computed energy values are divided by 2 binning factor before building the histogram. Binning Factor is usually set to 1, but for low count rates and wide peaks, it might be useful to set it to a larger value to obtain a spectrum with fewer bins, but more counts per bin.

Emin is reserved for a future function to subtract a constant “minimum energy” from the computed energy value before binning to essentially cut off the lower end of the spectrum.

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To compute energy filter response offline, the ideal settings are:

- Total trace length > $2(2\text{EFRT} + \text{EFFT})$
- Pre-trigger trace (Trace-delay) length > $(3 * \text{EFRT} + \text{EFFT})$

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- 16 channel gate input
 - 1 module gate input
 - 1 not use current
-

J151-J155

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single-ended TTL external input signals

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© Hongyi Wu *updated: 2018-09-17 13:33:27*

Logic

- Module Fast Trigger (for trigger)
- Module Validation Trigger (for control logic)
- Channel Validation Trigger(for trigger/control logic)
- Veto
- System FPGA (coincidence/multiplicity)

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Channel Validation Trigger(for trigger/control logic)

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System FPGA (coincidence/multiplicity)

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Control logic (module/channel validation trigger)

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© Hongyi Wu updated: 2018-09-17 13:32:55

Multiple Modules Synchronously

- Individual Clock Mode
- PXI Clock Mode
- Daisy-chained Clock Mode
- Multi-Crate Clock Mode
 - Installation of Pixie-16 Modules
 - Clock Jumper (JP101) Settings on the Pixie-16 Modules
 - Cable Connections for Pixie-16 Rear I/O Trigger Modules
 - Jumper Settings on the Pixie-16 Rear I/O Trigger Modules

When many Pixie-16 modules are operated together as a system, it may be required to synchronize clocks and timers between them and to distribute triggers across modules. It will also be necessary to ensure that runs are started and stopped synchronously in all modules. All these signals are distributed through the PXI backplane of the Pixie-16 crate.

Figure: JP101

In a multi-module system there will be one clock master and a number of clock slaves or repeaters. The clock function of a module can be selected by setting shunts on Jumper JP101 near the bottom right corner of the board. The 10-pin Jumper JP101 is shown in the picture on the top with those pins labelled in red color. Shunts are provided to connect pins that are appropriate for each chosen clock distribution mode. Four clock distribution modes, individual clock mode, PXI clock mode, daisy-chained clock mode, and multi- crate clock mode, are described below.

[warning] Please Note In 250 MHz or 500 MHz Pixie-16 modules, the frequency of signal processing clock in the FPGA has been divided down to either 125 MHz or 100 MHz, respectively, for more practical implementation of the design. That division might result in different clock phase and thus different timestamp offset for each channel within a given 250 MHz or 500 MHz Pixie-16 module whenever the module is reinitialized. Calibration might be needed to quantify the different timestamp offset for each channel.

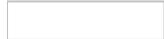
Individual Clock Mode

Figure: individual clock mode

If only one Pixie-16 module is used in the system, or if clocks between modules do not need to be synchronized, the module(s) should be set into individual clock mode. Connect pin 7 of JP101 (the clock input) with a shunt to pin 8 (loc – IN). This will use the 50 MHz local crystal oscillator of the Pixie-16 module as the clock source.

PXI Clock Mode

*Figure: PXI clock mode*

The preferred way to distribute clocks among multiple Pixie-16 modules is to use the PXI clock distributed on the backplane. This clock is by default generated on the backplane and is a 10MHz clock signal, which is then repeated by a fan out buffer and connected to each crate slot by a dedicated line with minimum skew(equal trace length to each slot). Although the 10MHz is too slow to be a useful clock for the Pixie-16, it can be overridden by a local clock signal from a Pixie-16 module that is installed in slot 2 through proper shunt settings on the JP101.

A Pixie-16 module can be configured to be the PXI clock master in slot 2 by connecting pins 6 and 8 (loc – BP) of the JP101. All modules, including the clock master, should be set to receive the PXI clock by connecting pin 1 and 3 on JP101 (PXI – IN). In this way, the 50 MHz clock from the Pixie-16 clock master is distributed to all Pixie-16 modules through the backplane with nearly identical clock phase.

One other advantage of the PXI clock mode over the daisy-chained clock mode, which will be discussed next, is that except for the Pixie-16 master module, which has to be installed in slot 2, other Pixie-16 slave modules can be installed in any other slot of the Pixie-16 crate. In contrast, when the daisy-chained clock mode is used, all Pixie-16 modules have to be installed next to each other, i.e. no gap is allowed between modules.

Daisy-chained Clock Mode

*Figure: daisy-chained clock mode*

A further option for clock distribution is to daisy-chain the clocks from one module to the other, with each module repeating the clock signal and transmitting it to the neighbor on the right. This requires one master module, located in the leftmost slot of the group of Pixie-16 modules. The master module uses its local crystal oscillator as the input and sends its output to the right (loc – IN, out – right). Other Pixie-16 modules in the crate should be configured as clock repeaters by using the signal from the left neighbor as the input and sending its output to the right (left – IN, out – right). However, as mentioned earlier, there must be no slot gap between modules.

Multi-Crate Clock Mode

*Figure: multi-crate clock mode*

In multi- crate systems, a global clock signal can be distributed among these crates using dedicated trigger and clock distribution cards, i.e. the Pixie-16 Rear I/O trigger modules, which are available from XIA.

An example of clock distribution between two crates is illustrated below.

Installation of Pixie-16 Modules

Multiple Pixie-16 modules can be installed in two 14-slot Pixie-16 crates, #1 and #2. For clock distribution purpose, crate #1 is called the Master crate, where the system-wide global clock for all Pixie-16 modules is originated, and crate #2 is called the Slave crate, which receives the global clock from the Master crate.

The Pixie-16 module installed in slot 2 of the Master crate is designated as the System Director Module, whose local 50 MHz crystal oscillator acts as the source of the system-wide global clock. The distribution of the clock signal from the System Director Module to all Pixie-16 modules in the 2-crate system is done through the Pixie-16 Rear I/O trigger modules.

The Pixie-16 module installed in slot 2 of the Slave crate is called the crate Master module, which is responsible for receiving the global clock from the Master crate and sending such clock to all modules in that crate through length-matched traces on the backplane. The System Director Module is also responsible for sending the global clock to all modules in the Master crate. Therefore, it is also a crate Master module. Other modules in these two crates are regular modules. Table shows the different types of modules in a 2-crate system.

Clock Jumper (JP101) Settings on the Pixie-16 Modules

Figure: Module Definitions in a 2-crate System

Cable Connections for Pixie-16 Rear I/O Trigger Modules

Figure: Clock Jumper JP101 Settings in a 2-crate System

Jumper Settings on the Pixie-16 Rear I/O Trigger Modules

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© Hongyi Wu *updated: 2018-05-23 08:18:06*

XIA API

It from **Programmer's Manual Digital Gamma Finder (DGF) PIXIE-16 Version 1.40, October 2009**

```
// Configure modules for communication in PXI chassis
// Use this function to configure the Pixie-16 modules in the PXI chassis.
// NumModules is the total number of Pixie-16 modules installed in the system. PXISlotMap is the pointer to an
array that must have at least as many entries as there are Pixie-16 modules in the chassis.
// PXISlotMap serves as a simple mapping of the logical module number and the physical slot number that the mod
ules reside in. The logical module number runs from 0. For instance, in a system with 5 Pixie-16 modules, these
5 modules may occupy slots 3 through 7. The user must fill PXISlotMap as follows: PXISlotMap = {3, 4, 5, 6, 7
...} since module number 0 resides in slot number 3, etc. To find out in which slot a module is located, any pi
ece of subsequent code can use the expression PXISlotMap[ModNum], where ModNum is the logic module number.
// OfflineMode is used to indicate to the API whether the system is running in OFFLINE mode (1) or ONLINE mode
(0). OFFLINE mode is useful for situations where no Pixie-16 modules are present but users can still test their
calls to the API functions in their application software.
// This function must be called as the first step in the boot process. It makes the modules known to the system
and "opens" each module for communication.
// The function relies on an initialization file (pxisys.ini) that contains information about the Host PC's PCI
buses, including the slot enumeration scheme. XIA's software distribution normally puts this file under the sa
me folder as Pixie-16 software installation folder. However, the user has the flexibility of putting it in othe
r folders by simply changing the definition of the string PCISysIniFile_Windows or PCISysIniFile_Linux in the h
eader part of the file pixie16sys.c, depending on which operating system is being used.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16InitSystem (
    unsigned short NumModules,      // total number of Pixie16 modules in the system
    unsigned short *PXISlotMap,     // an array containing the PXI slot number for each pixie16 module
    unsigned short OfflineMode ); // specify if the system is in offline mode
```

```
// Release user virtual addressees assigned to modules
// Use this function to release the user virtual addressees that are assigned to Pixie-16 modules when these mo
dules are initialized by function Pixie16InitSystem. This function should be called before a user's application
exits.
// If ModNum is set to less than the total number of modules in the system, only the module specified by ModNum
will be closed. But if ModNum is equal to the total number of modules in the system, e.g. there are 5 modules
in the chassis and ModNum = 5, then all modules in the system will be closed altogether. Note that the modules
are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ExitSystem (
    unsigned short ModNum );      // module number
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadModuleInfo (
    unsigned short ModNum,        // module number
    unsigned short *ModRev,       // returned module revision
    unsigned int *ModSerNum,      // returned module serial number
    unsigned short *ModADCBits,   // returned module ADC bits
    unsigned short *ModADCMSPS ); // returned module ADC sampling rate
```

```
// Boot modules so that they can be set up for data taking
// Use this function to boot Pixie-16 modules so that they can be set up for data taking. The function download
s to the Pixie-16 modules the communication FPGA configurations, signal processing FPGA configurations, trigger
FPGA configurations (Revision A modules only), executable code for the digital signal processor (DSP), and DSP
parameters.
// The FPGA configurations consist of a fixed number of words depending on the hardware mounted on the modules;
the DSP codes have a length which depends on the actual compiled code; and the set of DSP parameters always co
nsists of 1280 32-bit words for each module. The host software has to make the names of those boot data files o
n the hard disk available to the boot function.
// If ModNum is set to be less than the total number of modules in the system, only the module specified by Mod
Num will be booted. But if ModNum is equal to the total number of modules in the system, e.g. there are 5 modul
es in the chassis and ModNum = 5, then all modules in the system will be booted.
```

```
// The boot pattern is a bit mask (shown below) indicating which on-board chip will be booted. Under normal circumstances, all on-board chips should be booted, i.e. the boot pattern would be 0x7F. For Rev-B, C, D modules, bit 1, i.e., "Boot trigger FPGA", will be ignored even if that bit is set to 1.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16BootModule (
    char *ComFPGAConfigFile,           // name of communications FPGA configuration file
    char *SPFFPGAConfigFile,          // name of signal processing FPGA configuration file
    char *TrigFPGAConfigFile,         // name of trigger FPGA configuration file
    char *DSPCodeFile,                // name of executable code file for digital signal processor (DSP)
    char *DSPParFile,                 // name of DSP parameter file
    char *DSPVarFile,                 // name of DSP variable names file
    unsigned short ModNum,             // pixie module number
    unsigned short BootPattern ); // boot pattern bit mask
```

```
// Acquire ADC traces in single or multiple modules
// Use this function to acquire ADC traces from Pixie-16 modules. Specify the module using ModNum. If ModNum is set to be less than the total number of modules in the system, only the module specified by ModNum will have its ADC traces acquired. But if ModNum is equal to the total number of modules in the system, then all modules in the system will have their ADC traces acquired.
// After the successful return of this function, the DSP's internal memory will be filled with ADC trace data. A user's application software should then call another function Pixie16ReadSglChanADCTrace to read the ADC trace data out to the host computer, channel by channel.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16AcquireADCTrace (
    unsigned short ModNum );        // module number
```

```
// Read ADC trace data from a channel in a module
// Use this function to read ADC trace data from a Pixie-16 module. Before calling this function, another function Pixie16AcquireADCTrace should be called to fill the DSP internal memory first. Also, the host code should allocate appropriate amount of memory to store the trace data. The ADC trace data length for each channel is 8192. Since the trace data are 16-bit unsigned integers (actually only the lower 14-bit contains real data due to the on-board 14-bit ADC), two consecutive 16-bit words are packed into one 32-bit word in the DSP internal memory. So for each channel, 4096 32-bit words are read out first from the DSP, and then each 32-bit word is unpacked to form two 16-bit words.
// Specify the module using ModNum and the channel on the module using ChanNum. Note that both the modules and channels are counted starting at 0.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglChanADCTrace (
    unsigned short *Trace_Buffer, // trace data
    unsigned int Trace_Length,   // trace length
    unsigned short ModNum,       // module number
    unsigned short ChanNum ); // channel number
```

```
// Transfer data between host and DSP internal memory
// Use this function to directly transfer data between the host and the DSP internal memory of a Pixie-16 module. The DSP internal memory is split into two blocks with address range 0x40000 to 0x4FFFF for the first block and address range 0x50000 to 0x5FFFF for the second block. Within the first block, address range 0x40000 to 0x49FFF is reserved for program memory and shouldn't be accessed directly by the host. Address range 0x4A000 to 0x4A4FF is used by the DSP I/O parameters which are stored in the configuration files (.set files) in the host. Within this range, 0x4A000 to 0x4A33F can be both read and written, but 0x4A340 to 0x4A4FF can only be read but not written. The remaining address range (0x4A500 to 4FFFF) in the first block and the entire second block (0x50000 to 0x5FFFF) should only be read but not written by the host. Use Direction = 1 for read and Direction = 0 for write.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16IMbufferIO (
    unsigned int *Buffer,           // buffer data
    unsigned int NumWords,         // number of buffer words to read or write
    unsigned int Address,          // buffer address
    unsigned short Direction,      // I/O direction
    unsigned short ModNum );       // module number
```

```
// Transfer data between host and DSP external memory
// Use this function to directly read data from or write data to the on-board external memory of a Pixie-16 module. The valid memory address is from 0x0 to 0x7FFFF (32-bit wide). Use Direction = 1 for read and Direction = 0 for write.
// The external memory is used to store the histogram data accumulated for each of the 16 channels of a Pixie-16 module. Each channel has a fixed histogram length of 32768 words(32-bit wide), and the placement of the histo
```

```

gram data in the memory is in the same order of the channel number, i.e. channel 0 occupies memory address 0x0
to 0xFFFF, channel 1 occupies 0x8000 to 0xFFFF, and so on.
// NOTE: another function Pixie16ReadHistogramFromModule can also be used to read out the histograms except tha
t it needs to be called channel by channel.
// In Rev-A modules, part of the external memory is also used to store the list mode data in ping-pong bufferin
g mode. This function can be used to read list mode data from the buffers.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16EMbufferIO (
    unsigned int *Buffer,           // buffer data
    unsigned int NumWords,          // number of buffer words to read or write
    unsigned int Address,           // buffer address
    unsigned short Direction,       // I/O direction
    unsigned short ModNum );        // module number

// Start a list mode data acquisition run
// Use this function to start a list mode data acquisition run in Pixie-16 modules. List mode run is used to co
llect data on an event-by-event basis, gathering energies, timestamps, pulse shape analysis values, and wavefor
ms, for each event. Runs will continue until a preset number of events are reached or the user terminates the r
un by calling function Pixie16EndRun. Once the run is progress, if the run is set to terminate after a given nu
mber of events have been accumulated, another function, Pixie16CheckRunStatus, should be called to check if the
run has finished. To start the data acquisition this function has to be called for every Pixie-16 module in th
e system. If all modules are to run synchronously, The last module addressed will release all others and the ac
quisition starts then. The first module to end the run will immediately stop the run in all other modules.
// Use mode=NEW_RUN (=1) to erase histograms and statistics information before launching the new run. Note that
this will cause a start up delay of up to 1 millisecond. Use mode=RESUME_RUN (=0) to resume an earlier run. Th
is mode has a start up delay of only a few microseconds.
// For Rev-A modules, currently there are 4 list mode run types supported. They are 0x100 (general purpose run)
, 0x101 (without waveforms), 0x102 (without auxiliary data) and 0x103 (energy and timestamp only).
// For Rev-B, C, D modules, there are only one list mode run type supported, that is, 0x100. However, different
output data options can be chosen by enabling or disabling different CHANCSRA bits.
// Histograms and statistics data are updated incrementally from run to run provided RESUME_RUN mode is used.
// ModNum is the module number which starts counting at 0. If ModNum is set to be less than the total number of
modules in the system, only the module specified by ModNum will have its list mode run started. But if ModNum
is set to equal to the total number of modules in the system, then all modules in the system will have their ru
ns started together.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16StartListModeRun (
    unsigned short ModNum,          // module number
    unsigned short RunType,         // run type
    unsigned short mode );          // run mode

// Start a MCA histogram mode data acquisition run
// Use this function to begin a data acquisition run that accumulates energy histograms, one for each channel.
It launches a data acquisition run in which only energy information is preserved and histogrammed locally to ea
ch channel.
// Call this function for each Pixie-16 module in the system. The last module addressed will allow the actual d
ata acquisition to begin. Histogram run can be self-terminating when the elapsed run time exceeds the preset ru
n time, or the user can prematurely terminate the run by calling Pixie16EndRun. On completion, final histogram
and statistics data will be available.
// Use mode=NEW_RUN (=1) to erase histograms and statistics information before launching the new run. Use mode=
RESUME_RUN (=0) to resume an earlier run.
// ModNum is the module number which starts counting at 0. If ModNum is set to be less than the total number of
modules in the system, only the module specified by ModNum will have its histogram run started. But if ModNum
is set to be equal to the total number of modules in the system, then all modules in the system will have their
runs started together.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16StartHistogramRun (
    unsigned short ModNum,          // module number
    unsigned short mode );          // run mode

// Check status of a data acquisition run
// Use this function to check the run status of a Pixie-16 module while a list mode data acquisition run is in
progress. If the run is still in progress continue polling.
// If the return code of this function indicates the run has finished, there might still be some data in the ex
ternal memory (Rev-A modules) or external FIFO (Rev-B, C, D modules) that need to be read out to the host. In a
ddition, final run statistics and histogram data are available for reading out too.
// In MCA histogram run mode, this function can also be called to check if the run is still in progress even th

```

```

ough it is normally self-terminating.
// ModNum is the module number which starts counting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16CheckRunStatus (
    unsigned short ModNum );      // Pixie module number

// Stop a data acquisition run
// Use this function to end a histogram run, or to force the end of a list mode run. In a multi-module system,
if all modules are running synchronously, only one module needs to be addressed this way. It will immediately s
top the run in all other module in the system.
// ModNum is the module number which starts counting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16EndRun (
    unsigned short ModNum );      // Pixie module number

// Compute input count rate
// Use this function to calculate the input count rate on one channel of a Pixie-16 module. This function does
not communicate with Pixie-16 modules. Before calling this function, another function, Pixie16ReadStatisticsFr
mModule, should be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statist
ics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. Mo
dNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeInputCountRate (
    unsigned int *Statistics,
    unsigned short ModNum,
    unsigned short ChanNum );

// Compute output count rate of a channel
// Use this function to calculate the output count rate on one channel of a Pixie-16 module. This function does
not communicate with Pixie-16 modules. Before calling this function, another function, Pixie16ReadStatisticsFr
omModule, should be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statist
ics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. Mo
dNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeOutputCountRate (
    unsigned int *Statistics,
    unsigned short ModNum,
    unsigned short ChanNum );

// Compute live time that a channel accumulated in a run
// Use this function to calculate the live time that one channel of a Pixie-16 module has spent on data acquisi
tion. This function does not communicate with Pixie-16 modules. Before calling this function, another function,
Pixie16ReadStatisticsFromModule, should be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statist
ics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. Mo
dNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeLiveTime (
    unsigned int *Statistics,
    unsigned short ModNum,
    unsigned short ChanNum );

// Compute number of events processed by a channel
// Use this function to calculate the number of events that have been processed by a Pixie-16 module during a d
ata acquisition run. This function is only used by Rev-A modules. This function does not communicate with Pixie
-16 modules. Before calling this function, another function, Pixie16ReadStatisticsFromModule, should be called
to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statist
ics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. Mo
dNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeProcessedEvents (
    unsigned int *Statistics,
    unsigned short ModNum );

```

```

// Compute real time that a channel accumulated in a run
// Use this function to calculate the real time that a Pixie-16 module has spent on data acquisition. This function does not communicate with Pixie-16 modules. Before calling this function, another function, Pixie16ReadStatisticsFromModule, should be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statistics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. ModNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeRealTime (
    unsigned int    *Statistics,
    unsigned short  ModNum );

```

```

PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16complexFFT (
    double *data,
    unsigned int length );

// Test one bit of a 16-bit unsigned integer
PIXIE16APP_EXPORT unsigned short PIXIE16APP_API APP16_TstBit (
    unsigned short bit,
    unsigned short value );

// Set one bit of a 16-bit unsigned integer
PIXIE16APP_EXPORT unsigned short PIXIE16APP_API APP16_SetBit (
    unsigned short bit,
    unsigned short value );

// Clear one bit of a 16-bit unsigned integer
PIXIE16APP_EXPORT unsigned short PIXIE16APP_API APP16_ClrBit (
    unsigned short bit,
    unsigned short value );

// Set one bit of a 32-bit unsigned integer
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API APP32_SetBit (
    unsigned short bit,
    unsigned int   value );

// Clear one bit of a 32-bit unsigned integer
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API APP32_ClrBit (
    unsigned short bit,
    unsigned int   value );

```

```

// Test one bit of a 32-bit unsigned integer
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API APP32_TstBit (
    unsigned short bit,
    unsigned int   value );

```

```

// Program on-board DACs
// Use this function to reprogram the on-board digital to analog converters (DAC) of the Pixie-16 modules. In this operation the DSP uses data from the DSP parameters that were previously downloaded.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SetDACs (
    unsigned short ModNum );

```

```

// Program on-board signal processing FPGAs
// Use this function to program the on-board signal processing FPGAs of the Pixie-16 modules. After the host computer has written the DSP parameters to the DSP memory, the DSP needs to write some of these parameters to the FPGAs. This function makes the DSP perform that action.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ProgramFippi (
    unsigned short ModNum );

```

```

// Adjust DC-offsets in single or multiple modules
// Use this function to adjust the DC-offsets of Pixie-16 modules. Specify the module using ModNum. If ModNum is set to be less than the total number of modules in the system, only the module specified by ModNum will have its DC-offsets adjusted. But if ModNum is set to be equal to the total number of modules in the system, then al

```

```

1 modules in the system will have their DC-offsets adjusted.
// After the DC-offset levels have been adjusted, the baseline level of the digitized input signals will be determined by the DSP parameter BaselinePercent. For instance, if BaselinePercent is set to 10(%), the baseline level of the input signals will be ~ 1638 on the 14-bit ADC scale (minimum: 0; maximum: 16383).
// The main purpose of this function is to ensure the input signals fall within the voltage range of the ADCs to ensure all input signals can be digitized by the ADCs properly.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16AdjustOffsets (
    unsigned short ModNum );

```

```

// Acquire baselines from a module
// Use this function to acquire baselines from Pixie-16 modules. Specify the module using ModNum. If ModNum is set to be less than the total number of modules in the system, only the module specified by ModNum will have its baselines acquired. But if ModNum is set to be equal to the total number of modules in the system, then all modules in the system will have their baselines acquired.
// After the successful return of this function, the DSP's internal memory will be filled with baselines data. Users should then call another function Pixie16ReadSglChanBaselines to read the baselines data out to the host computer, channel by channel.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16AcquireBaselines (
    unsigned short ModNum );           // module number

```

```

// Read baselines from a channel in a module
// Use this function to read baselines data from a Pixie-16 module. Before calling this function, another function Pixie16AcquireBaselines should be called to fill the DSP internal memory first. Also, the host code should allocate appropriate amount of memory to store the baseline data. The baselines data length for each channel is 3640. In the DSP internal memory, each baseline data is a 32-bit IEEE floating point number. After being read out to the host, this function will convert each baseline data to a decimal number. In addition to baseline values, timestamps corresponding to each baseline were also returned after this function call.
// Specify the module using ModNum and the channel on the module using ChanNum. Note that the modules and channels are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglChanBaselines (
    double *Baselines,                // returned baselines values
    double *TimeStamps,               // time stamp for each baseline value
    unsigned short NumBases,          // number of baseline values to read
    unsigned short ModNum,             // module number
    unsigned short ChanNum );         // channel number

```

```

// Ramp Offset DACs of a module and record the baselines
// Use this function to execute the RAMP_OFFSETDACS control task run. Each Offset DAC has 65536 steps, and the RAMP_OFFSETDACS control task ramps the DAC from 0 to 65335 with a step size of 64, i.e., a total of 1024 steps. At each DAC step, the control task computes the baseline value as the representation of the signal baseline and stores it in the DSP memory. After the control task is finished, the stored baseline values are read out to the host computer and saved to a binary file called "rampdacs.bin" in the form of IEEE 32-bit floating point numbers. Users can then plot the baseline values vs. DAC steps to determine the appropriate DAC value to be set in the DSP in order to bring the input signals into the voltage range of the ADCs. However, this function is no longer needed due to the introduction of function Pixie16AdjustOffsets.
// If ModNum is set to less than the total number of modules in the system, only the module specified by ModNum will start the RAMP_OFFSETDACS control task run. But if ModNum is equal to the total number of modules in the system, e.g. there are 5 modules in the chassis and ModNum = 5, then all modules in the system will start the RAMP_OFFSETDACS control task run. Note that the modules are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16RampOffsetDacs (
    double *DCValues,                 // returned DC offset values
    unsigned short NumbCvals,          // number of DC values to read
    unsigned short ModNum );

```

```

// Execute special control tasks
// Use this function to call special control tasks. This may include programming the Fippi or setting the DACs after downloading DSP parameters.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ControlTaskRun (
    unsigned short ModNum,             // Pixie module number
    unsigned short ControlTask,        // Control task number
    unsigned int Max_Poll );          // Timeout control in unit of ms for control task run

```

```

// Find the Baseline Cut values of a module
// Use this function to find the Baseline Cut value for one channel of a Pixie-16 module. The baseline cut value is then downloaded to the DSP, where baselines are captured and averaged over time. The cut value would prevent a bad baseline value from being used in the averaging process, i.e., if a baseline value is outside the baseline cut range, it will not be used for computing the baseline average. Averaging baselines over time improves energy resolution measurement.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16BLcutFinder (
    unsigned short ModNum,           // Pixie module number
    unsigned short ChanNum,          // Pixie channel number
    unsigned int *BLcut );          // BLcut return value

// Find the exponential decay time of a channel
// Use this function to find the exponential decay time constant (Tau value) of the detector or preamplifier signal that is connected to one channel of a Pixie-16 module. The found Tau value is returned via pointer *Tau.

PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16TauFinder (
    unsigned short ModNum,           // Pixie module number
    double *Tau );                 // 16 returned Tau values, in #s

// Write a MODULE level parameter to a module
// Use this function to write a module parameter to a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16WriteSglModPar (
    char *ModParName,              // the name of the module parameter
    unsigned int *ModParData,       // the module parameter value to be written to the module
    unsigned short ModNum );       // module number

// Read a MODULE level parameter from a module
// Use this function to read a module parameter from a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglModPar (
    char *ModParName,              // the name of the module parameter
    unsigned int *ModParData,       // the module parameter value to be read from the module
    unsigned short ModNum );       // module number

// Write a CHANNEL level parameter to a module
// Use this function to write a channel parameter to a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16WriteSglChanPar (
    char *ChanParName,             // the name of the channel parameter
    double ChanParData,            // the channel parameter value to be written to the module
    unsigned short ModNum,         // module number
    unsigned short ChanNum );      // channel number

// Read a CHANNEL level parameter from a module
// Use this function to read a channel parameter from a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglChanPar (
    char *ChanParName,             // the name of the channel parameter
    double *ChanParData,            // the channel parameter value to be read from the module
    unsigned short ModNum,         // module number
    unsigned short ChanNum );      // channel number

// Read histogram data from a module
// Use this function to read out the histogram data from a Pixie-16 module's histogram memory. Before calling this function, the host code should allocate appropriate amount of memory to store the histogram data. The default histogram length is 32768. Histogram data are 32-bit unsigned integers.
// Specify the module using ModNum and the channel on the module using ChanNum. Note that both the modules and channels are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadHistogramFromModule (
    unsigned int *Histogram,        // histogram data
    unsigned int NumWords,          // number of words to be read out
    unsigned short ModNum,          // module number
    unsigned short ChanNum);       // channel number

```

```
// Read run statistics data from a module
// Use this function to read out statistics data from a Pixie-16 module. Before calling this function, the host
// code should allocate appropriate amount of memory to store the statistics data. The number of statistics data
// for each module is fixed at 448. Statistics data are 32-bit unsigned integers.
// Specify the module using ModNum. Note that the modules are counted starting at 0.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadStatisticsFromModule (
    unsigned int *Statistics, // run statistics data
    unsigned short ModNum ); // module number
```

```
// Read histogram data from a module and save to a file
// Use this function to read histogram data from a Pixie-16 module and save the data to a file. New data will b
e appended to the end of the file. So the same file name can be used for multiple modules and the data from eac
h module will be stored in the order that this function is called.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SaveHistogramToFile (
    char *FileName, // histogram data file name
    unsigned short ModNum ); // module number
```

```
// Parse a list mode data file to get events information
// Use this function to parse the list mode events in the list mode data file. The number of events for each mo
dule will be reported.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16GetModuleEvents (
    char *FileName, // the list mode data file name (with complete path)
    unsigned int *ModuleEvents ); // receives number of events for each module
```

```
// Get detailed events information from a data file
// Use this function to retrieve the detailed information of each event in the list mode data file for the desi
gnated module. Before calling this function to get the individual events information, another function Pixie16G
etModuleEvents should be called first to determine the number of events that have been recorded for each module
. If the number of events for a given module is nEvents, a memory block *EventInformation should be allocated w
ith a length of (nEvents*68):
```

```
// EventInformation = (unsigned long *)malloc(sizeof(unsigned long) * nEvents * 68);
// where 68 is the length of the information records of each event (energy, timestamps, etc.) and has the follo
wing structure.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16GetEventsInfo (
    char *FileName, // the list mode data file name (with complete path)
    unsigned int *EventInformation, // to hold event information
    unsigned short ModuleNumber); // the module whose events are to be retrieved
```

```
// Read trace data from a list mode data file
// Use this function to retrieve list mode trace from a list mode data file. It uses the trace length and file
location information obtained from function Pixie16GetEventsInfo for the selected event.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadListModeTrace (
    char *FileName, // list mode data file name
    unsigned short *Trace_Data, // list mode trace data (16-bit words)
    unsigned short NumWords, // number of 16-bit words to be read out
    unsigned int FileLocation); // the location of the trace in the file
```

```
// Read histogram data from a histogram data file
// Use this function to read histogram data from a histogram data file. Before calling this function, the host
// code should allocate appropriate amount of memory to store the histogram data. The default histogram length is
32768. Histogram data are 32-bit unsigned integers.
```

```
// Specify the module using ModNum and the channel on the module using ChanNum. Note that both the modules and
channels are counted starting at 0.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadHistogramFromFile (
    char *FileName, // the histogram data file name (with complete path)
    unsigned int *Histogram, // histogram data
    unsigned int NumWords, // number of words to be read out
    unsigned short ModNum, // module number
    unsigned short ChanNum); // channel number
```

```
// Read DSP parameters from modules and save to a file
// Use this function to save DSP parameters to a settings file. It will first read the values of DSP parameters
on each Pixie-16 module and then write them to the settings file. Each module has exactly 1280 DSP parameter v
alues (32-bit unsigned integers), and depending on the value of PRESET_MAX_MODULES (defined in pixie16app_defs.
h), the settings file should have exactly (1280 * PRESET_MAX_MODULES * 4) bytes when stored on the computer har
d drive.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SaveDSPParametersToFile (
    char *FileName ); // the DSP parameters file name (with complete path)
```

```
// Load DSP parameters to modules from a file
// Use this function to read DSP parameters from a settings file and then download the settings to Pixie-16 mod
ules that are installed in the system. Each module has exactly 1280 DSP parameter values (32-bit unsigned integ
ers), and depending on the value of PRESET_MAX_MODULES (defined in pixie16app_defs.h), the settings file should
have exactly (1280 * PRESET_MAX_MODULES * 4) bytes when stored on the computer hard drive.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16LoadDSPParametersFromFile (
    char *FileName ); // the DSP parameters file name (with complete path)
```

```
// Copy DSP parameters from a module to others
// Use this function to copy DSP parameters from one module to the others that are installed in the system.
// BitMask is bit pattern which designates which items should be copied from the source module to the destinati
on module(s).
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16CopyDSPParameters (
    unsigned short BitMask, // copy items bit mask
    unsigned short SourceModule, // source module
    unsigned short SourceChannel, // source channel
    unsigned short *DestinationMask ); // the destination module and channel bit mask
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadMSGFile (
    char *ReturnMsgStr );
```

```
// Convert a decimal into IEEE 32-bit floating point number
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API Decimal2IEEEFloating(double DecimalNumber);

// Convert an IEEE 32-bit floating point number to a decimal
PIXIE16APP_EXPORT double PIXIE16APP_API IEEEFloating2Decimal(unsigned int IEEEFloatingNumber);
```

```
// Read data from external FIFO and save to a file
// Use this function to read data from the external FIFO of a module. This function can only be used for Pixie-
16 Revision-B, C, and D modules.
// This function first checks the status of the external FIFO of a Pixie-16 module, and if there are data in th
e external FIFO, this function then reads list mode data (32-bit unsigned integers) from the external FIFO. So
this function essentially encapsulates both functions Pixie16CheckExternalFIFOStatus and Pixie16ReadDataFromExt
ernalFIFO within one function. The number of words that are read from the external FIFO is recorded in variable
*nFIFOWords.
// The function also expects setting the value of a variable called "EndOfRunRead" to indicate whether this rea
d is at the end of a run (1) or during the run (0). This is necessary since the external FIFO needs special tre
atment when the host reads the last few words from the external FIFO due to its pipelined structure.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SaveExternalFIFODataToFile (
    char *FileName, // list mode data file name
    unsigned int *nFIFOWords, // number of words read from external FIFO
    unsigned short ModNum, // module number
    unsigned short EndOfRunRead); // indicator whether this is the end of run read
```

```
// Read from or write to registers on a module
// Use this function to read data from or write data to a register in a Pixie-16 module.
// Specify the module using ModNum. Note that the modules are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16RegisterIO (
    unsigned short ModNum, // the Pixie module to communicate to
    unsigned int address, // register address
```

```
    unsigned short direction,      // either MOD_READ or MOD_WRITE
    unsigned int   *value );      // holds or receives the data
```

```
// Read Control & Status Register value from a module
// Use this function to read the host Control & Status Register (CSR) value.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadCSR (
    unsigned short ModNum,
    unsigned int   *CSR );
```

```
// Write to Control & Status Register in a module
// Use this function to write a value to the host Control & Status Register (CSR).
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16WriteCSR (
    unsigned short ModNum,
    unsigned int   CSR );
```

```
// Check status of external FIFO of a module
// Use this function to check the status of the external FIFO of a Pixie-16 module while a list mode data acquisition run is in progress. The function returns the number of words (32-bit) that the external FIFO currently has. If the number of words is greater than a user-set threshold, function Pixie16ReadDataFromExternalFIFO can then be used to read the data from the external FIFO. The threshold can be set by the user to either minimize reading overhead or to read data out of the FIFO as quickly as possible.
// *nFIFOWords returns the number of 32-bit words that the external FIFO currently has.
// ModNum is the module number which starts counting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16CheckExternalFIFOStatus (
    unsigned int   *nFIFOWords,
    unsigned short ModNum );
```

```
// Read data from external FIFO of a module
// Use this function to read data from the external FIFO of a module. This function can only be used for Pixie-16 Revision-B, C, and D modules.
// This function reads list mode data from the external FIFO of a Pixie-16 module. The data are 32-bit unsigned integers. Normally, function Pixie16CheckExternalFIFOStatus is called first to see how many words the external FIFO currently has, then this function is called to read the data from the FIFO.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadDataFromExternalFIFO (
    unsigned int   *ExtFIFO_Data, // To receive the external FIFO data
    unsigned int   nFIFOWords,   // number of words to read from external FIFO
    unsigned short ModNum );    // module number
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ComputeFastFiltersOffline (
    char          *FileName,           // the list mode data file name (with complete path)
    unsigned short ModuleNumber,       // the module whose events are to be analyzed
    unsigned short ChannelNumber,      // the channel whose events are to be analyzed
    unsigned int   FileLocation,        // the location of the trace in the file
    unsigned short RcdTraceLength,     // recorded trace length
    unsigned short *RcdTrace,          // recorded trace
    double         *fastfilter,        // fast filter response
    double         *cfd );             // cfd response
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ComputeSlowFiltersOffline (
    char          *FileName,           // the list mode data file name (with complete path)
    unsigned short ModuleNumber,       // the module whose events are to be analyzed
    unsigned short ChannelNumber,      // the channel whose events are to be analyzed
    unsigned int   FileLocation,        // the location of the trace in the file
    unsigned short RcdTraceLength,     // recorded trace length
    unsigned short *RcdTrace,          // recorded trace
    double         *slowfilter );      // slow filter response
```

```
// Add by Hongyi Wu
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API HongyiWuPixie16ComputeSlowFiltersOffline (
    char          *FileName,           // the list mode data file name (with complete path)
    unsigned short ModuleNumber,      // the module whose events are to be analyzed
    unsigned short ChannelNumber,     // the channel whose events are to be analyzed
    unsigned int   FileLocation,      // the location of the trace in the file
    unsigned short RcdTraceLength,    // recorded trace length
    unsigned short *RcdTrace,         // recorded trace
    double        *slowfilter,        // slow filter response
    unsigned int   bl,
    double        sl,
    double        sg,
    double        tau,
    int           sfr,
    int           pointtobl );
```

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