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---

## Output Data Structures

## Event header as the first 4 words RevD(12-bit,100MHz),RevF(14-bit,100MHz)

Index	Data					
0	[31]	[30:17]	[16:12]	[11:8]	[7:4]	[3:0]
	Finish Code	Event Length	Header Length	CratID	SlotID	Chan#
1	[31:0] EVTTIME_LO[31:0]					
2	[31]	[30:16] CFD Fractional Time[14:0] x 32768			[15:0] EVTTIME_HI[15:0]	
3	[31]	[30:16] Trace Length		[15:0] Event Energy		

(EVTTIME\_LO[31:0]+EVTTIME\_HI[15:0]x2<sup>32</sup>+(CFD\_Fractional\_Time[14:0]/32768))x10ns

Finish Code: 0-good event,1- pileup event

CFD forced trigger bit: 0- valid ,1-invalid (Threshold was set too high)

Trace Out-of-Range Flag: 0- trace in range, 1- trace is out of range

## Event header as the first 4 words RevF(12/14/16-bit,250MHz)

Index	Data					
0	[31]	[30:17]	[16:12]	[11:8]	[7:4]	[3:0]
	Finish Code	Event Length	Header Length	CratID	SlotID	Chan#
1	[31:0] EVTTIME_LO[31:0]					
2	[31]	[30]	[29:16] CFDFractionalTime[13:0]x16384		[15:0] EVTTIME_HI[15:0]	
3	[31]	[30:16] Trace Length		[15:0] Event Energy		

((EVTTIME\_LO[31:0]+EVTTIME\_HI[15:0]x2<sup>32</sup>)x2 - CFD trigger source bit + (CFD\_Fractional\_Time[13:0]/16384))x4ns

Finish Code: 0-good event,1- pileup event

CFD forced trigger bit: 0- valid ,1-invalid (Threshold was set too high)

CFD trigger source bit:

Trace Out-of-Range Flag: 0- trace in range, 1- trace is out of range

## Event header as the first 4 words RevF(12/14-bit,500MHz)

Index	Data					
0	[31]	[30:17]	[16:12]	[11:8]	[7:4]	[3:0]
	Finish Code	Event Length	Header Length	CratID	SlotID	Chan#
1	[31:0] EVTTIME_LO[31:0]					
2	[31:29]	[28:16] CFD Fractional Time[12:0] x 8192		[15:0] EVTTIME_HI[15:0]		
3	[31]	[30:16] Trace Length		[15:0] Event Energy		

(EVTTIME\_LO[31:0]+EVTTIME\_HI[15:0]x2<sup>32</sup>)x10 + ((CFD\_Fractional\_Time[12:0]/8192)+ CFD trigger source bits[2:0]-1 )x2ns  
 Finish Code: 0-good event, 1-pileup event  
 CFD trigger source bits:  
 Trace Out-of-Range Flag: 0-trace in range, 1-trace is out of range

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The screenshot shows a Mozilla Firefox browser window with the title "UEP6000/PL500 - Mozilla Firefox". The address bar shows the URL "UEP6000/PL500" and the IP "222.29.111.26". The page content includes:

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The screenshot shows a Mozilla Firefox browser window titled "UEP6000/PL500 - Mozilla Firefox". The address bar shows the URL "UEP6000/PL500" and the IP address "222.29.111.26". The page content includes:

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## 滤波步长(Filter Range)

To accommodate a wide range of energy filter rise times from tens of nanoseconds to tens of microseconds, the filters are implemented in the FPGA with different clock decimations(filter ranges). The ADC sampling rate is either 2ns, 4ns, or 10ns depending on the ADC variant that is used, but in higher clock decimations, several ADC samples are averaged before entering the energy filtering logic. In filter range 1,  $2^1$  samples are averaged,  $2^2$  samples in filter range 2, and so on. Since the sum of rise time and flat top is limited to 127 decimated clock cycles, filter time granularity and filter time are limited to the values listed in Table .

<b>Filter range</b>	<b>Filter granularity</b>	<b>max. <math>T_{rise}+T_{flat}</math></b>	<b>min. <math>T_{rise}</math></b>	<b>min. <math>T_{flat}</math></b>
1	$0.02\mu s$	$2.54\mu s$	$0.04\mu s$	$0.06\mu s$
2	$0.04\mu s$	$5.08\mu s$	$0.08\mu s$	$0.12\mu s$
3	$0.08\mu s$	$10.16\mu s$	$0.16\mu s$	$0.24\mu s$
4	$0.16\mu s$	$20.32\mu s$	$0.32\mu s$	$0.48\mu s$
5	$0.32\mu s$	$40.64\mu s$	$0.64\mu s$	$0.96\mu s$
6	$0.64\mu s$	$81.28\mu s$	$1.28\mu s$	$1.92\mu s$

Figure: Filter clock decimations and filter time granularity for 100 MHz or 500 MHz

Filter range	Filter granularity	max. $T_{rise}+T_{flat}$	min. $T_{rise}$	min. $T_{flat}$
1	0.016μs	2.032μs	0.032μs	0.048μs
2	0.032μs	4.064μs	0.064μs	0.096μs
3	0.064μs	8.128μs	0.128μs	0.192μs
4	0.128μs	16.256μs	0.256μs	0.384μs
5	0.256μs	32.512μs	0.512μs	0.768μs
6	0.512μs	65.024μs	1.024μs	1.536μs

Figure: Filter clock decimations and filter time granularity for 250 MHz

## 梯形濾波(Trapezoidal Filtering)

From this point onward, only trapezoidal filtering will be considered as it is implemented in a Pixie-16 module according to Equation  $\square$ . The result of applying such a filter with Length  $L=1$  us and Gap  $G=0.4$  us to a gamma-ray event is shown in Figure. The filter output is clearly trapezoidal in shape and has a rise time equal to  $L$ , a flattop equal to  $G$ , and a symmetrical fall time equal to  $L$ . The basewidth, which is a first-order measure of the filter's noise reduction properties, is thus  $2L+G$ .

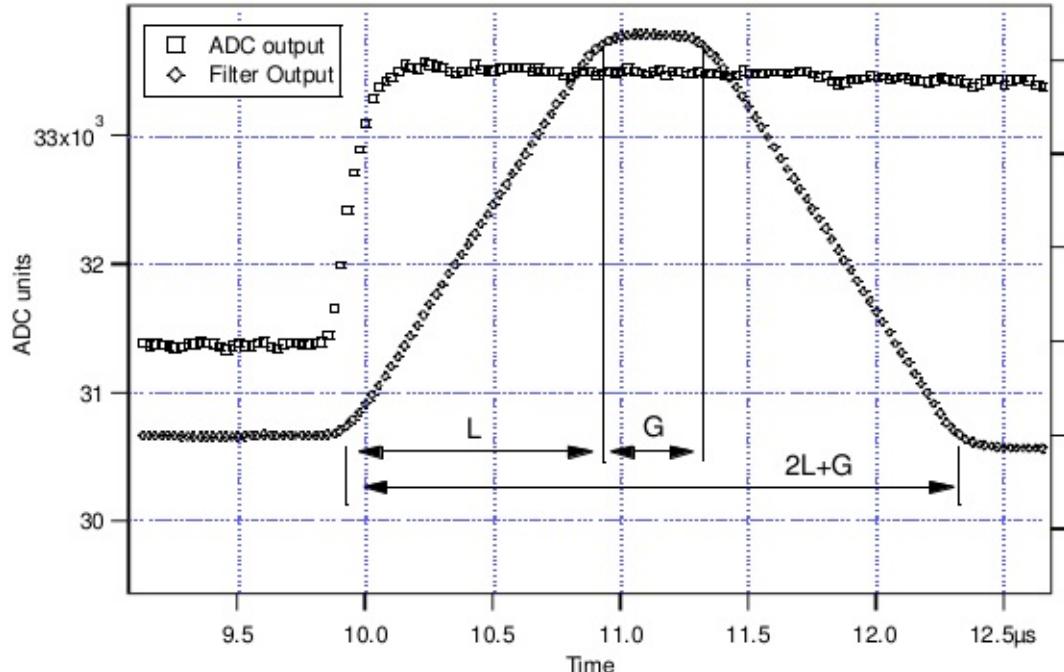
Figure: Trapezoidal filtering of a preamplifier step with  $L=1\mu s$  and  $G=0.4\mu s$

This raises several important points in comparing the noise performance of the Pixie-16 module to analog filtering amplifiers.

- First, semi-Gaussian filters are usually specified by a shaping time.
  - Their rise time is typically twice this and their pulses are not symmetric so that the basewidth is about 5.6 times the shaping time or 2.8 times their rise time.
- Thus a semi-Gaussian filter typically has a slightly better energy resolution than a triangular filter of the same rise time because it has a longer filtering time.

- This is typically accommodated in amplifiers offering both triangular and semi-Gaussian filtering by stretching the triangular rise time a bit, so that the true triangular rise time is typically 1.2 times the selected semi-Gaussian rise time.
- This also leads to an apparent advantage for the analog system when its energy resolution is compared to a digital system with the same nominal rise time.

One important characteristic of a digitally shaped trapezoidal pulse is its extremely sharp termination on completion of the basewidth  $2L+G$ . This may be compared to analog filtered pulses whose tails may persist up to 40% of the rise time, a phenomenon due to the finite bandwidth of the analog filter. As can be seen below, this sharp termination gives the digital filter a definite rate advantage in pileup free throughput.

## 基线与前放衰减时间(Baselines and Preamp. Decay Times)

Figure shows an event over a longer time interval and how the filter treats the preamplifier noise in regions when no gamma-ray pulses are present.

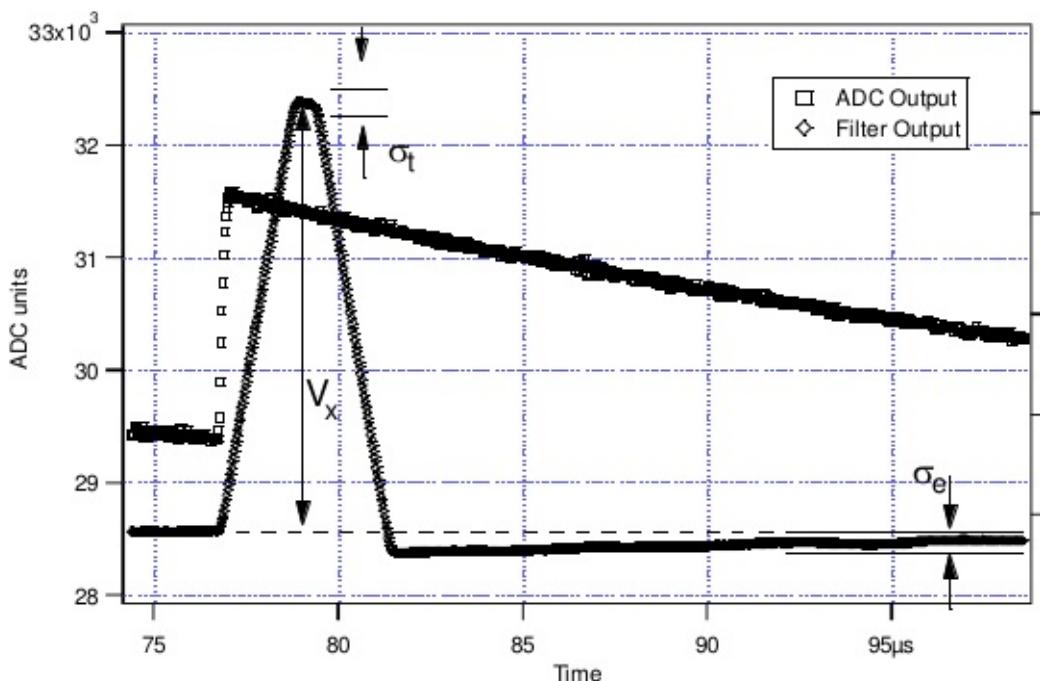
Figure: A gamma-ray event displayed over a longer time period to show baseline noise and the effect of preamplifier decay time

As may be seen the effect of the filter is both to reduce the amplitude of the fluctuations and reduce their high frequency content. This region is called the baseline because it establishes the reference level from which the gamma-ray peak amplitude  $V_x$  is to be measured. The fluctuations in the baseline have a standard deviation  $\sigma_e$  which is referred to as the electronic noise of the system, a number which depends on the rise time of the filter used. Riding on top of this noise, the gamma-ray peaks contribute an additional noise term, the Fano noise, which arises from statistical fluctuations in the amount of charge  $Q_x$  produced when the gamma-ray is absorbed in the detector. This Fano noise  $\square$  adds in quadrature with the electronic noise, so that the total noise  $\square$  in measuring  $V_x$  is found from:

The Fano noise is only a property of the detector material. The electronic noise, on the other hand, may have contributions from both the preamplifier and the amplifier. When the preamplifier and amplifier are both well designed and well matched, however, the amplifier's noise contribution should be essentially negligible. Achieving this in the mixed analog-digital environment of a digital pulse processor is a non-trivial task, however.

With a RC-type preamplifier, the slope of the preamplifier is rarely zero. Every step decays exponentially back to the DC level of the preamplifier. During such a decay, the baselines are obviously not zero. This can be seen in Figure, where the filter output during the exponential decay after the pulse is below the initial level. Note also that the flat top region is sloped downwards.

Using the decay constant  $\square$ , the baselines can be mapped back to the DC level. This allows precise determination of gamma-ray energies, even if the pulse sits on the falling slope of a previous pulse. The value of  $\square$ , being a characteristic of the preamplifier, has to be determined by the user and host software and downloaded to the module.

---

## 堆积检测(Pileup Inspection)


As noted above, the goal is to capture a value of Vx for each gamma-ray detected and use these values to construct a spectrum.

### [info] info

This process is also significantly different between digital and analog systems. In the analog system the peak value must be “captured” into an analog storage device, usually a capacitor, and “held” until it is digitized. Then the digital value is used to update a memory location to build the desired spectrum. During this analog to digital conversion process the system is dead to other events, which can severely reduce system throughput. Even single channel analyzer systems introduce significant deadtime at this stage since they must wait some period (typically a few microseconds) to determine whether or not the window condition is satisfied.

Digital systems are much more efficient in this regard, since the values output by the filter are already digital values. All that is required is to take the filter sums, reconstruct the energy Vx , and add it to the spectrum. In the Pixie-16, the filter sums are continuously updated in the FPGA, and are captured into event buffers. Reconstructing the energy and incrementing the spectrum is done by the DSP, so that the FPGA is ready to take new data immediately (unless the buffers are full). This is a significant source of the enhanced throughput found in digital systems.

The peak detection and sampling in a Pixie-16 module is handled as indicated in Figure *Peak detection and sampling*. Two trapezoidal filters are implemented, a fast filter and a slow filter. The fast filter is used to detect the arrival of gamma-rays, the slow filter is used for the measurement of Vx , with reduced noise at longer filter rise times. The fast filter has a filter length Lf = 0.1us and a gap Gf = 0.1us. The slow filter has Ls = 1.2us and Gs = 0.35us.

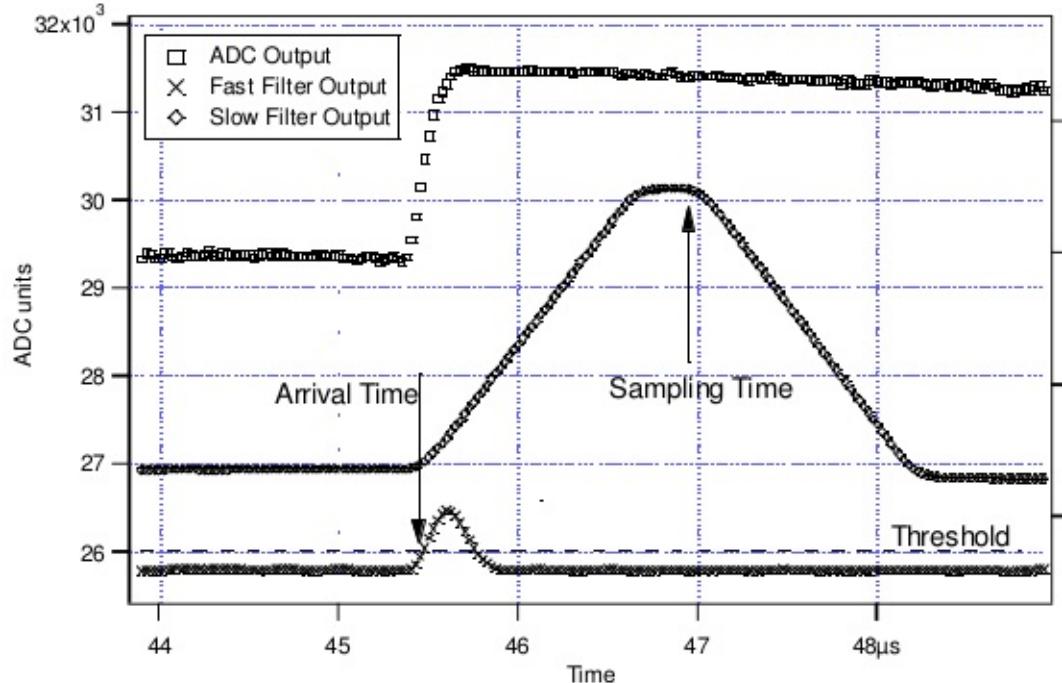*Figure: Peak detection and sampling*

The arrival of the gamma-ray step(in the preamplifier output) is detected by digitally comparing the fast filter output to **THRESHOLD**, a digital constant set by the user. Crossing the threshold starts a delay line to wait **PEAKSAMP** clock cycles to arrive at the appropriate time to sample the value of the slow filter. Because the digital filtering processes are deterministic, **PEAKSAMP** depends only on the values of the fast and slow filter constants. The slow filter value captured following **PEAKSAMP** is then the slow digital filter's estimate of  $V_x$ . Using a delay line allows to stage sampling of multiple pulses even within a **PEAKSAMP** interval (though the filter values themselves are then not correct representations of a single pulse's height).

The value  $V_x$  captured will only be a valid measure of the associated gamma-ray's energy provided that the filtered pulse is sufficiently well separated in time from its preceding and succeeding neighbor pulses so that their peak amplitudes are not distorted by the action of the trapezoidal filter. That is, if the pulse is not piled up. The relevant issues may be understood by reference to Figure, which shows 3 gamma-rays arriving separated by various intervals. The fast filter has a filter length  $L_f = 0.1\mu s$  and a gap  $G_f = 0.1\mu s$ . The slow filter has  $L_s = 1.2\mu s$  and  $G_s = 0.35\mu s$ .

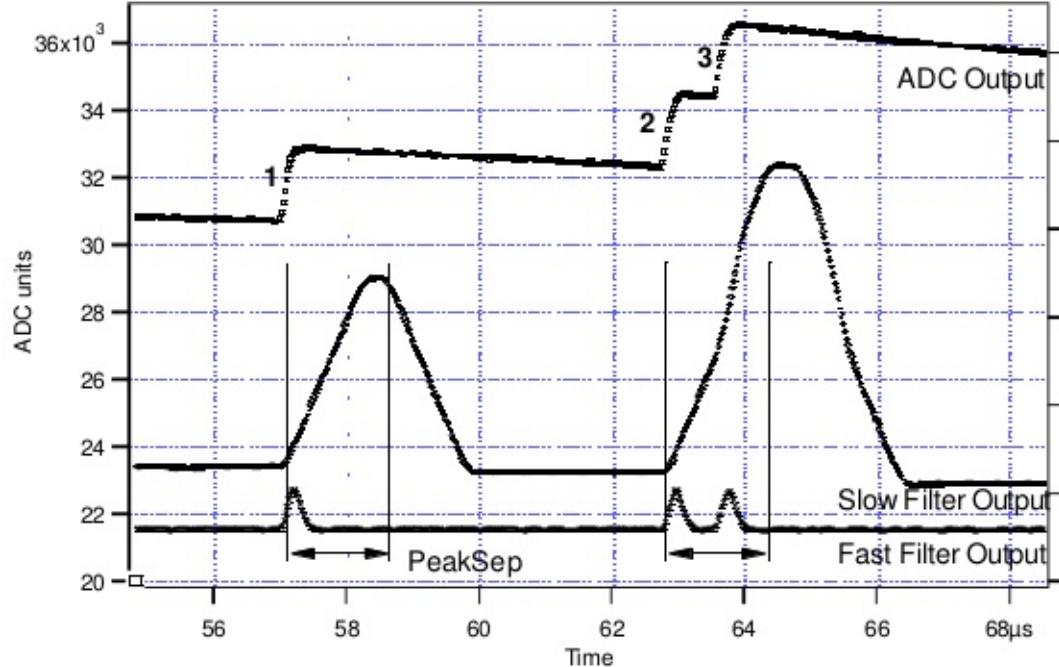
Figure: A sequence of 3 gamma-ray pulses separated by various intervals to show the origin of pileup and demonstrate how it is detected by the Pixie module

Because the trapezoidal filter is a linear filter, its output for a series of pulses is the linear sum of its outputs for the individual members in the series. Pileup occurs when the rising edge of one pulse lies under the peak (specifically the sampling point) of its neighbor. Thus, in Figure , peaks 1 an 2 are sufficiently well separated so that the leading edge of peak 2 falls after the peak of pulse 1. Because the trapezoidal filter function is symmetrical, this also means that pulse 1's trailing edge also does not fall under the peak of pulse 2. For this to be true, the two pulses must be separated by at least an interval of  $L+G$ . Peaks 2 and 3, which are separated by less than 1.0 us, are thus seen to pileup in the present example with a 1.2 us rise time.

This leads to an important point: whether pulses suffer slow pileup depends critically on the rise time of the filter being used. The amount of pileup which occurs at a given average signal rate will increase with longer rise times.

Because the fast filter rise time is only 0.1 us, these gamma-ray pulses do not pileup in the fast filter channel. The Pixie-16 module can therefore test for slow channel pileup by measuring the fast filter for the interval PEAKSEP after a pulse arrival time. If no second pulse occurs in this interval, then there is no trailing edge pileup and the pulse is validated for acquisition. **PEAKSEP** is usually set to a value close to  $L+G+1$ . Pulse 1 passes this test, as shown in Figure. Pulse 2, however, fails the **PEAKSEP** test because pulse 3 follows less than 1.0 us. Notice, by the symmetry of the trapezoidal filter, if pulse 2 is rejected because of pulse 3, then pulse 3 is similarly rejected because of pulse 2.

## CFD

ch #	CFD Delay[us]	CFD Frac[0-7]	CFD Thre
0	0.10	2.00	120.00
1	0.10	2.00	120.00
2	0.10	2.00	120.00
3	0.10	2.00	120.00
4	0.10	2.00	120.00
5	0.10	2.00	120.00
6	0.10	2.00	120.00
7	0.10	2.00	120.00
8	0.10	2.00	120.00
9	0.10	2.00	120.00
10	0.10	2.00	120.00
11	0.10	2.00	120.00
12	0.10	2.00	120.00
13	0.10	2.00	120.00
14	0.10	2.00	120.00
15	0.10	2.00	120.00

Status: **Ready**  
 Module #:

Select channel #:

Figure: CFD

### TODO

## 100 MHz and 250 MHz modules

The following CFD algorithm is implemented in the signal processing FPGA of the 100 MHz(Rev. B, C, D and F) and 250 MHz(Rev. F) Pixie-16 modules.

Assume the digitized waveform can be represented by data series  $\text{Trace}[i]$ ,  $i = 0, 1, 2, \dots$ . First the fast filter response(FF) of the digitized waveform is computed as follows:

Where FL is called the fast length and FG is called the fast gap of the digital trapezoidal filter. Then the CFD is computed as follows:

$$CFD[i + D] = FF[i + D] \times (1 - w/8) - FF[i]$$

Where D is called the CFD delay length and w is called the CFD scaling factor( $w=0, 1, \dots, 7$ ).

The CFD zero crossing point(ZCP) is then determined when  $\boxed{\quad}$  and  $CFD[i + 1] < 0$ . The timestamp is latched at Trace point  $i$ , and the fraction time  $f$  is given by the ratio of the two CFD response amplitudes right before and after the ZCP.

Where CFDout1 is the CFD response amplitude right before the ZCP, and CFDout2 is the CFD response amplitude right after the ZCP(subtraction is used in the denominator since CFDout2 is negative). The Pixie-16 DSP computes the CFD final value as follows and stores it in the output data stream for online or offline analysis.

Where N is scaling factor, which equals to 32768 for 100 MHz modules and 16384 for 250 MHz modules, respectively.

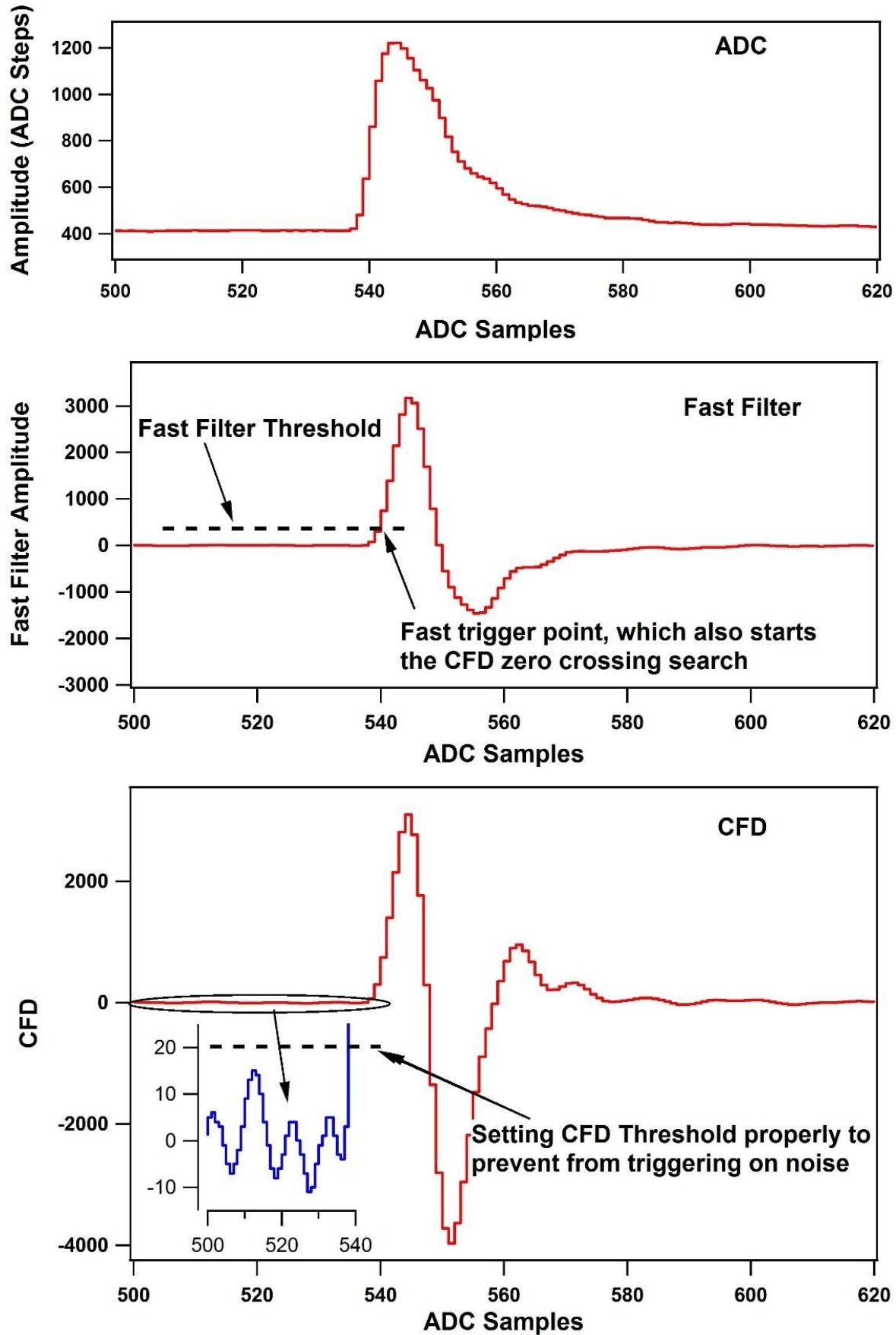
Figure shows a sample ADC trace, its fast filter response and its CFD response, respectively.

The top figure shows a raw ADC trace. After computing the fast filter response on the raw ADC trace using Equation  $FF[i]$ , the fast filter response is compared against the fast filter threshold as shown in the middle figure. The ADC sample where the fast filter response crosses the fast filter threshold is called the fast trigger point, which also starts the search for the CFD zero crossing point.

The CFD response is computed using Equation  $CFD[i + D]$  and is shown in the bottom figure (for actual implementation in the firmware, the fast filter response  $FF$  is delayed slightly before being used for computing the CFD response so that there are sufficient number of CFD response points to look for the zero crossing point after the fast trigger). To prevent premature CFD trigger as a result of the noise in the CFD response before the actual trigger, a DSP parameter called `CFDThresh` is used to suppress those noise-caused zero crossing. However, if a zero crossing point cannot be found within a certain period after the fast trigger (typically 32 clock cycles), e.g., due to unnecessarily high `CFDThresh`, a forced CFD Trigger will be issued and a flag will be set in an event header word to indicate that the recorded CFD time for this event is invalid.

However, the event will still have a valid timestamp which is latched by the fast filter trigger when fast filter crosses over the trigger threshold. The aforementioned CFD parameters correspond to the following DSP parameters.

CFD Parameters	DSP Parameters
FL	FastLength
FG	FastGap
Fast Filter Threshold	FastThresh
D	CFDDelay
W	CFDScale (valid values: 0, 1, 2, ... and 7)
CFD Threshold	CFDThresh

Figure: Corresponding DSP Parameters for the CFD Parameters

#### [info] 250 MHz

In the 250 MHz Pixie-16 modules, the event timestamp is counted with 125 MHz clock ticks, i.e., 8 ns intervals, and two consecutive 250 MHz ADC samples are captured in one 8 ns interval as well.

The CFD trigger also runs at 125 MHz, but the CFD zero crossing point is still reported as a fractional time between two neighboring 250 MHz ADC samples, which are processed by the FPGA in one 125 MHz clock cycle.

However, the CFD zero crossing point could be in either the odd or even clock cycle of the captured 250 MHz ADC waveforms.

Therefore, the firmware outputs a "CFD trigger source" bit in the output data stream to indicate whether the CFD zero crossing point is in the odd or even clock cycle of the captured 250 MHz ADC waveforms.

#### [info] 100 MHz

In the 100 MHz Pixie-16 modules, event timestamp, CFD trigger, and ADC waveform capture are all carried out with the same 100 MHz clock. So there is no need to report "CFD trigger source" for the 100 MHz Pixie-16 modules.

## 500 MHz modules

The CFD algorithm discussed in the previous section for the 100 MHz and 250 MHz Pixie-16 modules can also be written in the following format:

Where  $a(i)$  is the ADC trace data,  $k$  is the index, and  $w$ ,  $B$ ,  $D$ , and  $L$  are CFD parameters.

The CFD algorithm implemented in the 500 MHz Pixie-16 modules is special when compared to the one implemented in the 100 MHz and 250 MHz Pixie-16 modules in terms of the ability to adjust parameters w, B, D, and L.

The reason for this is that in the 500 MHz Pixie-16 modules, ADC data that come into the FPGA at the speed of 500 MHz is first slowed down with a ratio of 1:5, in other words, the FPGA captures 5 ADC samples at the rate of 100 MHz, i.e., every 10 ns. The FPGA then tries to find the CFD trigger point between any two adjacent 2-ns ADC samples within that 10 ns by first building sums of ADC samples and then calculating differences between delayed and non-delayed sums until the zero crossing point is found. However, in the 500 MHz Pixie-16 modules, the FPGA does not have enough resources to build sums for 5 ADC samples in parallel with variable delays. Therefore, the CFD algorithm for the 500 MHz modules was implemented using a set of fixed CFD parameters as shown in Table *Fixed CFD Parameter Values for 500 MHz Pixie-16 Modules*. Tests show these fixed parameters give best performance for LaBr<sub>3</sub>(Ce) detectors.

CFD Parameters	Fixed Values for 500 MHz Modules
w	1
B	5
D	5
L	1

Figure: Fixed CFD Parameter Values for 500 MHz Pixie-16 Modules

The CFD time given by the 500 MHz Pixie-16 modules consists of two parts: a shift within the 5 ADC samples and a fractional time between two ADC samples where the CFD zero crossing occurred. The shift within the 5 ADC samples is reported as the 3-bit CFD trigger source[2:0] is defined as follows.

CFD Trigger Source [2:0]	Zero Crossing Point (ZCP) Location
000	ZCP occurred between the 5th ADC sample of the previous 5-sample group and the 1st ADC sample of the current 5-sample group
001	ZCP occurred between the 1th ADC sample of the current 5-sample group and the 2nd ADC sample of the current 5-sample group
010	ZCP occurred between the 2nd ADC sample of the current 5-sample group and the 3rd ADC sample of the current 5-sample group
011	ZCP occurred between the 3rd ADC sample of the current 5-sample group and the 4th ADC sample of the current 5-sample group
100	ZCP occurred between the 4th ADC sample of the current 5-sample group and the 5th ADC sample of the current 5-sample group
101	Not used
110	Not used
111	CFD trigger is forced, so CFD time is invalid

Figure: Meanings of the CFD Trigger Source for 500 MHz Pixie-16 Modules

The CFD fractional time is given as follows:

# QDC

ch #	QDC len0[us]	QDC len1[us]	QDC len2[us]	QDC len3[us]	QDC len4[us]	QDC len5[us]	QDC len6[us]	QDC len7[us]
0	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
1	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
2	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
3	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
4	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
5	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
6	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
7	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
8	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
9	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
10	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
11	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
12	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
13	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
14	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13
15	0.30	0.63	0.88	1.13	1.38	1.63	1.88	2.13

Status: Ready

Module #

Select channel #

Figure: QDC

Eight QDC sums, each of which can have different lengths, are computed in the Signal Processing FPGA of a Pixie-16 module for each channel and the sums are written to the list mode output data stream if the user requests so.

The recording of QDC sums starts at the waveform point which is *Pre-trigger Trace Length* or *Trace Delay* earlier than the trigger point, which is either the CFD trigger or channel fast trigger depending on whether or not CFD trigger mode is enabled.

The eight QDC sums are computed one by one continuously, but they are not overlapping. The recording of QDC sums ends when the eight intervals have all passed.

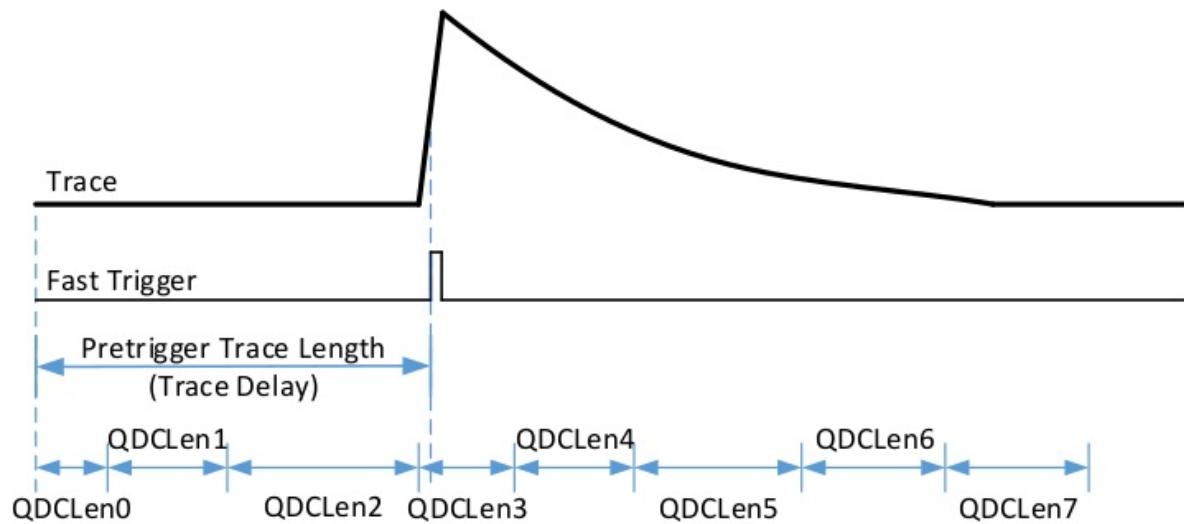
Figure: The 8 QDC sums of a triggered event

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## Module Variables

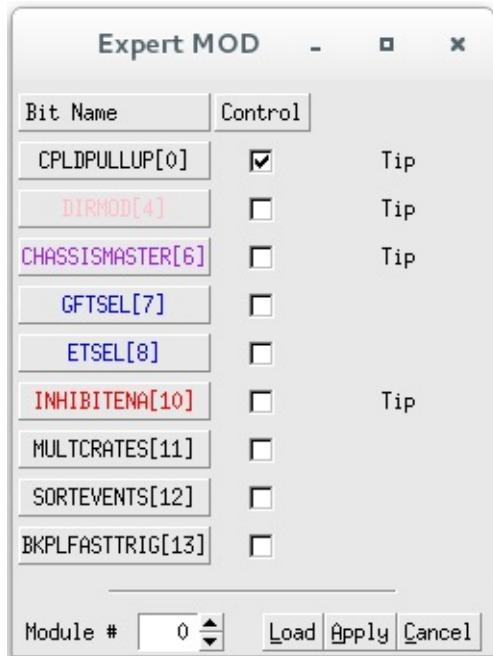
Figure: Module Variables

In addition to distributing the global clock signal, the Pixie-16 rear I/O trigger module can also share global triggers and run synchronization signals. The global trigger signals include the global validation trigger and global fast trigger, plus the Pixie-16 FPGA data storage buffers' full flag signal. The run synchronization signals include synchronous run start and stop signals that can be shared among multiple crates.

In order to enable the distribution of such global triggers and run synchronization signals, certain Pixie-16 parameters have to be set properly. The parameter that controls the trigger distribution and run synchronization is the Module Control Register B (ModCSRB).

ModCSRB is a 32-bit parameter with each of 32 bits controlling different operation modes of the Pixie-16 module.

### [info] Trigger Distribution and Run Synchronization

For the System Director module that is installed in the Master crate, bits 0, 4, 6 and 11 of ModCSRB should be set to 1 (checked & enabled).

For the Crate Master module that is installed in the Slave crate, bits 0, 6 and 11 of ModCSRB should be set to 1 (checked & enabled).

For the General modules that are installed in both the Slave crate and Master crate, bit 11 of ModCSRB should be set to 1 (checked & enabled).

## ## Register definition

Module Control Register B affecting the module as a whole.

- bit 0 - MODCSRB\_CPLDPULLUP

- Enable pullups for PXI trigger lines on the backplane through an onboard CPLD.
- With the pullups, those PXI trigger lines default to logic high state.
- Only when one module actively pulls a line to logic low state will such a line be in the low state.
- Therefore signals transmitted over those PXI trigger lines are actively low signals.
- **Note: enable this bit only for one module per crate (e.g. the crate master module)**
- bit 4 - MODCSRB\_DIRMOD
  - Set this module as the Director module so that it can send triggers, trace and header DPM full signal and run synchronization signal to all crates through the rear I/O trigger modules.
  - Here triggers include fast trigger and validation trigger
  - **Note: enable this bit only for one module among all crates (e.g. the system director module in multi-crate configuration)**
- bit 6 - MODCSRB\_CHASSISMASTER
  - Set this module as the chassis master module so that it can send triggers, trace and header DPM full signal and run synchronization signal to the backplane of the local crate.
  - Here triggers include fast trigger and validation trigger
  - **Note: enable this bit only for one module per crate(e.g. the crate master module)**
- bit 7 - MODCSRB\_GFTSEL
  - Select external fast trigger source(=1: external validation trigger, =0: external fast trigger, in case these two signals are swapped at the Pixie-16 front panel input connectors)
- bit 8 - MODCSRB\_ETSEL
  - Select external validation trigger source(=1: external fast trigger, =0: external validation trigger, in case these two signals are swapped at the Pixie-16 front panel input connectors)
- bit 10 - MODCSRB\_INHIBITENA
  - Enable(=1) or disable(=0) the use of external INHIBIT signal.
  - When enabled, the external INHIBIT signal in the logic high state will prevent the run from starting until this external INHIBIT signal goes to logic low state.
- bit 11 - MODCSRB\_MULTCRATES
  - Set this module to run in the multi-crate mode(=1) or in the local-crate mode(=0).
  - If the module is running in multi-crate mode, it will use the trace and header DPM full signal and run synchronization signal that are generated and distributed among multiple crates.
  - If the module is running in local-crate mode, it will use the trace and header DPM full signal and run synchronization signal generated in the local crate.
- bit 12 - MODCSRB\_SORTEVENTS
  - Sort(=1) or do not sort(=0) events from all 16 channels of a Pixie-16 module based on the timestamps of the events, before storing the events in the external FIFO.
  - Note: all 16 channels must have the same DAQ parameters setting to use this feature
- bit 13 - MODCSRB\_BKPLFASTTRIG
  - Enable(=1) or disable(=0) the sending of 16 local fast triggers to the 16 lines on the backplane of the crate.
  - **Note: only one module can enable this option in each PCI bus segment of a crate(not limited to the crate master module, e.g. any module in each PCI bus segment)**

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## Register definition

Channel Control Register A affecting each channel individually

- bit 0 - CCSRA\_FTRIGSEL
  - Channel fast trigger selection(=1: module fast trigger from the System FPGA; =0: the selection depends on the value of another bit CCSRA\_GROUPTRIGSEL – if CCSRA\_GROUPTRIGSEL = 1, select the channel validation trigger from the System FPGA, and if CCSRA\_GROUPTRIGSEL = 0, select this channel’s local fast trigger)
- bit 1 - CCSRA\_EXTTRIGSEL
  - Module validation trigger selection(=1: module gate input from the Pixie-16 front panel Module Gate LVDS connector; =0: module validation trigger from the System FPGA)
- bit 2 - CCSRA\_GOOD
  - Set this channel as a Good channel(=1) or a not Good channel(=0).
  - **When a channel is set to be a not Good channel, it still generates local fast triggers, which could be used in multiplicity computation, etc., but this channel will not record list mode data or MCA data, and will not update its baseline value**
- bit 3 - CCSRA\_CHANTRIGSEL
  - Channel validation signal selection(=1: channel gate input from the Pixie-16 front panel Channel Gate LVDS connector; =0:channel validation trigger from the System FPGA)
- bit 4 - CCSRA\_SYNCDATAACQ
  - Choose the level of synchronous data acquisition for this channel(=1: stops taking data when the trace or header DPM for any channel of any Pixie-16 module in the system is full; =0: stops taking data only when the trace or header DPM for this channel of this Pixie-16 module is full)
- bit 5 - CCSRA\_POLARITY
  - Choose this channel’s input signal polarity(=1: invert input signal’s polarity; =0: do not invert input signal’s polarity).
  - **Please note in Pixie-16, signal processing requires positive rising input signal. So if input signal has a negative falling edge, it should be inverted by setting this CCSRA\_POLARITY bit to 1**
- bit 6 - CCSRA\_VETOENA
  - Enable(=1) or disable(=0) this channel’s veto.
  - If veto is enabled, this channel’s fast trigger will be vetoed by either the module veto signal(see bit 20 CCSRA\_MODVETOSEL below) or channel veto signal(see bit 19 CCSRA\_CHANVETOSEL below).
  - But if veto is disabled, this channel’s fast trigger will not be vetoed by either veto signal, even if either veto signal is present
- bit 7 - CCSRA\_HISTOE
  - Enable(=1) or disable(=0) the histogramming of pulse energy values in the onboard MCA memory.
  - However, the current Pixie-16 firmware always histograms pulse energy values in the onboard MCA memory.
  - So this CCSRA\_HISTOE is essentially not in use at the moment
- bit 8 - CCSRA\_TRACEENA
  - Enable(=1) or disable(=0) trace capture in the list mode run for this channel
- bit 9 - CCSRA\_QDCENA
  - Enable(=1) or disable(=0) QDC sums recording in the list mode run for this channel.
  - There are a total of 8 QDC sums for each event
- bit 10 - CCSRA\_CFDMODE
  - Enable(=1) or disable(=0) CFD trigger in the list mode run for this channel.
  - CFD trigger is used to latch sub-sample timing for the event time of arrival or timestamp
- bit 11 - CCSRA\_GLOBTRIG
  - Enable(=1) or disable(=0) the requirement of module validation trigger for this channel.
  - If enabled, only when module validation trigger overlaps the channel fast trigger will the events be recorded for this channel
- bit 12 - CCSRA\_ESUMSENA
  - Enable(=1) or disable(=0) the recording of raw energy sums and baseline values in the list mode run for this channel.
  - There are a total of three raw energy sums and one baseline value for each event.
  - **Please note the baseline value is stored in the format of 32-bit IEEE float point(IEEE 754)**

- bit 13 - CCSRA\_CHANTRIG
  - Enable(=1) or disable(=0) the requirement of channel validation trigger for this channel.
  - If enabled, only when channel validation trigger overlaps the channel fast trigger will the events be recorded for this channel
- bit 14 - CCSRA\_ENARELAY
  - Switch between two attenuations or gains for the input signal in this channel through an input relay(=1: close the input relay resulting in no input signal attenuation; =0: open the input relay resulting in a 1/4 input signal attenuation)
- bit 15/16 - CCSRA\_PILEUPCTRL/CCSRA\_INVERSEPILEUP
  - Control normal pileup rejection(bit 15) and inverse pileup rejection(bit 16) for list mode runs:
  - 1) Bits [16:15] = 00, record all events
  - 2) Bits [16:15] = 01, only record single events, i.e., reject piled up events
  - 3) Bits [16:15] = 10, record everything for piled up events, but will not record trace for single events even if trace recording is enabled, i.e., only record event header
  - 4) Bits [16:15] = 11, only record piled up events, i.e., reject single events
  - In all cases, if the event is piled up, no energy will be computed for such event
- bit 17 - CCSRA\_ENAENERGYCUT
  - Enable(=1) or disable(=0) the “no traces for large pulses” feature.
  - If enabled, trace will not be recorded if the event energy is larger than the value set in DSP parameter EnergyLow
- bit 18 - CCSRA\_GROUPTRIGSEL
  - Select channel fast trigger – this bit works together with the CCSRA\_FTRIGSEL bit(bit 0): if CCSRA\_FTRIGSEL=1, this CCSRA\_GROUPTRIGSEL bit has no effect; if CCSRA\_FTRIGSEL=0, then if CCSRA\_GROUPTRIGSEL=1, select the channel validation trigger from the System FPGA, and if CCSRA\_GROUPTRIGSEL=0, select this channel’s local fast trigger
- bit 19 - CCSRA\_CHANVETOSEL
  - Channel veto signal selection(=1: channel validation trigger from the System FPGA; =0: channel gate input from the Pixie-16 front panel Channel Gate LVDS connector)
- bit 20 - CCSRA\_MODVETOSEL
  - Module veto signal selection(=1: module validation trigger from the System FPGA; =0: module gate input from the Pixie-16 front panel Module Gate LVDS connector)
- bit 21 - CCSRA\_EXTTSENA
  - Enable(=1) or disable(=0) the recording of the 48-bit external clock timestamp in the event header during list mode run for this channel

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## Hist & XDT

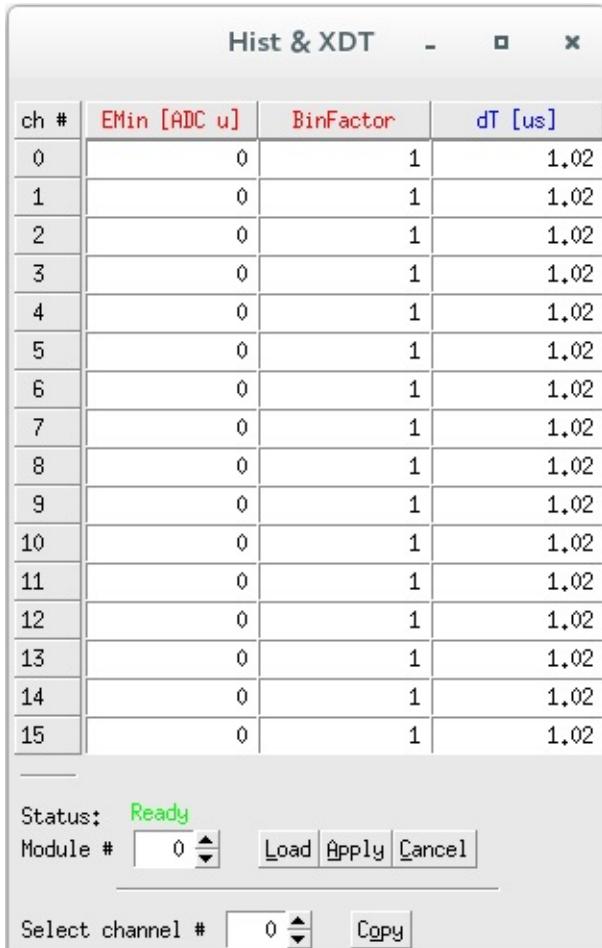
Figure: Hist &amp; XDT

The binning factor controls the number of MCA bins in the spectrum. Energies are computed as 16 bit numbers, allowing in principle 64K MCA bins.

However, spectrum memory for each channel is limited to 32K bins, so computed energy values are divided by \$2^{binning factor}\$.  $E_{min}$  is reserved for a future function to subtract a constant “minimum energy” from the computed energy value before binning to essentially cut off the lower end of the spectrum.

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In the figure, EFRT stands for Energy Filter Rise Time and EFFT stands for Energy Filter Flat Top.

To compute energy filter response offline, the ideal settings are:

- Total trace length >  $2(2 \text{EFRT} + \text{EFFT})$
- Pre-trigger trace (Trace-delay) length >  $(3 \times \text{EFRT} + \text{EFFT})$

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the ADC trace display also includes the option to view a FFT of the acquired trace. This is useful to diagnose noise contributions.

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```

s
created the directory /home/wuhongyi/data/0320.
open: /home/wuhongyi/data/0320/data_R0320_M00.bin
open: /home/wuhongyi/data/0320/data_R0320_M01.bin
open: /home/wuhongyi/data/0320/data_R0320_M02.bin
open: /home/wuhongyi/data/0320/data_R0320_M03.bin
RUN START
SHM Opend!
Running No. 320
Auto Run Mode: CLOSE
Send Shared Memory Online: OPEN
Write Data To Harddisk: OPEN
a
The times for each run: 180 s
h

[ q]   Quit
[ s]   Start/Stop Acquisition
[ o]   Open/Close Send Shared Memory Online
[ w]   Open/Close Write Data To Harddisk
[ a]   Open/Close Auto Run Mode
[ e]   Update Online Energy Spectrum
[ c]   Clean Screen
[ h]   Print Run Status
-----
Running No. 320
Auto Run Mode: OPEN - 180 s per run
Send Shared Memory Online: OPEN
Write Data To Harddisk: OPEN

```

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```

# Only use in NOGUI, unit: second
AutoRunModeTimes 180

```

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## Front Panel

- Analog Signal Input Connectors(all revisions)
- LVDS I/O Port (all revisions)
- Digital I/O Connectors(Rev. F only)
- Front Panel LEDs (all revisions)
- 3.3V I/O Connector(Rev. D only)
- GATE Inputs(Rev. D only)
- 3.3V I/O Connector(Rev. B and C only)
- Digital Signals in Standard Firmware(all revisions)

On the front panel of each Pixie-16 module, there are 16 analog signal input connectors, one LVDS I/O port, five digital I/O connectors as well as three LEDs near the bottom of the front panel. In addition, a sticker showing Pixie-16 model number (e.g., P16L-250-14, meaning the 14-bit, 250 MHz variant of the Pixie-16) is affixed to the top handle of the front panel, and another sticker indicating the serial number of the Pixie-16 module (e.g., S/N 1100) is placed at the bottom handle of the front panel.

## Analog Signal Input Connectors(all revisions)

- Connector Labels
  - 0 to 15 for 16 channels
- Connector Type
  - SMB Jack

Each Pixie-16 module accepts 16 analog input signals, and each input connector is a SMB Jack (male contact) connector.

## LVDS I/O Port (all revisions)

Connected Signals	4 LVDS pairs ( $F_{O1p}/F_{O1n}$ , $F_{I1p}/F_{I1n}$ , $F_{I5p}/F_{I5n}$ , $F_{O5p}/F_{O5n}$ , see below for pin layout)
Signal Direction	Input or Output, software configurable
Cable Type	Cat 5 or Cat 6 (the same ones used for Ethernet)

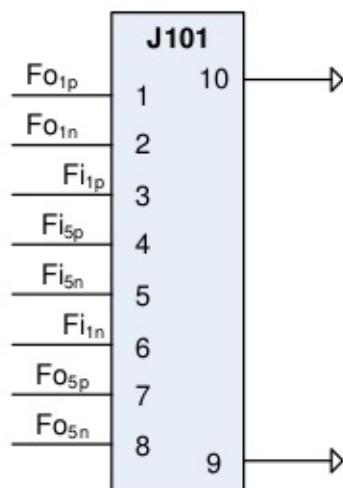
Figure: LVDS I/O Port

Each Pixie-16 module is equipped with one LVDS I/O port on its front panel. LVDS stands for low voltage differential signaling. The LVDS I/O connector is a RJ45 connector, which implies that the same Cat 5 or Cat 6 Ethernet cables can be used to connect signals to or from this I/O port. However, no Ethernet connectivity is available through this Pixie-16 I/O port.

Four differential signal pairs, i.e., between pin-pairs  $Fo_{1p}/Fo_{1n}$ ,  $Fi_{1p}/Fi_{1n}$ ,  $Fi_{5p}/Fi_{5n}$ , and  $Fo_{5p}/Fo_{5n}$ , are available from this I/O port. Each pair can be configured as either an input or output signal.

## Digital I/O Connectors(Rev. F only)

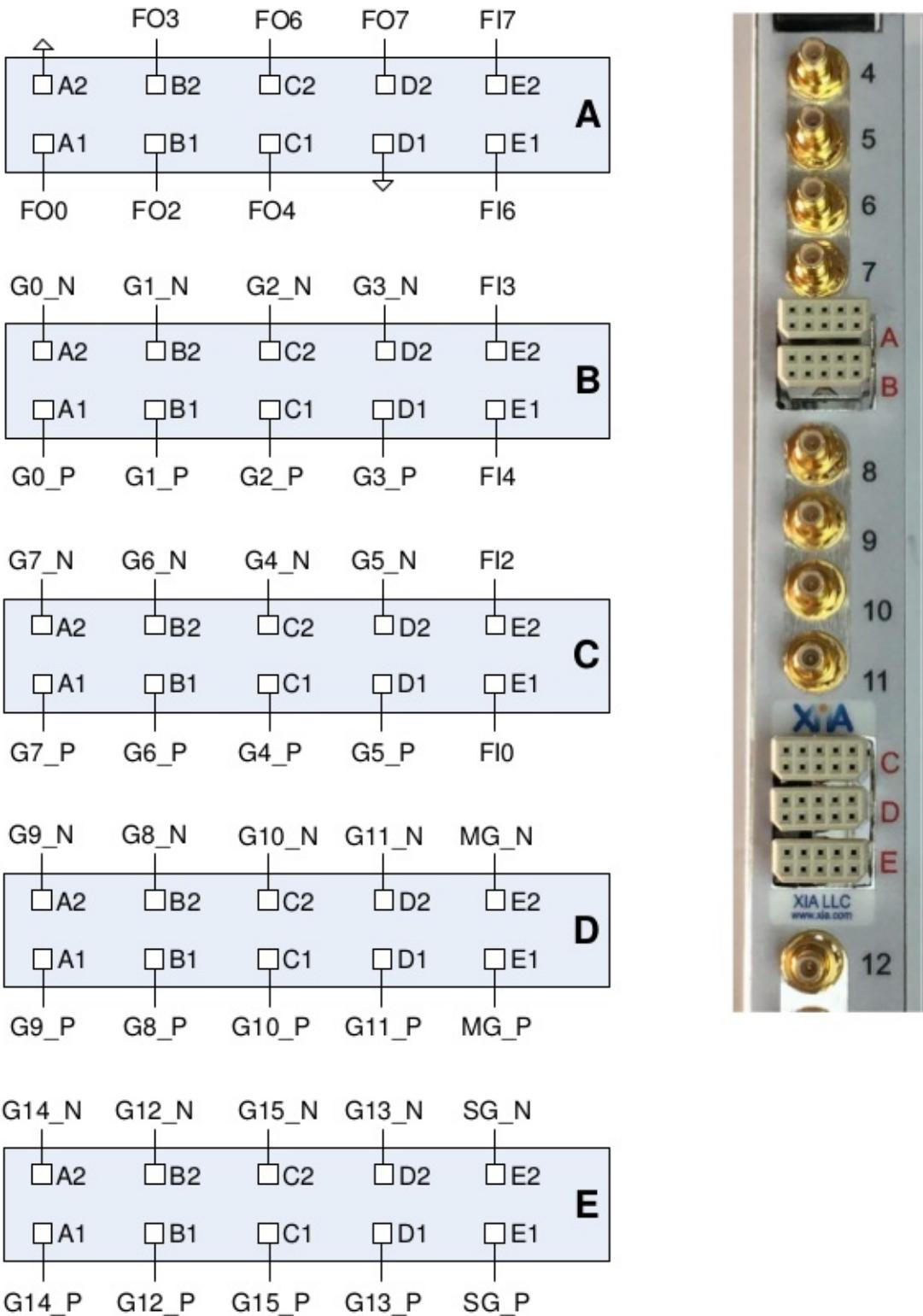
Figure: Digital I/O Connectors A, B, C, D and E

The Pixie-16 Rev. F modules are equipped with five har-link® connectors on its front panel which act as digital I/O connectors. The 2mm pitch har-link® connector from HARTING is designed for high speed data transfer with rates up to 2 Gbit/s. Its EMI shielding, shown in Figure 2mm pitch har-link® connector from HARTING, guarantees excellent performance in EM-polluted environment.

Figure: 2mm pitch har-link® connector from HARTING

Each har-link® connector has 2 rows with 5 pins on each row, and is labelled using one of the five letters in red color font, from A to E. The signals connected to each pin of these five connectors are shown in Table Rev. F Module's Digital I/O Connectors.

Connector Type	har-link® (HARTING, 2mm pin spacing)
FI <sub>0</sub> , FI <sub>2</sub> , FI <sub>3</sub> , FI <sub>4</sub> , FI <sub>6</sub> , FI <sub>7</sub>	TTL digital input signals (max. 5V)
FO <sub>0</sub> , FO <sub>2</sub> , FO <sub>3</sub> , FO <sub>4</sub> , FO <sub>6</sub> , FO <sub>7</sub>	Digital outputs for test/debug purpose (TTL 5V)
Gx_P/Gx_N (x=0-15)	Channel Gate Inputs (0-15 for 16 channels) (LVDS format)
MG_P/MG_N	Module Gate Input (LVDS format)
SG_P/SG_N	Spare Gate Input (LVDS format)

Figure: Rev. F Module's Digital I/O Connectors

Among them, FI<sub>0</sub>, FI<sub>2</sub>, FI<sub>3</sub>, FI<sub>4</sub>, FI<sub>6</sub>, FI<sub>7</sub> are six TTL digital input signals. They can be signals like global fast trigger, global validation trigger, external clock, run inhibit, etc. The specific usage of each input pin is determined by the specific firmware that is downloaded to the Pixie-16 module (see Table 1-9 for input signals supported by the standard firmware). The six digital output signals, FO<sub>0</sub>, FO<sub>2</sub>, FO<sub>3</sub>, FO<sub>4</sub>, FO<sub>6</sub>, FO<sub>7</sub>, which are connected to six test output pins on the System FPGA of the Pixie-16, can be used to assist a user in the process of system setup. These test pins are connected to various internal signals of the Pixie-16 to provide insight of the current status of the system.

The Channel Gate Inputs (0-15 for 16 channels) are LVDS format input signals which independently gate the data acquisition of each of the 16 channels of a Pixie-16 module.

The Channel Gate signal is level sensitive signal, i.e., when the level of the Channel Gate Signal is logic high(1), the gate signal is effective; when the level of the Channel Gate Signal is logic low(0), the gate signal is not in use. In normal cases, the Channel Gate Signal is set up to veto the data acquisition in a given channel, i.e., at the time of the arrival of fast trigger in that channel, if the Channel Gate Signal is logic high(1), that fast trigger is discarded since it is vetoed. However, this type of logic can be

reversed through setting corresponding registers in the FPGA via software. In such cases, the Channel Gate Signal is set up to validate the data acquisition in a given channel, i.e., at the time of fast trigger in that channel, only if the Channel Gate Signal is logic high(1) will that fast trigger be accepted to have the event recorded.

The Module Gate Input is a LVDS format signal that gates the data acquisition in all 16 channels of a Pixie-16 module. It is also a level sensitive signal, i.e., when the level of the Module Gate Signal is logic high(1), the gate signal is effective; when the level of the Module Gate Signal is logic low(0), the gate signal is not in use. In normal cases, the Module Gate Signal is set up to veto the data acquisition in all 16 channels, i.e., at the time of the arrival of fast trigger in any of the 16 channels, if the Module Gate Signal is logic high(1), that fast trigger of that channel is discarded since it is vetoed. However, this type of logic can be reversed through setting corresponding registers in the FPGA via software.

In such cases, the Module Gate Signal is set up to validate the data acquisition in all 16 channels, i.e., at the time of fast trigger in any of the 16 channel, only if the Module Gate Signal is logic high(1) will that fast trigger of that channel be accepted to have the event recorded.

The Spare Gate Input is a LVDS format signal that is reserved for special applications.

Such applications typically require development of custom firmware to support special functionalities of the Pixie-16 system.

Figure: frontpanel

## Front Panel LEDs (all revisions)

Near the bottom of the Pixie-16 front panel, there are three LEDs. They are labelled as RUN, I/O, and ERR, from left to right. They correspond to three different colors, **green, yellow, and red**, respectively.

LED Name	Color	Function
RUN	Green	ON when run is in progress, and OFF if run is stopped or not started yet
I/O	Yellow	Flashing when there is I/O activity on the PCI bus between the Pixie-16 module and host computer
ERR	Red	ON when there is no more space in the External FIFO for storage of list mode event data, and OFF when there is sufficient space to store at least one more list mode event data <b>(ON does not indicate any actual error condition. Rather, it simply indicates the External FIFO's FULL condition)</b>

Figure: Front Panel LEDs for the Pixie-16 Modules

The **RUN LED** will be turned on when a run in the Pixie-16 module is in progress, and will be turned off when the run is stopped or not started yet.

The **I/O LED** will blink when there is I/O activity on the PCI bus between the Pixie-16 module and host computer.

The **ERR LED** is, in fact, not to signal any error condition in the Pixie-16 module. Instead, it is used to indicate whether or not the External FIFO of the Pixie-16 module is full. It will be ON when there is no more space in the External FIFO for storage of list mode event data, and OFF when there is sufficient space to store at least one more list mode event data. When the External FIFO is full, no more list mode event data can be written into it until the host software reads out part of the data in the External FIFO through the PCI bus.

### 3.3V I/O Connector(Rev. D only)

Connector Type	Single-ended, 2mm pin spacing
FO <sub>0</sub> , FO <sub>2</sub> , FO <sub>3</sub> , FO <sub>4</sub> , FO <sub>6</sub> , FO <sub>7</sub>	Digital outputs for test/debug purpose (3.3V)

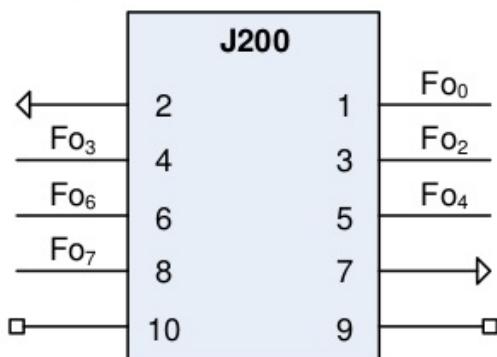
Figure: Rev. D Module's 3.3V I/O Connector

On Rev. D Pixie-16 modules, between analog input SMB connectors for channel 7 and channel 8, respectively, is the 3.3V I/O Connector(J200). It has 10 single-ended pins with 2mm spacing. Pins #1, 3, 4, 5, 6, and 8 are connected to six digital output signals from the System FPGA of the Pixie-16 module, i.e. FO0, FO2, FO3, FO4, FO6, FO7, mainly for the purpose of testing and debugging. Pins #2 and 7 are ground pins, and pins #9 and 10 are not in use.

## GATE Inputs(Rev. D only)

Connector Type	Amphenol FCI® 55 Position Header, 2mm pin spacing
FI <sub>0</sub> , FI <sub>2</sub> , FI <sub>3</sub> , FI <sub>4</sub> , FI <sub>6</sub> , FI <sub>7</sub>	TTL digital input signals (max. 5V)
Gx <sub>in+</sub> /Gx <sub>in-</sub> (x=0-15)	Channel Gate Inputs (0-15 for 16 channels) (LVDS format)
MG <sub>in+</sub> /MG <sub>in-</sub>	Module Gate Input (LVDS format)
SG <sub>in+</sub> /SG <sub>in-</sub>	Spare Gate Input (LVDS format)

Figure: Rev. D Module's GATE Inputs

On Rev. D Pixie-16 modules, between analog input SMB connectors for channel 11 and channel 12, respectively, is the GATE INPUTS connector. This connector is an Amphenol FCI® 55 Position Header with 2mm pin spacing. The layout of these 55 pins is shown in Figure *Rev. D Module's GATE INPUTS Connector*. The 11 pins from the middle pin column (J150C) are all tied to the Ground. Among the first 8 rows of the GATE INPUTS connector, each differential pair of pins from the A/B columns(J150A/J150B) or the D/E columns (J150D/J150E) corresponds to one channel's GATE INPUT, which has the LVDS format, e.g. Gx<sub>in+</sub>/Gx<sub>in-</sub>(x=0-15). Differential pair of pins at J150A3/J150B3 is the Module Gate Input signal, MG<sub>in+</sub>/MG<sub>in-</sub>. Channel Gate Input signal can be used to veto or validate that given channel's own trigger. Module Gate Input signal works on the whole module level, i.e. it can be used to veto or validate all 16 channels' own trigger of that given module. Differential pair of pins at J150D3/J150E3 is the Spare Gate Input signal, SG<sub>in+</sub>/SG<sub>in-</sub>. Spare Gate Input signal can be used for special applications which require a custom firmware.

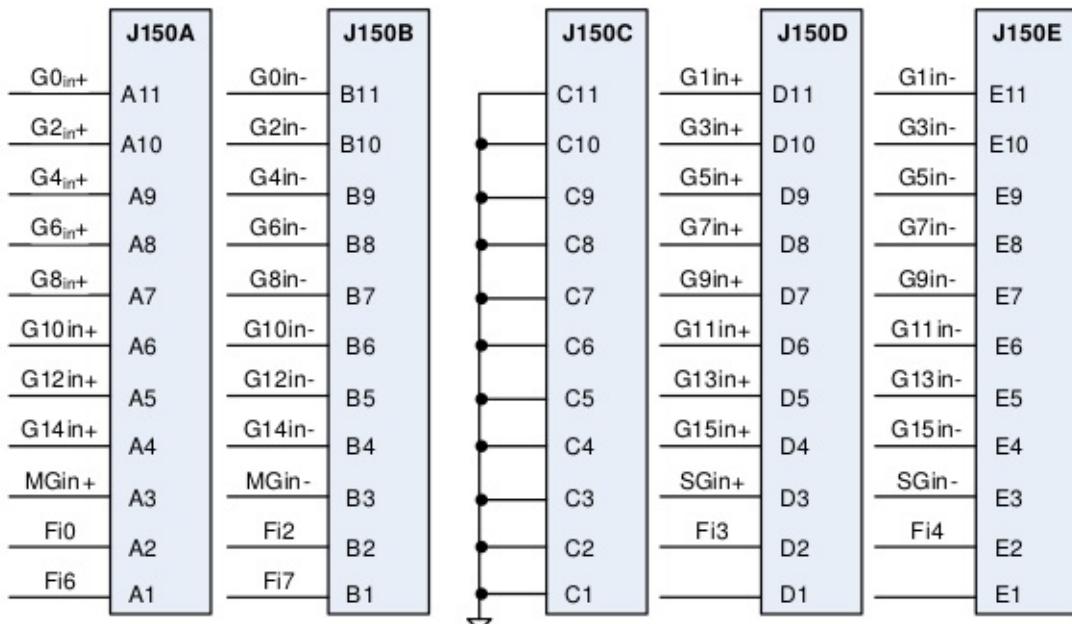
Figure: Rev. D Module's GATE INPUTS Connector

On Rev. D Pixie-16 modules, the TTL digital input signals (max. 5V), i.e. FI<sub>0</sub>, FI<sub>2</sub>, FI<sub>3</sub>, FI<sub>4</sub>, FI<sub>6</sub>, FI<sub>7</sub>, are distributed among the bottom two rows of the GATE INPUTS Connector, as illustrated in Figure Rev. D Module's GATE INPUTS Connector.

### 3.3V I/O Connector(Rev. B and C only)

Connector Type	Single-ended, 2mm pin spacing
FI <sub>0</sub> , FI <sub>2</sub> , FI <sub>3</sub> , FI <sub>4</sub> , FI <sub>6</sub> , FI <sub>7</sub>	TTL digital input signals (max. 5V)
FO <sub>0</sub> , FO <sub>2</sub> , FO <sub>3</sub> , FO <sub>4</sub> , FO <sub>6</sub> , FO <sub>7</sub>	Digital outputs for test/debug purpose (3.3V)

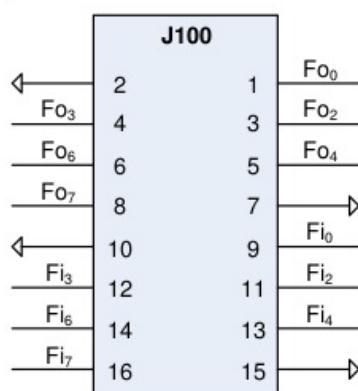
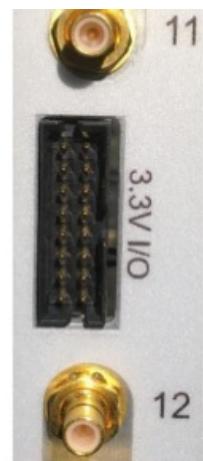
Figure: Rev. B and C Module's 3.3V I/O Connector

On Rev. B and C Pixie-16 modules, between analog input SMB connectors for channel 11 and channel 12, respectively, is the 3.3V I/O Connector(J100). It has 16 single-ended pins with 2mm spacing. Pins #1, 3, 4, 5, 6, and 8 are connected to six digital output signals from the System FPGA of the Pixie-16 module, i.e. FO<sub>0</sub>, FO<sub>2</sub>, FO<sub>3</sub>, FO<sub>4</sub>, FO<sub>6</sub>, FO<sub>7</sub>, mainly for the purpose of

testing and debugging. Pins #2, 7, 10 and 15 are ground pins. Pins #9, 11, 12, 13, 14 and 16 are connected to the six TTL digital input signals (max. 5V), i.e. FI0, FI2, FI3, FI4, FI6, FI7.

## Digital Signals in Standard Firmware(all revisions)

The standard firmware of the Pixie-16 supports input and output of digital signals through its front panel I/O connectors, which were discussed earlier.

TTL digital input signals	Connected signals in standard firmware	Direction	Description
FI <sub>0</sub>	EXT_FASTTRIG	Input	External fast trigger signal
FI <sub>2</sub>	INHIBIT	Input	Run inhibit signal
FI <sub>3</sub>	EXT_TS_CLK	Input	External timestamp clock signal
FI <sub>4</sub>	EXT_VALIDTRIG	Input	External validation signal
FI <sub>6</sub>	not used		
FI <sub>7</sub>	EXT_TS_CLR	Input	External timestamp clear signal

*Figure: TTL Digital Input Signals*

Table *TTL Digital Input Signals* shows the five TTL digital input signals supported by the Pixie-16 standard firmware.

Among them, the signals **EXT\_TS\_CLK** and **EXT\_TS\_CLR** are used for external timestamping in the Pixie-16, i.e. the Pixie-16 accepting an external clock signal(the frequency of this external clock is not recommended to exceed about 20 MHz in order to avoid the clock signal integrity issue), counting such clock signal with a 48-bit counter, and outputting such counter value to the list mode data stream when an event trigger occurs.

The external timestamping is useful for synchronizing the Pixie-16 data acquisition system with another data acquisition system through correlating the external timestamps of the events recorded by both systems.

The **INHIBIT** signal is used by an external system to inhibit the data acquisition run in a Pixie-16 system when synchronization requirement is enabled in the Pixie-16 modules. It is a level sensitive signal, i.e. when the **INHIBIT** signal is at the logic high level, the run in the Pixie-16 won't start. Only when the **INHIBIT** signal goes to the logic low level will the run start in the Pixie-16. During the run, if the **INHIBIT** signal returns to the logic high level, the run will be aborted.

The **EXT\_FASTTRIG** signal is the external fast trigger signal, which can be used to replace the local fast trigger for recording events in the Pixie-16 modules. The **EXT\_VALIDTRIG** signal is the external validation signal, which can be used to validate events in the Pixie-16 modules.

Connector J101 Pins	Connected signals in standard firmware	Direction	Description
F <sub>01p</sub> /F <sub>01n</sub>	not used		
F <sub>i1p</sub> /F <sub>i1n</sub>	LVDS_VALIDTRIG	Input	External validation trigger signal in LVDS format
F <sub>i5p</sub> /F <sub>i5n</sub>	LVDS_FASTTRIG	Input	External fast trigger signal in LVDS format
F <sub>o5p</sub> /F <sub>o5n</sub>	SYNC_LVDS_FP	Output	Pixie-16 synchronization output signal in LVDS format (to synchronize with other DAQ systems)

Figure: Connector J101 LVDS I/O Port Signals

Table *Connector J101 LVDS I/O Port Signals* shows the Pixie-16 connector J101 LVDS I/O port signals. This J101 LVDS I/O port can use the regular Ethernet cable for connection but it does not have Ethernet connectivity. Among the four LVDS pairs available from this J101 port, one pair is currently not in use, two pairs are used for input and one pair is used for output. The **LVDS\_VALIDTRIG** is the external validation trigger signal in LVDS format, and the **LVDS\_FASTTRIG** is the external fast trigger signal in LVDS format. The **SYNC\_LVDS\_FP** is an output signal from the Pixie-16 module to indicate to external data acquisition systems the synchronization status of the Pixie-16 system so that both data acquisition systems can be synchronized.

TTL digital output signals	Connected signals in standard firmware		Direction	Description	
F <sub>O<sub>0</sub></sub>	FTRIG_DELAY	FTRIG_DELAY	Output	Delayed local fast trigger of one of the 16 channels	Delayed local fast trigger of one of the 16 channels
F <sub>O<sub>2</sub></sub>	FTRIG_VAL	VETO_CE	Output	Validated, delayed local fast trigger one of the 16 channels	Stretched veto trigger of one of the 16 channels
F <sub>O<sub>3</sub></sub>	ETRIG_CE	LDPMFULL	Output	Stretched external global validation trigger of one of the 16 channels	Module level dual port memory (DPM) full status flag
F <sub>O<sub>4</sub></sub>	CHANTRIG_CE	SDPMFULL	Output	Stretched channel validation trigger of one of the 16 channels	System level dual port memory (DPM) full status flag
F <sub>O<sub>6</sub></sub>	FTIN_OR	FTIN_OR	Output	OR of 16 local fast triggers	OR of 16 local fast triggers
F <sub>O<sub>7</sub></sub>	TEST_SEL	TEST_SEL	Output	Selected test signal	Selected test signal

Figure: TTL Digital Output Signals

Table *TTL Digital Output Signals* lists the six Pixie-16 TTL digital output signals. Two groups of six output signals can be chosen through software settings (see bits [14:12] and [19:16] of TrigConfig0). The last output signal TEST\_SEL can be further selected through software settings. More details about these signals will be provided in later sections of this manual.

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© Hongyi Wu      *updated: 2018-11-03 15:08:17*

## Logic

- Module Fast Trigger (for trigger)
- Module Validation Trigger (for control logic)
- Channel Validation Trigger(for trigger/control logic)
- Veto
- System FPGA (coincidence/multiplicity)

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Control logic :

- Module validation trigger
- Channel validation trigger
- Veto
- Pileup
- ...

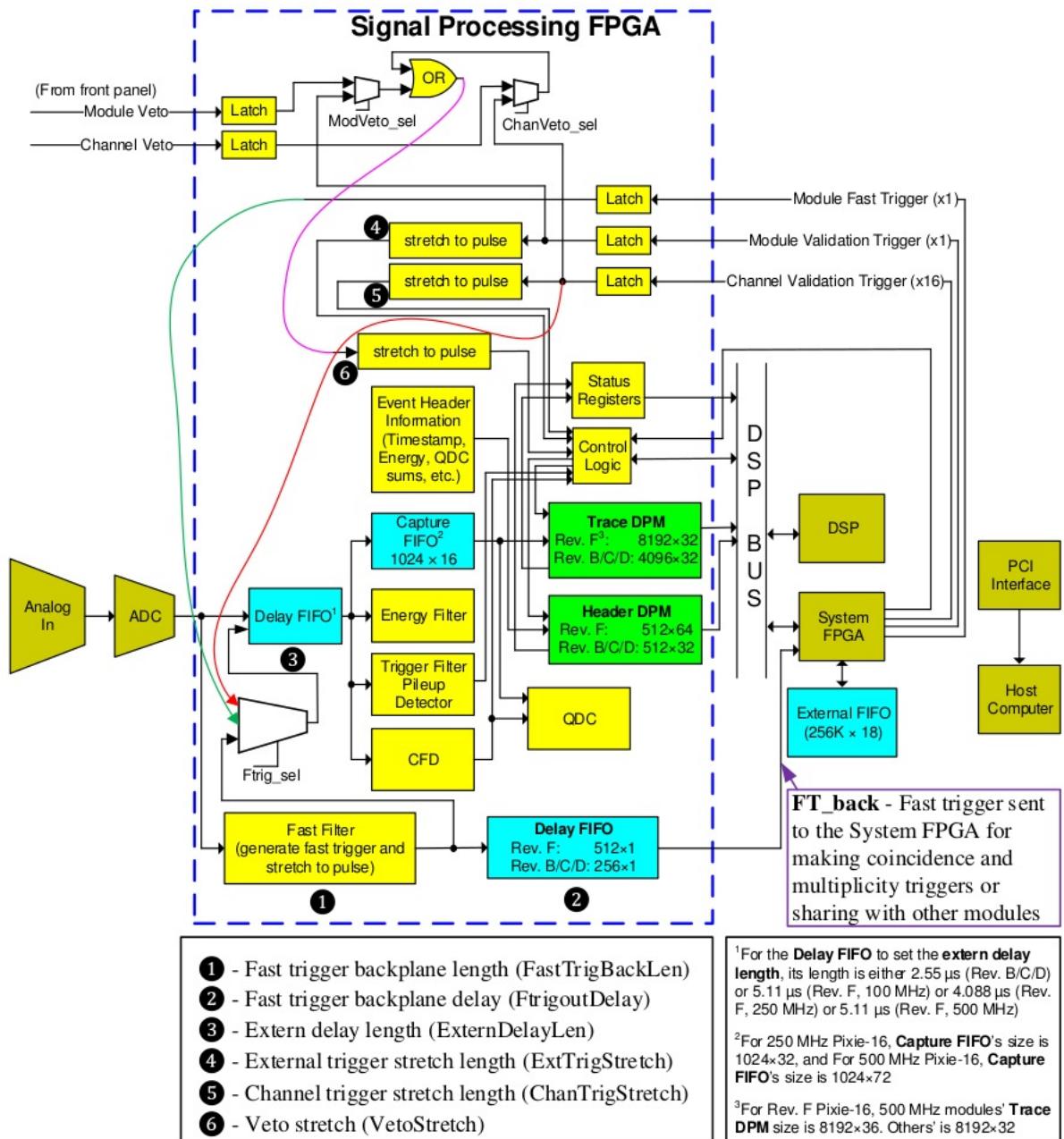
Figure: illustrates the signal processing in the Pixie-16 modules

As shown in Figure, the incoming analog pulse is first digitized by the ADC and then enters the signal processing circuitries in the Signal Processing FPGA, each of which processes ADC data from 4 channels of a Pixie-16 module.

The digitized data stream is first fed into two branches: a fast filter generating fast triggers to be sent to the System FPGA and a Delay FIFO which could be used to compensate for the delay between fast triggers and the external triggers.

The digitized data stream passing through the Delay FIFO is then branched into four parts:

- 1) energy filter which samples energy running sums at the PeakSample time;
- 2) trigger filter which detects pulse and performs pileup inspection;
- 3) capture FIFO which delays the ADC data according to the trace delay parameter value before the ADC data is streamed into the Trace Dual Port Memory (DPM) when a valid pulse is detected;
- 4) CFD circuitry where a CFD trigger is generated to trigger the computation of QDC sums, latch timestamps and record

traces.

The Control Logic in the signal processing FPGA utilizes the local fast trigger, CFD trigger, veto and external triggers to determine whether and when to stream waveform data into the Trace DPM and to write event information into the Header DPM. The DSP polls the status of the DPMs through the Status Registers and moves event data into the External FIFO through the DSP bus and the System FPGA.

### Trigger Stretch Lengths

- **External trigger stretch** is used to stretch the module validation trigger pulse.
- **Channel trigger stretch** is used to stretch the channel validation trigger pulse.
- **Veto stretch** is used to stretch the veto pulse for this channel.
- **Fast trigger backplane length** is used to stretch the fast trigger pulse to be sent to the System FPGA, where this fast trigger can be sent to the backplane to be shared with other modules, or can be used for making coincidence or multiplicity triggers.

### FIFO Delays

- **External delay length** is used to delay the incoming ADC waveform and the local fast trigger in order to compensate for the delayed arrival of the external trigger pulses, e.g., module validation trigger, channel validation trigger, etc.
- **Fast trigger backplane delay** is used to delay the fast trigger pulse before it is sent to the System FPGA for sharing with other modules through the backplane or making coincidence or multiplicity triggers.

Parameters	Range	
	100 or 500 MHz	250 MHz
<i>External trigger stretch</i>	0.01 – 40.95 µs	0.008 – 32.76 µs
<i>Channel trigger stretch</i>	0.01 – 40.95 µs	0.008 – 32.76 µs
<i>Veto stretch</i>	0.01 – 40.95 µs	0.008 – 32.76 µs
<i>Fast trigger backplane length</i>	0.01 – 40.95 µs	0.008 – 32.76 µs
<i>External delay length</i>	0 – 2.55 µs (Rev. B/C/D) 0 – 5.11 µs (Rev. F)	0 – 4.088 µs
<i>Fast trigger backplane delay</i>	0 – 2.55 µs (Rev. B/C/D) 0 – 5.11 µs (Rev. F)	0 – 4.088 µs

Figure: Range for Trigger Stretch Lengths and FIFO Delays in Pixie-16 Modules

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## Channel Validation Trigger(for trigger/control logic)

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## System FPGA (coincidence/multiplicity)

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Control logic (module/channel validation trigger)

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© Hongyi Wu      updated: 2018-11-03 15:08:17

## Multiple Modules Synchronously

- Individual Clock Mode
- PXI Clock Mode
- Daisy-chained Clock Mode
- Multi-Crate Clock Mode
  - Installation of Pixie-16 Modules
  - Clock Jumper (JP101) Settings on the Pixie-16 Modules
  - Cable Connections for Pixie-16 Rear I/O Trigger Modules
  - Jumper Settings on the Pixie-16 Rear I/O Trigger Modules

When many Pixie-16 modules are operated together as a system, it may be required to synchronize clocks and timers between them and to distribute triggers across modules. It will also be necessary to ensure that runs are started and stopped synchronously in all modules. All these signals are distributed through the PXI backplane of the Pixie-16 crate.

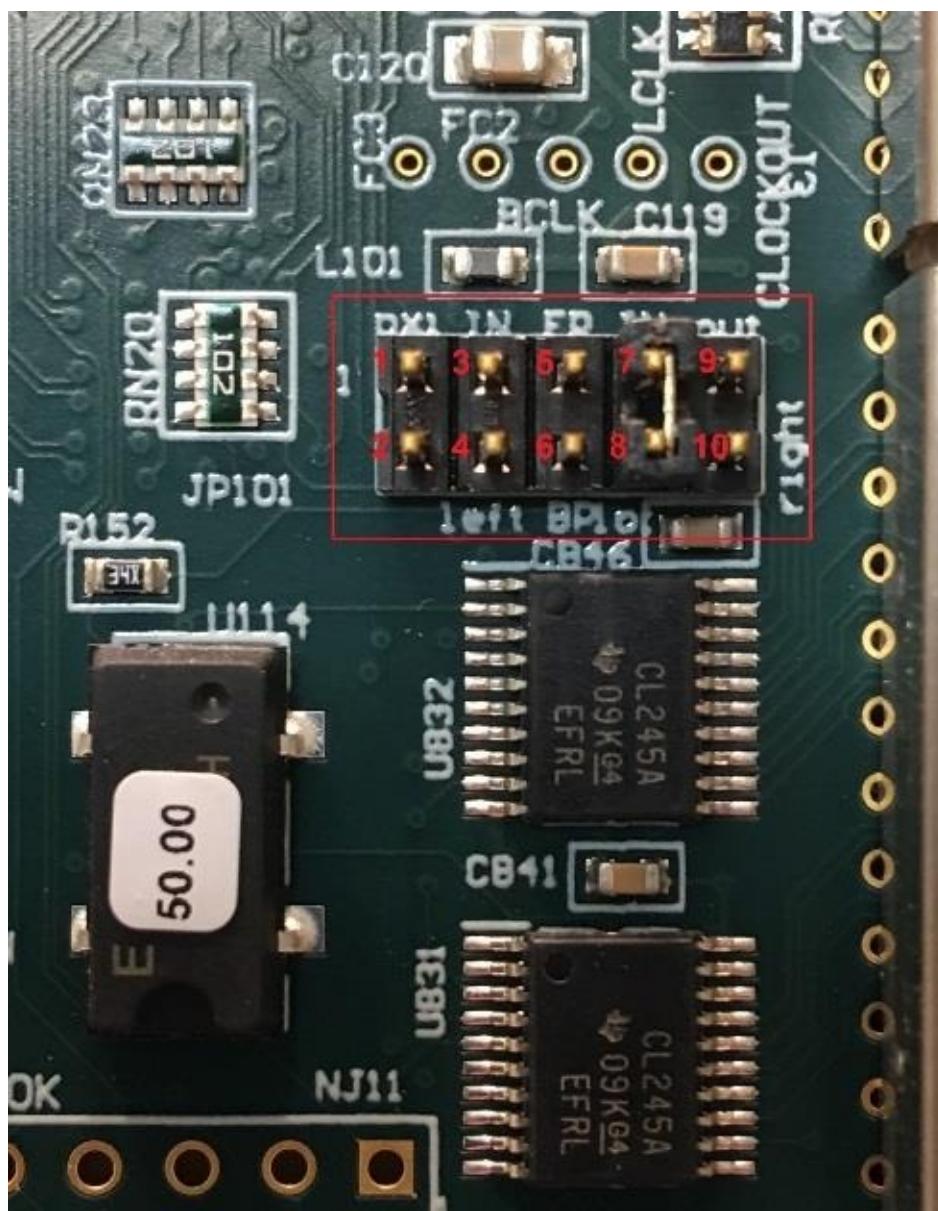
Figure: JP101

In a multi-module system there will be one clock master and a number of clock slaves or repeaters. The clock function of a module can be selected by setting shunts on Jumper JP101 near the bottom right corner of the board. The 10-pin Jumper JP101 is shown in the picture on the top with those pins labelled in red color. Shunts are provided to connect pins that are appropriate for each chosen clock distribution mode. Four clock distribution modes, individual clock mode, PXI clock mode, daisy-chained clock mode, and multi- crate clock mode, are described below.

**[warning] Please Note**

In 250 MHz or 500 MHz Pixie-16 modules, the frequency of signal processing clock in the FPGA has been divided down to either 125 MHz or 100 MHz, respectively, for more practical implementation of the design. That division might result in different clock phase and thus different timestamp offset for each channel within a given 250 MHz or 500 MHz Pixie-16 module whenever the module is reinitialized. Calibration might be needed to quantify the different timestamp offset for each channel.

## Individual Clock Mode


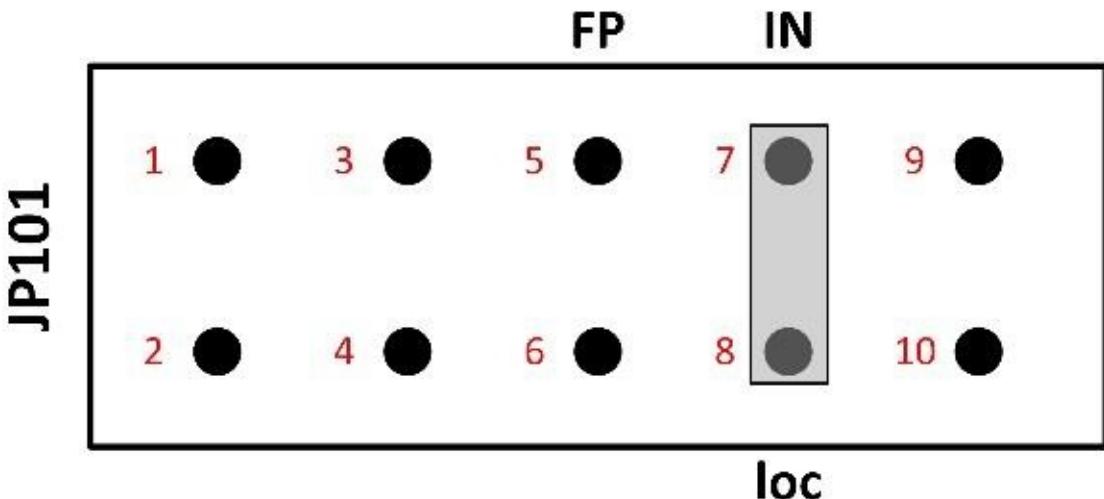
Figure: individual clock mode

If only one Pixie-16 module is used in the system, or if clocks between modules do not need to be synchronized, the module(s) should be set into individual clock mode. Connect pin 7 of JP101 (the clock input) with a shunt to pin 8 (loc – IN). This will use the 50 MHz local crystal oscillator of the Pixie-16 module as the clock source.

## PXI Clock Mode

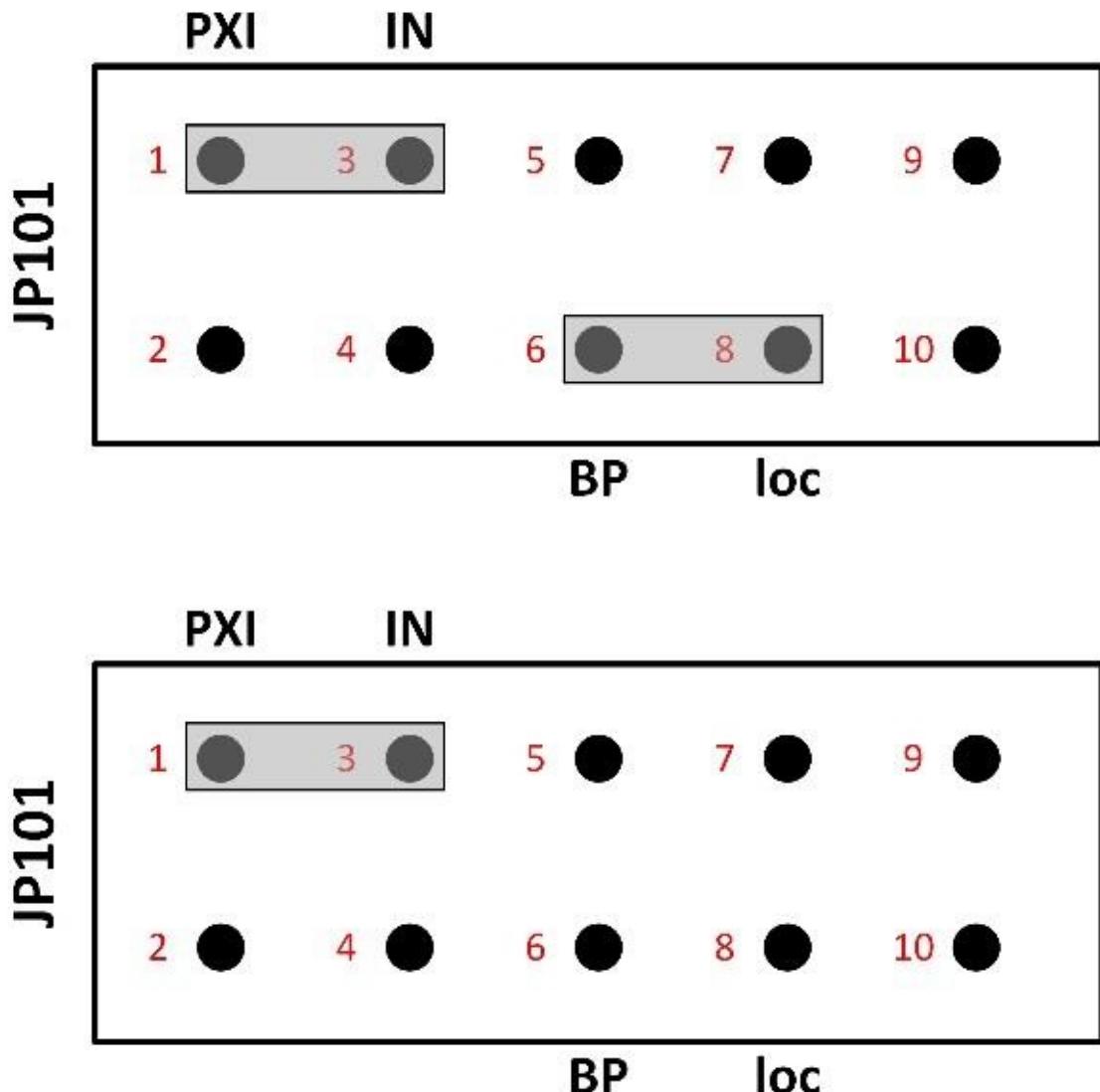Top: PXI clock master (slot 2)  
Bottom: PXI clock recipient (slots 3-14)

*Figure: PXI clock mode*


The preferred way to distribute clocks among multiple Pixie-16 modules is to use the PXI clock distributed on the backplane. This clock is by default generated on the backplane and is a 10MHz clock signal, which is then repeated by a fan out buffer and connected to each crate slot by a dedicated line with minimum skew(equal trace length to each slot). Although the 10MHz is too slow to be a useful clock for the Pixie-16, it can be overridden by a local clock signal from a Pixie-16 module that is installed in slot 2 through proper shunt settings on the JP101.

A Pixie-16 module can be configured to be the PXI clock master in slot 2 by connecting pins 6 and 8 (loc – BP) of the JP101. All modules, including the clock master, should be set to receive the PXI clock by connecting pin 1 and 3 on JP101 (PXI – IN). In this way, the 50 MHz clock from the Pixie-16 clock master is distributed to all Pixie-16 modules through the backplane with nearly

identical clock phase.

One other advantage of the PXI clock mode over the daisy-chained clock mode, which will be discussed next, is that except for the Pixie-16 master module, which has to be installed in slot 2, other Pixie-16 slave modules can be installed in any other slot of the Pixie-16 crate. In contrast, when the daisy-chained clock mode is used, all Pixie-16 modules have to be installed next to each other, i.e. no gap is allowed between modules.

## Daisy-chained Clock Mode



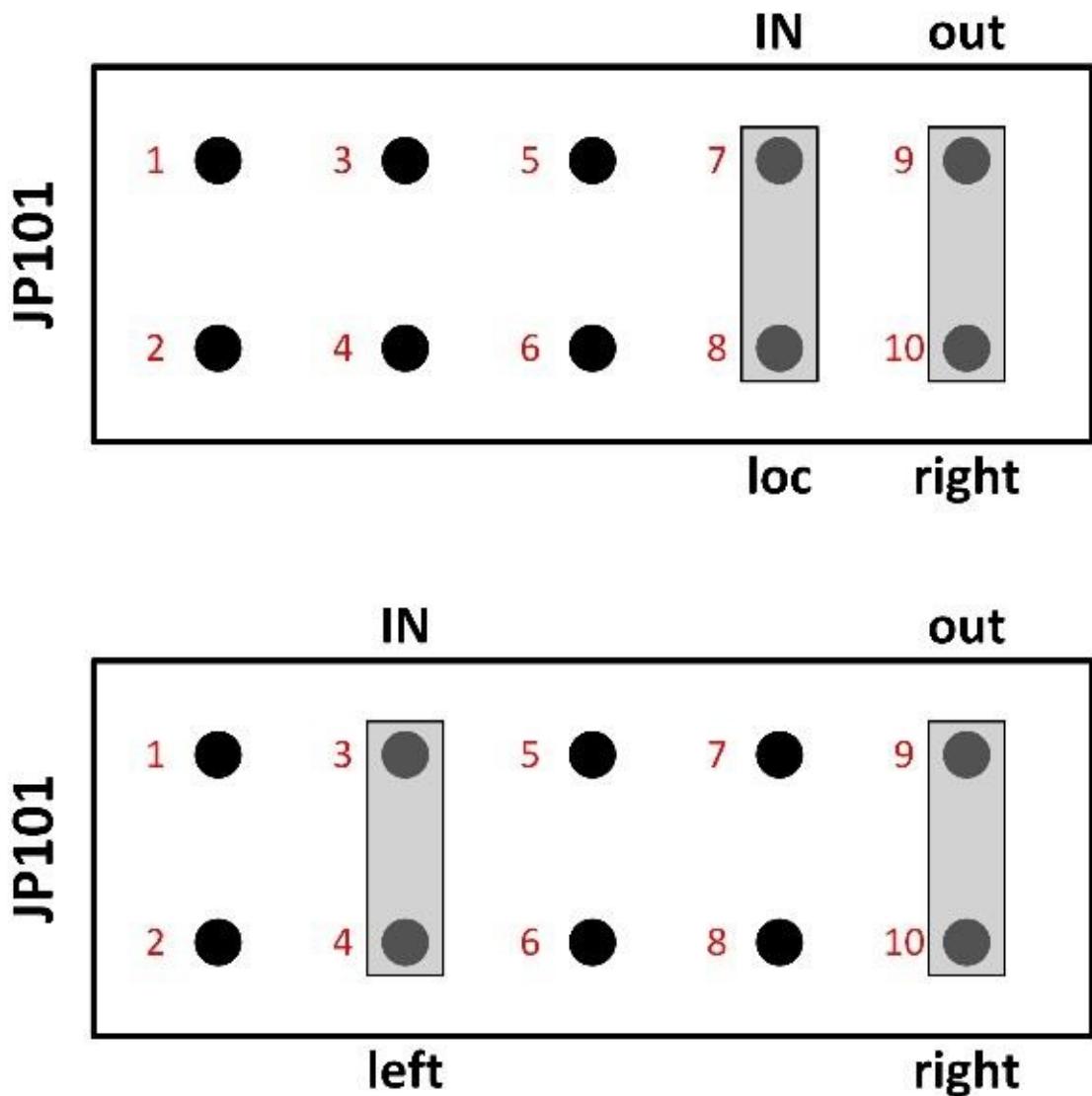
Top: Daisy-chain clock master (leftmost module)  
Bottom: Daisy-chain clock repeater

Figure: daisy-chained clock mode

A further option for clock distribution is to daisy-chain the clocks from one module to the other, with each module repeating the clock signal and transmitting it to the neighbor on the right. This requires one master module, located in the leftmost slot of the group of Pixie-16 modules. The master module uses its local crystal oscillator as the input and sends its output to the right (loc – IN, out – right). Other Pixie-16 modules in the crate should be configured as clock repeaters by using the signal from the left neighbor as the input and sending its output to the right (left – IN, out – right). However, as mentioned earlier, there must be no slot gap between modules.

---

## Multi-Crate Clock Mode

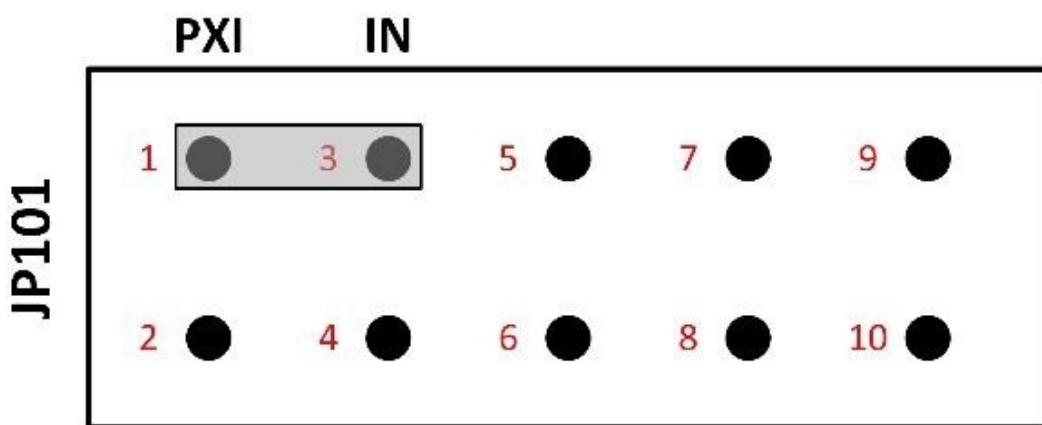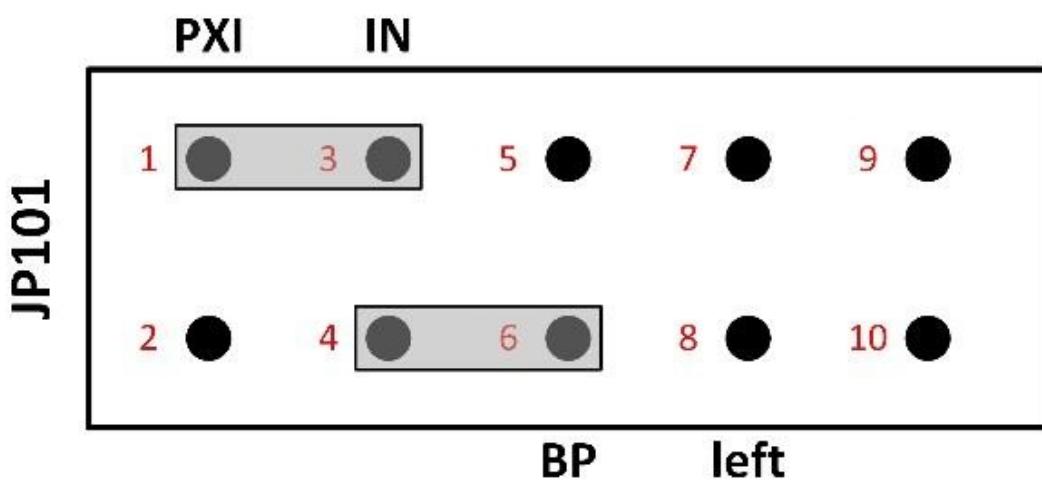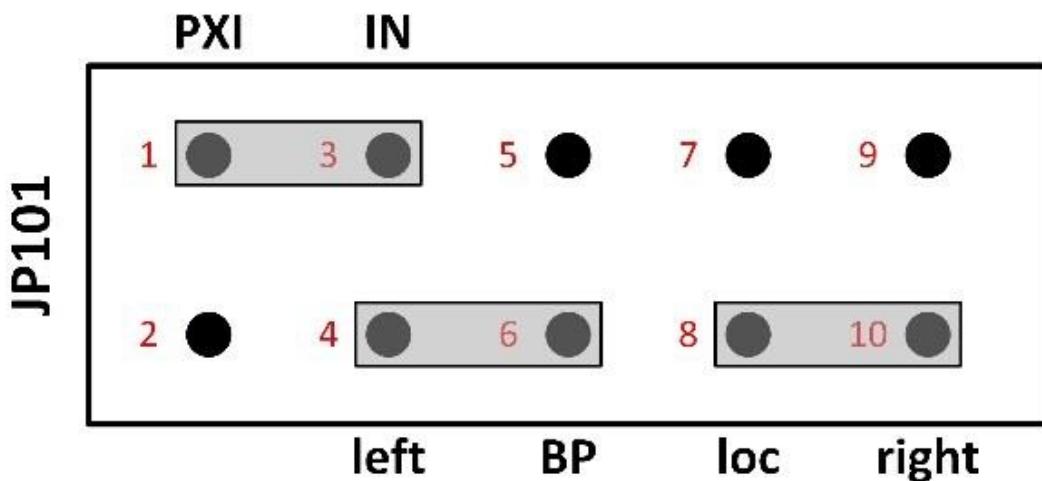Top: System Director Module

Middle: Crate Master Module

Bottom: Regular Module

*Figure: multi-crate clock mode*

In multi-crate systems, a global clock signal can be distributed among these crates using dedicated trigger and clock distribution cards, i.e. the Pixie-16 Rear I/O trigger modules, which are available from XIA.

*An example of clock distribution between two crates is illustrated below.*

## Installation of Pixie-16 Modules

Multiple Pixie-16 modules can be installed in two 14-slot Pixie-16 crates, #1 and #2. For clock distribution purpose, crate #1 is called the Master crate, where the system-wide global clock for all Pixie-16 modules is originated, and crate #2 is called the Slave crate, which receives the global clock from the Master crate.

The Pixie-16 module installed in slot 2 of the Master crate is designated as the System Director Module, whose local 50 MHz crystal oscillator acts as the source of the system-wide global clock. The distribution of the clock signal from the System Director Module to all Pixie-16 modules in the 2-crate system is done through the Pixie-16 Rear I/O trigger modules.

The Pixie-16 module installed in slot 2 of the Slave crate is called the crate Master module, which is responsible for receiving the global clock from the Master crate and sending such clock to all modules in that crate through length-matched traces on the backplane. The System Director Module is also responsible for sending the global clock to all modules in the Master crate. Therefore, it is also a crate Master module. Other modules in these two crates are regular modules. Table shows the different types of modules in a 2-crate system.

Crate #	1				
Slot #	2	3	...	13	14
Module	System Director Module / Crate Master Module	Regular Module	...	Regular Module	Regular Module
Crate #	2				
Slot #	2	3	...	13	14
Module	Crate Master Module	Regular Module	...	Regular Module	Regular Module

*Figure: Module Definitions in a 2-crate System*

## Clock Jumper (JP101) Settings on the Pixie-16 Modules

For all Pixie-16 modules in a 2-crate system to use the same global clock signal, the clock jumper (JP101) in all modules should be set according to Table *Clock Jumper JP101 Settings in a 2-crate System* and Figure *multi-crate clock mode*.


<b>System Director Module</b>	Connect pins 1 and 3, 4 and 6, 8 and 10.
<b>Crate Master Module</b>	Connect pins 1 and 3, 4 and 6.
<b>Regular Module</b>	Connect pins 1 and 3.

*Figure: Clock Jumper JP101 Settings in a 2-crate System*

## Cable Connections for Pixie-16 Rear I/O Trigger Modules

The Pixie-16 Rear I/O trigger modules are installed at the rear side of each crate where a 6U card cage is installed. Figure *rear I/O trigger modules* shows a Pixie-16 Rear I/O trigger module is installed directly behind either the Director or the Master module, respectively, to share clock, triggers, and run start or stop synchronization signals among multiple Pixie-16 crates. The rear of the backplane has connectors J3, J4 and J5, but it does not have J1 and J2, since it does not need to use CompactPCI or PXI communication.

Typically the first slot at the rear of the backplane with J3, J4, J5 connectors installed is the slot where the Pixie-16 Rear I/O trigger module should be installed. While installing the module, please ensure the alignment of top and bottom rails with the trigger module to avoid damage to the backplane pins.

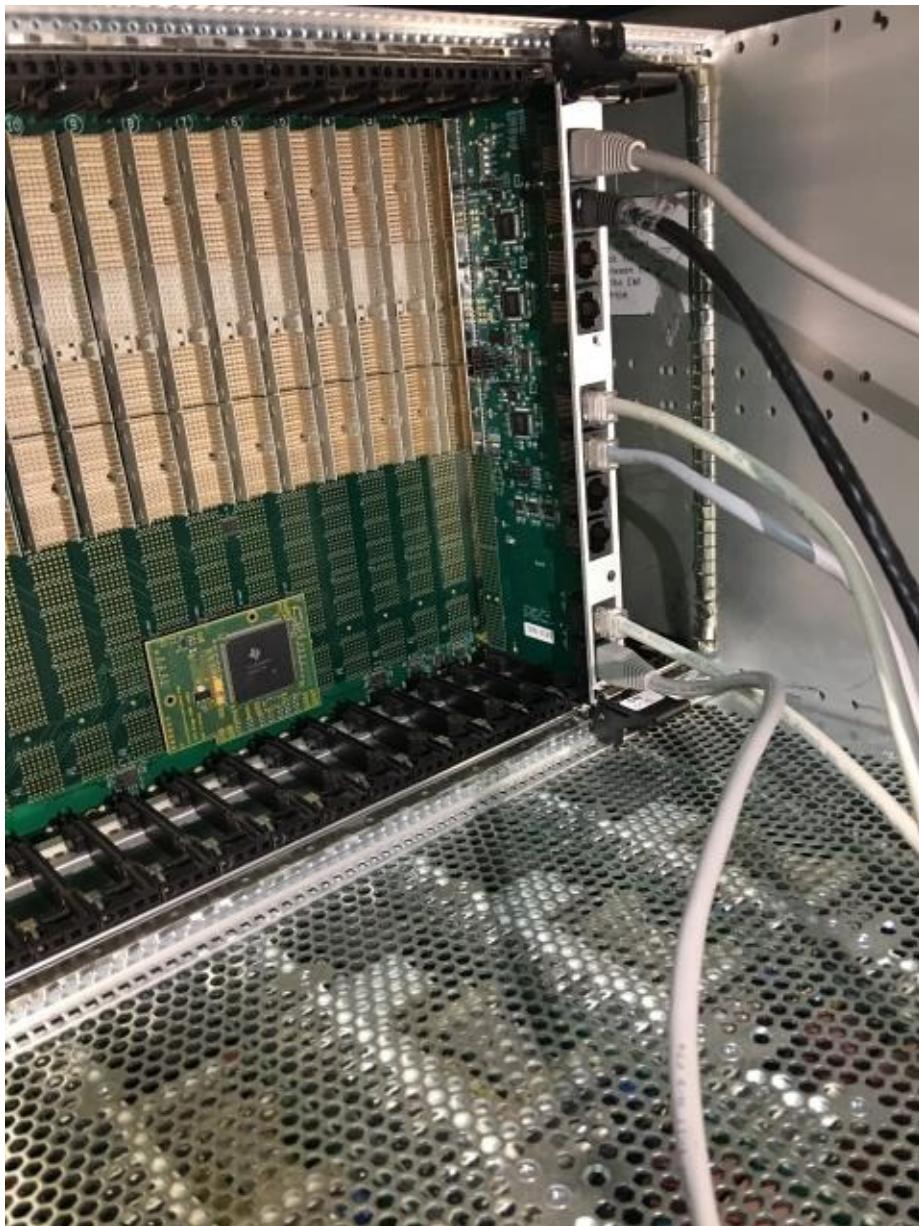*Figure: rear I/O trigger modules*

Figure *Cable connections between two Pixie-16 rear I/O trigger modules* shows the cable connections between two Pixie-16 rear I/O trigger modules that are installed in two separate crates. All connection cables are Category 5 or 6 Ethernet cables and shall have the same length to minimize clock phase difference between Pixie-16 modules in the two crates.

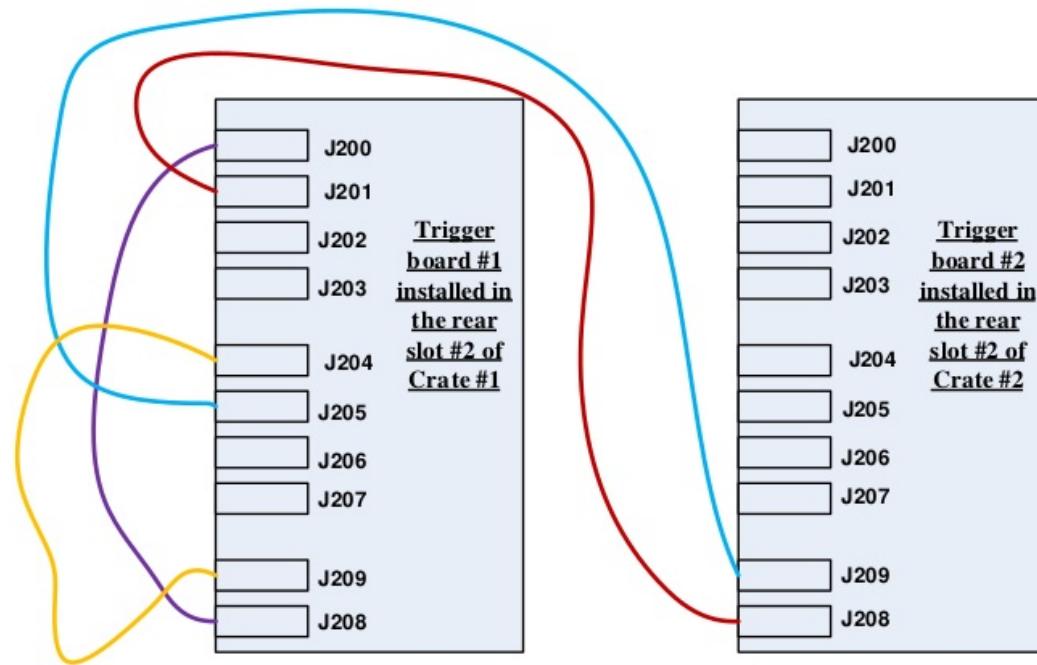
Figure: Cable connections between two Pixie-16 rear I/O trigger modules

## Jumper Settings on the Pixie-16 Rear I/O Trigger Modules

Trigger module #1 is installed in the rear slot #2 of crate #1. As mentioned earlier, the rear slot #2 is located at the back of the crate and is at the direct opposite side of the front slot #2 of the crate. Care should be taken when installing the trigger module into the rear slot #2 by avoiding bending any pins of the rear side of the backplane, since that could cause the 3.3V pin to be shorted to neighboring ground pin and thus damage the whole backplane.

Please note pin numbering for all jumpers on the trigger module is counted from right to left when facing the top side of the module, i.e. the backplane connectors J3 to J5 are on the left (only exception is JP1, which is in vertical orientation and should be counted from bottom to top). A tiny '1' label is painted on the right hand side of the jumpers, indicating pin 1. Figure *Pin numbering for the jumpers on the Pixie-16 rear I/O trigger module* shows the pin '1' in red boxes.

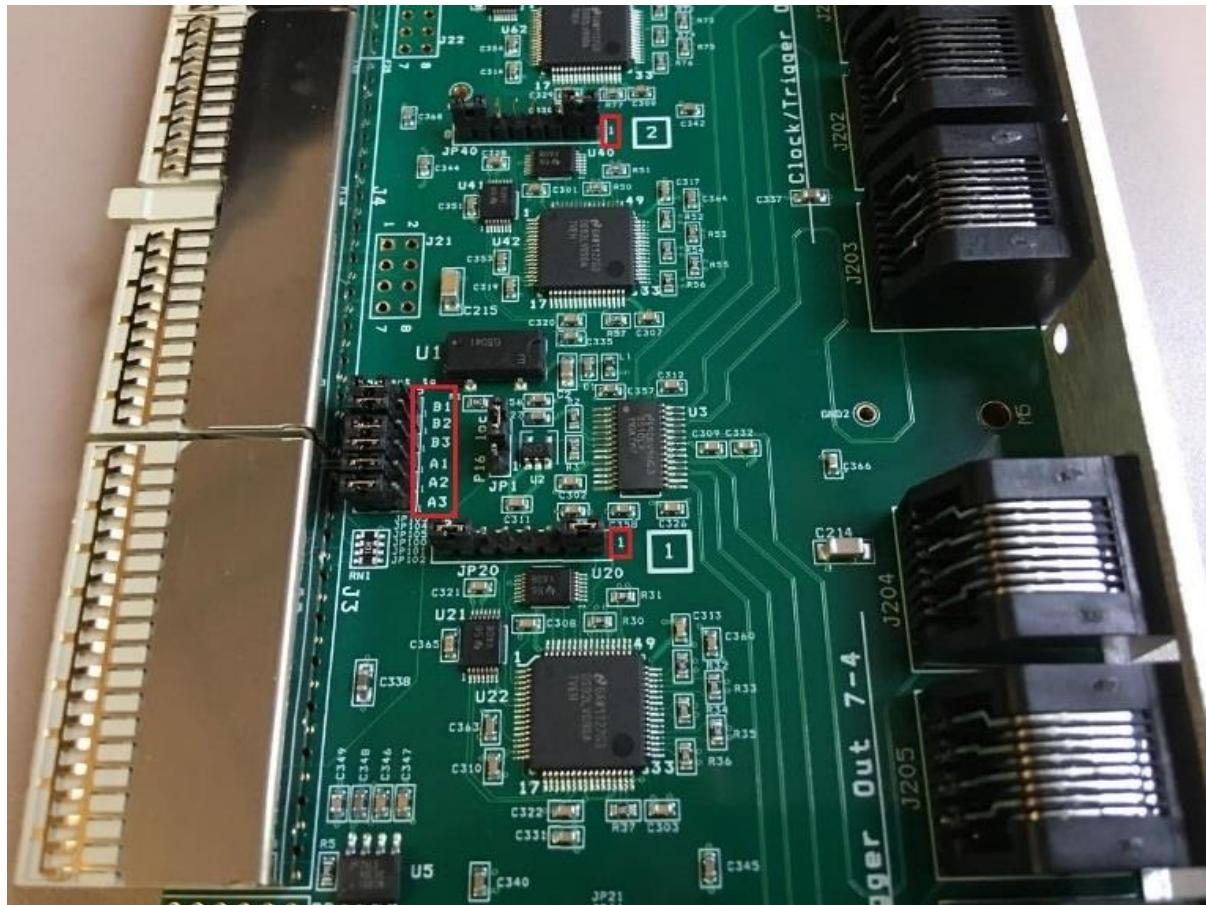
Figure: Pin numbering for the jumpers on the Pixie-16 rear I/O trigger module

Table Rear I/O Trigger Module #1's Jumper Settings shows the jumper settings of the Pixie-16 rear I/O trigger module #1 in a 2-crane system.

<b>JP1</b>	Connect pins 1 and 2 for "P16"
<b>JP20</b>	Connect pins 2 and 3, 6 and 7
<b>JP40</b>	Connect pins 2 and 3, 6 and 7
<b>JP60</b>	Connect pins 1 and 2, 7 and 8
<b>JP21</b>	Connect pins 2 and 3
<b>JP41</b>	Connect pins 2 and 3
<b>JP61</b>	Connect pins 1 and 2
<b>JP100</b>	Connect pins 2 and 3 (connect to J4)
<b>JP101</b>	Connect pins 2 and 3 (connect to J4)
<b>JP102</b>	Connect pins 2 and 3 (connect to J4)
<b>JP103</b>	Connect pins 2 and 3 (connect to J4)
<b>JP104</b>	Connect pins 2 and 3 (connect to J4)
<b>JP105</b>	Connect pins 2 and 3 (connect to J4)

Figure: Rear I/O Trigger Module #1's Jumper Settings

Trigger module #2 is installed in the rear slot #2 of crate #2. Table *Rear I/O Trigger Module #2's Jumper Settings* shows the jumper settings of the Pixie-16 rear I/O trigger module #2 in a 2-crate system.

<b>JP1</b>	Connect pins 2 and 3 for “loc”
<b>JP20</b>	Connect pins 1 and 2, 7 and 8
<b>JP40</b>	Connect pins 1 and 2, 7 and 8
<b>JP60</b>	Connect pins 2 and 3, 6 and 7
<b>JP21</b>	Don't connect any pin
<b>JP41</b>	Don't connect any pin
<b>JP61</b>	Connect pins 1 and 2
<b>JP100</b>	Connect pins 2 and 3 (connect to J4)
<b>JP101</b>	Connect pins 2 and 3 (connect to J4)
<b>JP102</b>	Connect pins 2 and 3 (connect to J4)
<b>JP103</b>	Connect pins 2 and 3 (connect to J4)
<b>JP104</b>	Connect pins 2 and 3 (connect to J4)
<b>JP105</b>	Connect pins 2 and 3 (connect to J4)

*Figure: Rear I/O Trigger Module #2's Jumper Settings*

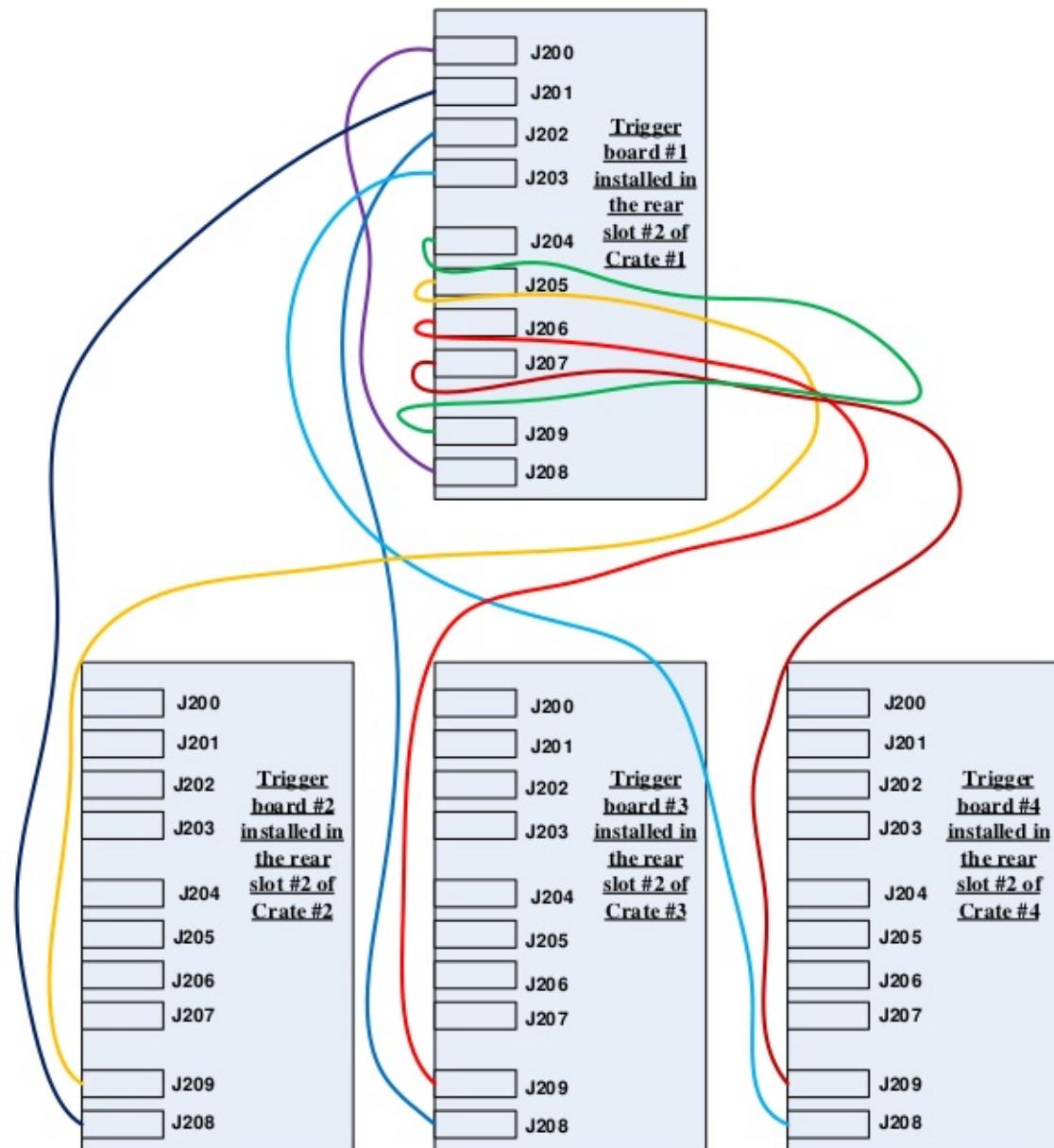
Figure: Cable connections among four Pixie-16 rear I/O trigger modules

Please note, if there are a total of four crates, the cable connections among those four Pixie-16 rear I/O trigger modules that are installed in those four separate crates should follow the connection methods shown in Figure 1-17. For the jumper settings on the Pixie-16 rear I/O trigger modules, trigger module #1 and #2 should use the same jumper settings as those in the trigger module #1 and #2 of the 2-crate system, respectively, whereas trigger module #3 and #4 should use the same jumper settings as those in trigger module #2.

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© Hongyi Wu      *updated: 2018-11-27 20:14:06*

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- 脉冲产生装置
- 100M module
  - 2 channel in one module
  - 2 channel in one crate
  - 2 channel in different crate
- 250M module
  - 2 channel in one module
  - 2 channel in one crate
  - 2 channel in different crate
- 100M & 250M module
  - 2 channel in one crate
  - 2 channel in different crate

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## LaBr3

### 250M

- fast filter
    - FL: 0.08
    - GF: 0.02
  - slow filter
    - SL: 0.14/0.13
    - SG: 0.05
    - Range: 1
    - TAU: 以实际测量为准(0.02)
  - cfd filter
    - dealy: 0.02
    - scale: 0
  - Trace
    - Delay: 1.496
    - Length: 3.504
  - QDC
    - 0.12,0.25,0.35,0.45,0.55,0.65,0.75,0.85
-

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© Hongyi Wu      *updated: 2018-11-03 15:08:17*

# XIA API

It from **Programmer's Manual Digital Gamma Finder (DGF) PIXIE-16 Version 1.40, October 2009**

```
// Configure modules for communication in PXI chassis
// Use this function to configure the Pixie-16 modules in the PXI chassis.
// NumModules is the total number of Pixie-16 modules installed in the system. PXISlotMap is the pointer to an
array that must have at least as many entries as there are Pixie-16 modules in the chassis.
// PXISlotMap serves as a simple mapping of the logical module number and the physical slot number that the mod
ules reside in. The logical module number runs from 0. For instance, in a system with 5 Pixie-16 modules, these
5 modules may occupy slots 3 through 7. The user must fill PXISlotMap as follows: PXISlotMap = {3, 4, 5, 6, 7
...} since module number 0 resides in slot number 3, etc. To find out in which slot a module is located, any pi
ece of subsequent code can use the expression PXISlotMap[ModNum], where ModNum is the logic module number.
// OfflineMode is used to indicate to the API whether the system is running in OFFLINE mode (1) or ONLINE mode
(0). OFFLINE mode is useful for situations where no Pixie-16 modules are present but users can still test their
calls to the API functions in their application software.
// This function must be called as the first step in the boot process. It makes the modules known to the system
and "opens" each module for communication.
// The function relies on an initialization file (pxisys.ini) that contains information about the Host PC's PCI
buses, including the slot enumeration scheme. XIA's software distribution normally puts this file under the sa
me folder as Pixie-16 software installation folder. However, the user has the flexibility of putting it in othe
r folders by simply changing the definition of the string PCISysIniFile_Windows or PCISysIniFile_Linux in the h
eader part of the file pixie16sys.c, depending on which operating system is being used.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16InitSystem (
    unsigned short NumModules,      // total number of Pixie16 modules in the system
    unsigned short *PXISlotMap,     // an array containing the PXI slot number for each pixie16 module
    unsigned short OfflineMode ); // specify if the system is in offline mode
```

```
// Release user virtual addressees assigned to modules
// Use this function to release the user virtual addressees that are assigned to Pixie-16 modules when these mo
dules are initialized by function Pixie16InitSystem. This function should be called before a user's application
exits.
// If ModNum is set to less than the total number of modules in the system, only the module specified by ModNum
will be closed. But if ModNum is equal to the total number of modules in the system, e.g. there are 5 modules
in the chassis and ModNum = 5, then all modules in the system will be closed altogether. Note that the modules
are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ExitSystem (
    unsigned short ModNum );      // module number
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadModuleInfo (
    unsigned short ModNum,        // module number
    unsigned short *ModRev,       // returned module revision
    unsigned int *ModSerNum,      // returned module serial number
    unsigned short *ModADCBits,   // returned module ADC bits
    unsigned short *ModADCMSPS ); // returned module ADC sampling rate
```

```
// Boot modules so that they can be set up for data taking
// Use this function to boot Pixie-16 modules so that they can be set up for data taking. The function download
s to the Pixie-16 modules the communication FPGA configurations, signal processing FPGA configurations, trigger
FPGA configurations (Revision A modules only), executable code for the digital signal processor (DSP), and DSP
parameters.
// The FPGA configurations consist of a fixed number of words depending on the hardware mounted on the modules;
the DSP codes have a length which depends on the actual compiled code; and the set of DSP parameters always co
nsists of 1280 32-bit words for each module. The host software has to make the names of those boot data files o
n the hard disk available to the boot function.
// If ModNum is set to be less than the total number of modules in the system, only the module specified by Mod
Num will be booted. But if ModNum is equal to the total number of modules in the system, e.g. there are 5 modul
es in the chassis and ModNum = 5, then all modules in the system will be booted.
```

```
// The boot pattern is a bit mask (shown below) indicating which on-board chip will be booted. Under normal circumstances, all on-board chips should be booted, i.e. the boot pattern would be 0x7F. For Rev-B, C, D modules, bit 1, i.e., "Boot trigger FPGA", will be ignored even if that bit is set to 1.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16BootModule (
    char *ComFPGAConfigFile,           // name of communications FPGA configuration file
    char *SPFFPGAConfigFile,          // name of signal processing FPGA configuration file
    char *TrigFPGAConfigFile,         // name of trigger FPGA configuration file
    char *DSPCodeFile,                // name of executable code file for digital signal processor (DSP)
    char *DSPParFile,                 // name of DSP parameter file
    char *DSPVarFile,                 // name of DSP variable names file
    unsigned short ModNum,            // pixie module number
    unsigned short BootPattern ); // boot pattern bit mask
```

```
// Acquire ADC traces in single or multiple modules
// Use this function to acquire ADC traces from Pixie-16 modules. Specify the module using ModNum. If ModNum is set to be less than the total number of modules in the system, only the module specified by ModNum will have its ADC traces acquired. But if ModNum is equal to the total number of modules in the system, then all modules in the system will have their ADC traces acquired.
// After the successful return of this function, the DSP's internal memory will be filled with ADC trace data. A user's application software should then call another function Pixie16ReadSglChanADCTrace to read the ADC trace data out to the host computer, channel by channel.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16AcquireADCTrace (
    unsigned short ModNum ); // module number
```

```
// Read ADC trace data from a channel in a module
// Use this function to read ADC trace data from a Pixie-16 module. Before calling this function, another function Pixie16AcquireADCTrace should be called to fill the DSP internal memory first. Also, the host code should allocate appropriate amount of memory to store the trace data. The ADC trace data length for each channel is 8192. Since the trace data are 16-bit unsigned integers (actually only the lower 14-bit contains real data due to the on-board 14-bit ADC), two consecutive 16-bit words are packed into one 32-bit word in the DSP internal memory. So for each channel, 4096 32-bit words are read out first from the DSP, and then each 32-bit word is unpacked to form two 16-bit words.
// Specify the module using ModNum and the channel on the module using ChanNum. Note that both the modules and channels are counted starting at 0.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglChanADCTrace (
    unsigned short *Trace_Buffer, // trace data
    unsigned int Trace_Length,   // trace length
    unsigned short ModNum,       // module number
    unsigned short ChanNum ); // channel number
```

```
// Transfer data between host and DSP internal memory
// Use this function to directly transfer data between the host and the DSP internal memory of a Pixie-16 module. The DSP internal memory is split into two blocks with address range 0x40000 to 0x4FFFF for the first block and address range 0x50000 to 0x5FFFF for the second block. Within the first block, address range 0x40000 to 0x49FFF is reserved for program memory and shouldn't be accessed directly by the host. Address range 0x4A000 to 0x4A4FF is used by the DSP I/O parameters which are stored in the configuration files (.set files) in the host. Within this range, 0x4A000 to 0x4A33F can be both read and written, but 0x4A340 to 0x4A4FF can only be read but not written. The remaining address range (0x4A500 to 4FFFF) in the first block and the entire second block (0x50000 to 0x5FFFF) should only be read but not written by the host. Use Direction = 1 for read and Direction = 0 for write.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16IMbufferIO (
    unsigned int *Buffer,             // buffer data
    unsigned int NumWords,            // number of buffer words to read or write
    unsigned int Address,             // buffer address
    unsigned short Direction,        // I/O direction
    unsigned short ModNum ); // module number
```

```
// Transfer data between host and DSP external memory
// Use this function to directly read data from or write data to the on-board external memory of a Pixie-16 module. The valid memory address is from 0x0 to 0x7FFFF (32-bit wide). Use Direction = 1 for read and Direction = 0 for write.
// The external memory is used to store the histogram data accumulated for each of the 16 channels of a Pixie-16 module. Each channel has a fixed histogram length of 32768 words(32-bit wide), and the placement of the histo
```

```

gram data in the memory is in the same order of the channel number, i.e. channel 0 occupies memory address 0x0
to 0xFFFF, channel 1 occupies 0x8000 to 0xFFFF, and so on.
// NOTE: another function Pixie16ReadHistogramFromModule can also be used to read out the histograms except tha
t it needs to be called channel by channel.
// In Rev-A modules, part of the external memory is also used to store the list mode data in ping-pong bufferin
g mode. This function can be used to read list mode data from the buffers.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16EMbufferIO (
    unsigned int *Buffer,           // buffer data
    unsigned int NumWords,          // number of buffer words to read or write
    unsigned int Address,           // buffer address
    unsigned short Direction,       // I/O direction
    unsigned short ModNum );        // module number

// Start a list mode data acquisition run
// Use this function to start a list mode data acquisition run in Pixie-16 modules. List mode run is used to co
llect data on an event-by-event basis, gathering energies, timestamps, pulse shape analysis values, and wavefor
ms, for each event. Runs will continue until a preset number of events are reached or the user terminates the r
un by calling function Pixie16EndRun. Once the run is progress, if the run is set to terminate after a given nu
mber of events have been accumulated, another function, Pixie16CheckRunStatus, should be called to check if the
run has finished. To start the data acquisition this function has to be called for every Pixie-16 module in th
e system. If all modules are to run synchronously, The last module addressed will release all others and the ac
quisition starts then. The first module to end the run will immediately stop the run in all other modules.
// Use mode=NEW_RUN (=1) to erase histograms and statistics information before launching the new run. Note that
this will cause a start up delay of up to 1 millisecond. Use mode=RESUME_RUN (=0) to resume an earlier run. Th
is mode has a start up delay of only a few microseconds.
// For Rev-A modules, currently there are 4 list mode run types supported. They are 0x100 (general purpose run)
, 0x101 (without waveforms), 0x102 (without auxiliary data) and 0x103 (energy and timestamp only).
// For Rev-B, C, D modules, there are only one list mode run type supported, that is, 0x100. However, different
output data options can be chosen by enabling or disabling different CHANCSRA bits.
// Histograms and statistics data are updated incrementally from run to run provided RESUME_RUN mode is used.
// ModNum is the module number which starts counting at 0. If ModNum is set to be less than the total number of
modules in the system, only the module specified by ModNum will have its list mode run started. But if ModNum
is set to equal to the total number of modules in the system, then all modules in the system will have their ru
ns started together.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16StartListModeRun (
    unsigned short ModNum,          // module number
    unsigned short RunType,         // run type
    unsigned short mode );          // run mode

// Start a MCA histogram mode data acquisition run
// Use this function to begin a data acquisition run that accumulates energy histograms, one for each channel.
It launches a data acquisition run in which only energy information is preserved and histogrammed locally to ea
ch channel.
// Call this function for each Pixie-16 module in the system. The last module addressed will allow the actual d
ata acquisition to begin. Histogram run can be self-terminating when the elapsed run time exceeds the preset ru
n time, or the user can prematurely terminate the run by calling Pixie16EndRun. On completion, final histogram
and statistics data will be available.
// Use mode=NEW_RUN (=1) to erase histograms and statistics information before launching the new run. Use mode=
RESUME_RUN (=0) to resume an earlier run.
// ModNum is the module number which starts counting at 0. If ModNum is set to be less than the total number of
modules in the system, only the module specified by ModNum will have its histogram run started. But if ModNum
is set to be equal to the total number of modules in the system, then all modules in the system will have their
runs started together.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16StartHistogramRun (
    unsigned short ModNum,          // module number
    unsigned short mode );          // run mode

// Check status of a data acquisition run
// Use this function to check the run status of a Pixie-16 module while a list mode data acquisition run is in
progress. If the run is still in progress continue polling.
// If the return code of this function indicates the run has finished, there might still be some data in the ex
ternal memory (Rev-A modules) or external FIFO (Rev-B, C, D modules) that need to be read out to the host. In a
ddition, final run statistics and histogram data are available for reading out too.
// In MCA histogram run mode, this function can also be called to check if the run is still in progress even th

```

```

ough it is normally self-terminating.
// ModNum is the module number which starts counting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16CheckRunStatus (
    unsigned short ModNum );           // Pixie module number

// Stop a data acquisition run
// Use this function to end a histogram run, or to force the end of a list mode run. In a multi-module system,
if all modules are running synchronously, only one module needs to be addressed this way. It will immediately s
top the run in all other module in the system.
// ModNum is the module number which starts counting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16EndRun (
    unsigned short ModNum );           // Pixie module number

// Compute input count rate
// Use this function to calculate the input count rate on one channel of a Pixie-16 module. This function does
not communicate with Pixie-16 modules. Before calling this function, another function, Pixie16ReadStatisticsFr
mModule, should be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statist
ics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. Mo
dNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeInputCountRate (
    unsigned int *Statistics,
    unsigned short ModNum,
    unsigned short ChanNum );

// Compute output count rate of a channel
// Use this function to calculate the output count rate on one channel of a Pixie-16 module. This function does
not communicate with Pixie-16 modules. Before calling this function, another function, Pixie16ReadStatisticsFr
omModule, should be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statist
ics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. Mo
dNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeOutputCountRate (
    unsigned int *Statistics,
    unsigned short ModNum,
    unsigned short ChanNum );

// Compute live time that a channel accumulated in a run
// Use this function to calculate the live time that one channel of a Pixie-16 module has spent on data acquisi
tion. This function does not communicate with Pixie-16 modules. Before calling this function, another function,
Pixie16ReadStatisticsFromModule, should be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statist
ics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. Mo
dNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeLiveTime (
    unsigned int *Statistics,
    unsigned short ModNum,
    unsigned short ChanNum );

// Compute number of events processed by a channel
// Use this function to calculate the number of events that have been processed by a Pixie-16 module during a d
ata acquisition run. This function is only used by Rev-A modules. This function does not communicate with Pixie
-16 modules. Before calling this function, another function, Pixie16ReadStatisticsFromModule, should be called
to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statist
ics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. Mo
dNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeProcessedEvents (
    unsigned int *Statistics,
    unsigned short ModNum );

```

```
// Compute real time that a channel accumulated in a run
// Use this function to calculate the real time that a Pixie-16 module has spent on data acquisition. This function does not communicate with Pixie-16 modules. Before calling this function, another function, Pixie16ReadStatisticsFromModule, should be called to read statistics data from the module first.
// *Statistics is a pointer to an array whose size is exactly 448 unsigned integer words (32-bit). The *Statistics array is filled with data from a Pixie-16 module after calling function Pixie16ReadStatisticsFromModule. ModNum is the module number which starts counting at 0. ChanNum is the channel number which starts counting at 0.
```

```
PIXIE16APP_EXPORT double PIXIE16APP_API Pixie16ComputeRealTime (
    unsigned int    *Statistics,
    unsigned short ModNum );
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16complexFFT (
    double *data,
    unsigned int length );

// Test one bit of a 16-bit unsigned integer
PIXIE16APP_EXPORT unsigned short PIXIE16APP_API APP16_TstBit (
    unsigned short bit,
    unsigned short value );

// Set one bit of a 16-bit unsigned integer
PIXIE16APP_EXPORT unsigned short PIXIE16APP_API APP16_SetBit (
    unsigned short bit,
    unsigned short value );

// Clear one bit of a 16-bit unsigned integer
PIXIE16APP_EXPORT unsigned short PIXIE16APP_API APP16_ClrBit (
    unsigned short bit,
    unsigned short value );

// Set one bit of a 32-bit unsigned integer
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API APP32_SetBit (
    unsigned short bit,
    unsigned int value );

// Clear one bit of a 32-bit unsigned integer
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API APP32_ClrBit (
    unsigned short bit,
    unsigned int value );

// Test one bit of a 32-bit unsigned integer
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API APP32_TstBit (
    unsigned short bit,
    unsigned int value );
```

```
// Program on-board DACs
// Use this function to reprogram the on-board digital to analog converters (DAC) of the Pixie-16 modules. In this operation the DSP uses data from the DSP parameters that were previously downloaded.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SetDACs (
    unsigned short ModNum );
```

```
// Program on-board signal processing FPGAs
// Use this function to program the on-board signal processing FPGAs of the Pixie-16 modules. After the host computer has written the DSP parameters to the DSP memory, the DSP needs to write some of these parameters to the FPGAs. This function makes the DSP perform that action.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ProgramFippi (
    unsigned short ModNum );
```

```
// Adjust DC-offsets in single or multiple modules
// Use this function to adjust the DC-offsets of Pixie-16 modules. Specify the module using ModNum. If ModNum is set to be less than the total number of modules in the system, only the module specified by ModNum will have its DC-offsets adjusted. But if ModNum is set to be equal to the total number of modules in the system, then al
```

```

1 modules in the system will have their DC-offsets adjusted.
// After the DC-offset levels have been adjusted, the baseline level of the digitized input signals will be determined by the DSP parameter BaselinePercent. For instance, if BaselinePercent is set to 10(%), the baseline level of the input signals will be ~ 1638 on the 14-bit ADC scale (minimum: 0; maximum: 16383).
// The main purpose of this function is to ensure the input signals fall within the voltage range of the ADCs to ensure all input signals can be digitized by the ADCs properly.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16AdjustOffsets (
    unsigned short ModNum );

```

```

// Acquire baselines from a module
// Use this function to acquire baselines from Pixie-16 modules. Specify the module using ModNum. If ModNum is set to be less than the total number of modules in the system, only the module specified by ModNum will have its baselines acquired. But if ModNum is set to be equal to the total number of modules in the system, then all modules in the system will have their baselines acquired.
// After the successful return of this function, the DSP's internal memory will be filled with baselines data. Users should then call another function Pixie16ReadSglChanBaselines to read the baselines data out to the host computer, channel by channel.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16AcquireBaselines (
    unsigned short ModNum );           // module number

```

```

// Read baselines from a channel in a module
// Use this function to read baselines data from a Pixie-16 module. Before calling this function, another function Pixie16AcquireBaselines should be called to fill the DSP internal memory first. Also, the host code should allocate appropriate amount of memory to store the baseline data. The baselines data length for each channel is 3640. In the DSP internal memory, each baseline data is a 32-bit IEEE floating point number. After being read out to the host, this function will convert each baseline data to a decimal number. In addition to baseline values, timestamps corresponding to each baseline were also returned after this function call.
// Specify the module using ModNum and the channel on the module using ChanNum. Note that the modules and channels are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglChanBaselines (
    double *Baselines,                // returned baselines values
    double *TimeStamps,               // time stamp for each baseline value
    unsigned short NumBases,          // number of baseline values to read
    unsigned short ModNum,             // module number
    unsigned short ChanNum );         // channel number

```

```

// Ramp Offset DACs of a module and record the baselines
// Use this function to execute the RAMP_OFFSETDACS control task run. Each Offset DAC has 65536 steps, and the RAMP_OFFSETDACS control task ramps the DAC from 0 to 65335 with a step size of 64, i.e., a total of 1024 steps. At each DAC step, the control task computes the baseline value as the representation of the signal baseline and stores it in the DSP memory. After the control task is finished, the stored baseline values are read out to the host computer and saved to a binary file called "rampdacs.bin" in the form of IEEE 32-bit floating point numbers. Users can then plot the baseline values vs. DAC steps to determine the appropriate DAC value to be set in the DSP in order to bring the input signals into the voltage range of the ADCs. However, this function is no longer needed due to the introduction of function Pixie16AdjustOffsets.
// If ModNum is set to less than the total number of modules in the system, only the module specified by ModNum will start the RAMP_OFFSETDACS control task run. But if ModNum is equal to the total number of modules in the system, e.g. there are 5 modules in the chassis and ModNum = 5, then all modules in the system will start the RAMP_OFFSETDACS control task run. Note that the modules are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16RampOffsetDacs (
    double *DCValues,                 // returned DC offset values
    unsigned short NumbCvals,          // number of DC values to read
    unsigned short ModNum );

```

```

// Execute special control tasks
// Use this function to call special control tasks. This may include programming the Fippi or setting the DACs after downloading DSP parameters.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ControlTaskRun (
    unsigned short ModNum,             // Pixie module number
    unsigned short ControlTask,        // Control task number
    unsigned int Max_Poll );          // Timeout control in unit of ms for control task run

```

```

// Find the Baseline Cut values of a module
// Use this function to find the Baseline Cut value for one channel of a Pixie-16 module. The baseline cut value is then downloaded to the DSP, where baselines are captured and averaged over time. The cut value would prevent a bad baseline value from being used in the averaging process, i.e., if a baseline value is outside the baseline cut range, it will not be used for computing the baseline average. Averaging baselines over time improves energy resolution measurement.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16BLcutFinder (
    unsigned short ModNum,           // Pixie module number
    unsigned short ChanNum,          // Pixie channel number
    unsigned int *BLcut );          // BLcut return value

// Find the exponential decay time of a channel
// Use this function to find the exponential decay time constant (Tau value) of the detector or preamplifier signal that is connected to one channel of a Pixie-16 module. The found Tau value is returned via pointer *Tau.

PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16TauFinder (
    unsigned short ModNum,           // Pixie module number
    double *Tau );                 // 16 returned Tau values, in #s

// Write a MODULE level parameter to a module
// Use this function to write a module parameter to a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16WriteSglModPar (
    char *ModParName,              // the name of the module parameter
    unsigned int *ModParData,       // the module parameter value to be written to the module
    unsigned short ModNum );       // module number

// Read a MODULE level parameter from a module
// Use this function to read a module parameter from a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglModPar (
    char *ModParName,              // the name of the module parameter
    unsigned int *ModParData,       // the module parameter value to be read from the module
    unsigned short ModNum );       // module number

// Write a CHANNEL level parameter to a module
// Use this function to write a channel parameter to a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16WriteSglChanPar (
    char *ChanParName,             // the name of the channel parameter
    double ChanParData,            // the channel parameter value to be written to the module
    unsigned short ModNum,         // module number
    unsigned short ChanNum );      // channel number

// Read a CHANNEL level parameter from a module
// Use this function to read a channel parameter from a Pixie-16 module.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadSglChanPar (
    char *ChanParName,             // the name of the channel parameter
    double *ChanParData,            // the channel parameter value to be read from the module
    unsigned short ModNum,         // module number
    unsigned short ChanNum );      // channel number

// Read histogram data from a module
// Use this function to read out the histogram data from a Pixie-16 module's histogram memory. Before calling this function, the host code should allocate appropriate amount of memory to store the histogram data. The default histogram length is 32768. Histogram data are 32-bit unsigned integers.
// Specify the module using ModNum and the channel on the module using ChanNum. Note that both the modules and channels are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadHistogramFromModule (
    unsigned int *Histogram,        // histogram data
    unsigned int NumWords,          // number of words to be read out
    unsigned short ModNum,          // module number
    unsigned short ChanNum);       // channel number

```

```
// Read run statistics data from a module
// Use this function to read out statistics data from a Pixie-16 module. Before calling this function, the host
// code should allocate appropriate amount of memory to store the statistics data. The number of statistics data
// for each module is fixed at 448. Statistics data are 32-bit unsigned integers.
// Specify the module using ModNum. Note that the modules are counted starting at 0.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadStatisticsFromModule (
    unsigned int *Statistics, // run statistics data
    unsigned short ModNum ); // module number
```

```
// Read histogram data from a module and save to a file
// Use this function to read histogram data from a Pixie-16 module and save the data to a file. New data will b
e appended to the end of the file. So the same file name can be used for multiple modules and the data from eac
h module will be stored in the order that this function is called.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SaveHistogramToFile (
    char *FileName, // histogram data file name
    unsigned short ModNum); // module number
```

```
// Parse a list mode data file to get events information
// Use this function to parse the list mode events in the list mode data file. The number of events for each mo
dule will be reported.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16GetModuleEvents (
    char *FileName, // the list mode data file name (with complete path)
    unsigned int *ModuleEvents ); // receives number of events for each module
```

```
// Get detailed events information from a data file
// Use this function to retrieve the detailed information of each event in the list mode data file for the desi
gnated module. Before calling this function to get the individual events information, another function Pixie16G
etModuleEvents should be called first to determine the number of events that have been recorded for each module
. If the number of events for a given module is nEvents, a memory block *EventInformation should be allocated w
ith a length of (nEvents*68):
```

```
// EventInformation = (unsigned long *)malloc(sizeof(unsigned long) * nEvents * 68);
// where 68 is the length of the information records of each event (energy, timestamps, etc.) and has the follo
wing structure.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16GetEventsInfo (
    char *FileName, // the list mode data file name (with complete path)
    unsigned int *EventInformation, // to hold event information
    unsigned short ModuleNumber); // the module whose events are to be retrieved
```

```
// Read trace data from a list mode data file
// Use this function to retrieve list mode trace from a list mode data file. It uses the trace length and file
location information obtained from function Pixie16GetEventsInfo for the selected event.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadListModeTrace (
    char *FileName, // list mode data file name
    unsigned short *Trace_Data, // list mode trace data (16-bit words)
    unsigned short NumWords, // number of 16-bit words to be read out
    unsigned int FileLocation); // the location of the trace in the file
```

```
// Read histogram data from a histogram data file
// Use this function to read histogram data from a histogram data file. Before calling this function, the host
// code should allocate appropriate amount of memory to store the histogram data. The default histogram length is
32768. Histogram data are 32-bit unsigned integers.
```

```
// Specify the module using ModNum and the channel on the module using ChanNum. Note that both the modules and
channels are counted starting at 0.
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadHistogramFromFile (
    char *FileName, // the histogram data file name (with complete path)
    unsigned int *Histogram, // histogram data
    unsigned int NumWords, // number of words to be read out
    unsigned short ModNum, // module number
    unsigned short ChanNum); // channel number
```

```
// Read DSP parameters from modules and save to a file
// Use this function to save DSP parameters to a settings file. It will first read the values of DSP parameters
on each Pixie-16 module and then write them to the settings file. Each module has exactly 1280 DSP parameter v
alues (32-bit unsigned integers), and depending on the value of PRESET_MAX_MODULES (defined in pixie16app_defs.
h), the settings file should have exactly (1280 * PRESET_MAX_MODULES * 4) bytes when stored on the computer har
d drive.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SaveDSPParametersToFile (
    char *FileName ); // the DSP parameters file name (with complete path)
```

```
// Load DSP parameters to modules from a file
// Use this function to read DSP parameters from a settings file and then download the settings to Pixie-16 mod
ules that are installed in the system. Each module has exactly 1280 DSP parameter values (32-bit unsigned integ
ers), and depending on the value of PRESET_MAX_MODULES (defined in pixie16app_defs.h), the settings file should
have exactly (1280 * PRESET_MAX_MODULES * 4) bytes when stored on the computer hard drive.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16LoadDSPParametersFromFile (
    char *FileName ); // the DSP parameters file name (with complete path)
```

```
// Copy DSP parameters from a module to others
// Use this function to copy DSP parameters from one module to the others that are installed in the system.
// BitMask is bit pattern which designates which items should be copied from the source module to the destinati
on module(s).
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16CopyDSPParameters (
    unsigned short BitMask, // copy items bit mask
    unsigned short SourceModule, // source module
    unsigned short SourceChannel, // source channel
    unsigned short *DestinationMask ); // the destination module and channel bit mask
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadMSGFile (
    char *ReturnMsgStr );
```

```
// Convert a decimal into IEEE 32-bit floating point number
PIXIE16APP_EXPORT unsigned int PIXIE16APP_API Decimal2IEEEFloating(double DecimalNumber);

// Convert an IEEE 32-bit floating point number to a decimal
PIXIE16APP_EXPORT double PIXIE16APP_API IEEEFloating2Decimal(unsigned int IEEEFloatingNumber);
```

```
// Read data from external FIFO and save to a file
// Use this function to read data from the external FIFO of a module. This function can only be used for Pixie-
16 Revision-B, C, and D modules.
// This function first checks the status of the external FIFO of a Pixie-16 module, and if there are data in th
e external FIFO, this function then reads list mode data (32-bit unsigned integers) from the external FIFO. So
this function essentially encapsulates both functions Pixie16CheckExternalFIFOStatus and Pixie16ReadDataFromExt
ernalFIFO within one function. The number of words that are read from the external FIFO is recorded in variable
*nFIFOWords.
// The function also expects setting the value of a variable called "EndOfRunRead" to indicate whether this rea
d is at the end of a run (1) or during the run (0). This is necessary since the external FIFO needs special tre
atment when the host reads the last few words from the external FIFO due to its pipelined structure.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16SaveExternalFIFODataToFile (
    char *FileName, // list mode data file name
    unsigned int *nFIFOWords, // number of words read from external FIFO
    unsigned short ModNum, // module number
    unsigned short EndOfRunRead); // indicator whether this is the end of run read
```

```
// Read from or write to registers on a module
// Use this function to read data from or write data to a register in a Pixie-16 module.
// Specify the module using ModNum. Note that the modules are counted starting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16RegisterIO (
    unsigned short ModNum, // the Pixie module to communicate to
    unsigned int address, // register address
```

```
    unsigned short direction,      // either MOD_READ or MOD_WRITE
    unsigned int   *value );      // holds or receives the data
```

```
// Read Control & Status Register value from a module
// Use this function to read the host Control & Status Register (CSR) value.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadCSR (
    unsigned short ModNum,
    unsigned int   *CSR );
```

```
// Write to Control & Status Register in a module
// Use this function to write a value to the host Control & Status Register (CSR).
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16WriteCSR (
    unsigned short ModNum,
    unsigned int   CSR );
```

```
// Check status of external FIFO of a module
// Use this function to check the status of the external FIFO of a Pixie-16 module while a list mode data acquisition run is in progress. The function returns the number of words (32-bit) that the external FIFO currently has. If the number of words is greater than a user-set threshold, function Pixie16ReadDataFromExternalFIFO can then be used to read the data from the external FIFO. The threshold can be set by the user to either minimize reading overhead or to read data out of the FIFO as quickly as possible.
// *nFIFOWords returns the number of 32-bit words that the external FIFO currently has.
// ModNum is the module number which starts counting at 0.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16CheckExternalFIFOStatus (
    unsigned int   *nFIFOWords,
    unsigned short ModNum );
```

```
// Read data from external FIFO of a module
// Use this function to read data from the external FIFO of a module. This function can only be used for Pixie-16 Revision-B, C, and D modules.
// This function reads list mode data from the external FIFO of a Pixie-16 module. The data are 32-bit unsigned integers. Normally, function Pixie16CheckExternalFIFOStatus is called first to see how many words the external FIFO currently has, then this function is called to read the data from the FIFO.
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ReadDataFromExternalFIFO (
    unsigned int   *ExtFIFO_Data, // To receive the external FIFO data
    unsigned int   nFIFOWords,   // number of words to read from external FIFO
    unsigned short ModNum );   // module number
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ComputeFastFiltersOffline (
    char          *FileName,           // the list mode data file name (with complete path)
    unsigned short ModuleNumber,       // the module whose events are to be analyzed
    unsigned short ChannelNumber,      // the channel whose events are to be analyzed
    unsigned int   FileLocation,        // the location of the trace in the file
    unsigned short RcdTraceLength,     // recorded trace length
    unsigned short *RcdTrace,          // recorded trace
    double         *fastfilter,        // fast filter response
    double         *cfd );            // cfd response
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API Pixie16ComputeSlowFiltersOffline (
    char          *FileName,           // the list mode data file name (with complete path)
    unsigned short ModuleNumber,       // the module whose events are to be analyzed
    unsigned short ChannelNumber,      // the channel whose events are to be analyzed
    unsigned int   FileLocation,        // the location of the trace in the file
    unsigned short RcdTraceLength,     // recorded trace length
    unsigned short *RcdTrace,          // recorded trace
    double         *slowfilter );      // slow filter response
```

```
// Add by Hongyi Wu
```

```
PIXIE16APP_EXPORT int PIXIE16APP_API HongyiWuPixie16ComputeSlowFiltersOffline (
    char          *FileName,           // the list mode data file name (with complete path)
    unsigned short ModuleNumber,      // the module whose events are to be analyzed
    unsigned short ChannelNumber,     // the channel whose events are to be analyzed
    unsigned int   FileLocation,      // the location of the trace in the file
    unsigned short RcdTraceLength,    // recorded trace length
    unsigned short *RcdTrace,         // recorded trace
    double        *slowfilter,        // slow filter response
    unsigned int   bl,
    double        sl,
    double        sg,
    double        tau,
    int           sfr,
    int           pointtbl );
```

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